



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification 5 : G01N 33/546, 35/02	A1	(11) International Publication Number: WO 93/20444 (43) International Publication Date: 14 October 1993 (14.10.93)
(21) International Application Number: PCT/US93/02776 (22) International Filing Date: 24 March 1993 (24.03.93) (30) Priority data: 07/859,218 27 March 1992 (27.03.92) US 07/916,737 20 July 1992 (20.07.92) US (71) Applicant: ABBOTT LABORATORIES [US/US]; CHAD 0377/AP6D-2, One Abbott Park Road, Abbott Park, IL 60064-3500 (US). (72) Inventors: RAYMOURE, William, J. ; 352 Briar Lane, Lake Luff, IL 60044 (US). CLARK, Frederic, L. ; 2712 Chamberlain Circle, Plano, TX 75023 (US). CLIFT, Gilbert ; 4514 Live Oak, Mesquite, TX 75150 (US). HENDRICK, Kendall, B. ; 1335 Forest Lane, Southlake, TX 76092 (US). KANEWSKE, William, J., III ; 1502 W. Colorado, Dallas, TX 75208 (US). LAGOCKI, Peter, A. ; 225 N. Hamilton Avenue, Park Ridge, IL 60068 (US). MARTIN, Richard, R. ; 8804 Saddlehorn #311, Irving, TX 75063 (US). MITCHELL, James, E. ; 184 River Road, Lake Barrington, IL 60010 (US). MOORE, Larry, E. ; 2713 Hunters Creek, Plano, TX 75075 (US). PENNINGTON, Charles, D. ; 980 Honey Lake Road, Lake Zurich, IL 60047 (US). WALKER, Edna, S. ; 3231 W.		Warner, Chicago, IL 60618 (US). SMITH, Jane, B. ; 26 Lindon Lane, Vernon Hills, IL 60061 (US). TAYL, Apparao ; 846 Langley Court, Grayslake, IL 60030 (US). VAUGHT, James, A. ; 908 Rosewood Court, Euless, TX 76039 (US). YOST, David, A. ; 19617 Selby Avenue, Poolesville, MD 20837 (US). (74) Agents: GORMAN, Edward, H., Jr. et al.; Abbott Laboratories, CHAD-0377/AP6D-2, One Abbott Park Road, Abbott Park, IL 60064-3500 (US). (81) Designated States: AU, CA, JP, KR, European patent (AT, BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE). Published <i>With international search report.</i>
(54) Title: METHODS FOR PROVIDING HOMOGENEOUS REAGENTS		
(57) Abstract A method for modifying a liquid assay reagent to provide prolonged homogeneity thereof, particularly where the liquid assay reagent comprises microparticles for performing a heterogeneous immunoassay, is provided wherein the addition of an inert material to a liquid assay reagent achieves neutral density to thereby prolong the homogeneity thereof for extended periods of time. A method for the automated agitation of assay reagents to maintain the homogeneity thereof with an automated, continuous and random access analytical instrument is also provided. The automated mixing is accomplished by a back and forth motion of a carousel onto which assay reagent containers or packs are mounted with asymmetric pauses which can be completed within a short period of time. The carousel acceleration, velocity, distance moved, and pause-asymmetry are optimized to provide rapid assay reagent resuspension without foaming or bubble formation. Accordingly, periodic removal of assay reagent packs by an operator in order to mix the reagents is not necessary.		

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	FR	France	MR	Mauritania
AU	Australia	GA	Gabon	MW	Malawi
BB	Barbados	GB	United Kingdom	NL	Netherlands
BE	Belgium	GN	Guinea	NO	Norway
BF	Burkina Faso	GR	Greece	NZ	New Zealand
BG	Bulgaria	HU	Hungary	PL	Poland
BJ	Benin	IE	Ireland	PT	Portugal
BR	Brazil	IT	Italy	RO	Romania
CA	Canada	JP	Japan	RU	Russian Federation
CF	Central African Republic	KP	Democratic People's Republic of Korea	SD	Sudan
CG	Congo	KR	Republic of Korea	SE	Sweden
CH	Switzerland	KZ	Kazakhstan	SK	Slovak Republic
CI	Côte d'Ivoire	LJ	Liechtenstein	SN	Senegal
CM	Cameroon	LK	Sri Lanka	SU	Soviet Union
CS	Czechoslovakia	LJ	Luxembourg	TD	Chad
CZ	Czech Republic	MC	Monaco	TG	Togo
DE	Germany	MG	Madagascar	UA	Ukraine
DK	Denmark	ML	Mali	US	United States of America
ES	Spain	MN	Mongolia	VN	Viet Nam
FI	Finland				

METHODS FOR PROVIDING HOMOGENEOUS REAGENTS

This application is a continuation-in-part application of U.S. Patent Application Serial No. 07/859,218, filed March 27, 1992, which is incorporated by reference herein.

Field of the Invention

The present invention relates to homogeneous reagents for the analysis of a test sample. In particular, the present invention relates to methods for maintaining the homogeneity of assay reagents and use thereof for the analysis of a test sample.

Background of the Invention

In order to meet the growing demands of the modern clinical laboratory to provide cost effective services, the use of automated clinical analyzers has increased. Reliance on automated clinical analyzers improves the efficiency of the laboratory procedures inasmuch as the technician has fewer tasks to performed. Automated clinical analyzers provide results much more rapidly while frequently avoiding operator or technician error, thus placing emphasis on accuracy and repeatability of a variety of tests.

Generally, analysis of a test sample involves the reaction of test samples with one or more reagents with respect to one or more analytes wherein it is frequently desired that the analysis be performed on a selective basis with respect to each test sample. Typically, such analysis involves forming a reaction mixture comprising the test sample and one or more assay reagents, and the reaction mixture is then analyzed by an apparatus as described above for one or more characteristics of the test sample.

Automated clinical analyzers are presently available for automatically performing such analyses of a test sample. Such analyzers typically include transport systems such as

conveyor systems and carousels designed to transport containers of sample liquids between various operating stations. For example, a reaction tube or cuvette containing a test sample may pass through a reagent filling station, mixing station, reaction forming station, detection stations, analysis stations, and the like. In particular, various automated immunoassay analyzers have been provided such as the Abbott IMx® analyzer and the Abbott TDx® analyzer (Abbott Laboratories, Abbott Park, Illinois, USA), often referred to as batch analyzers, which utilize procedures involving a variety of different assay steps which typically rely on detection and measurement of optical changes in a reaction mixture during the assay process. Random access analyzers have also been described which not only can analyze multiple test samples, but multiple analytes may be analyzed from each test sample. In addition, presently available sequential and random access analyzers include various reagents within the apparatus itself or placed near the apparatus for pipetting purposes. Liquid reagents, in bulk form, are selected for the various types of tests which are to be performed on the test sample, and are stored in or near the apparatus. The reagent delivery units, such as pumps and the like, along with valves, control and pipette mechanisms, are included in these automated analyzers so that different reagents can be mixed according to the type of test to be performed. Recently, apparatus and methods have been proposed for performing, selectively on the same sample, various homogeneous and heterogeneous assays concurrently in a random access fashion. Such apparatus and methods provide for the analysis of a plurality of liquid samples wherein each sample is analyzed with respect to at least one analyte utilizing both homogeneous and heterogeneous assay techniques.

Typically, such automated analyzers include assay reagent packs or containers from which assay reagents contained therein are removed during operation of the analyzer in order to carry out a particular analysis. Alternatively, a technician may be required to remove such assay reagents

from a container not associated with the analyzer for introduction thereof into the analyzer. In many instances, such assay reagents comprise one or more components which must be substantially homogeneous upon removal from the reagent container or pack in order to provide consistently accurate and reliable results when used for performing a particular assay. Since such assay reagents must be stored in assay reagent containers or packs for accessibility by a technician or the automated analyzer for extended periods of time, they tend to lack uniformity or homogeneity as a result of, for example, settling of assay reagent components, such as the settling of microparticles in an assay reagent for use in a heterogeneous immunoassay as described above.

In order to maintain homogeneity of such assay reagents, removal of the assay reagent container or pack and manual manipulation thereof by the operator, such as by inversion or agitation of the assay reagent container or pack, is required. Where the analysis of a test sample employing such assay reagents or portion thereof is performed manually, such as the addition of assay reagents prior to introduction into an automated analyzer, manual manipulation of an assay reagent container or pack by a technician is also nevertheless needed. In either instance, such manual manipulation is time consuming and cumbersome, and can lead to operator error, loss of assay reagents, and, when removal of an assay reagent container or pack from an automated analyzer is required, could result in damage to the automated analyzer. In addition, manual manipulation of an assay reagent container typically causes foaming and bubbles which result in inaccurate pipetting.

Summary of the Invention

According to the present invention, a method for modifying a liquid assay reagent comprising one or more particulate assay components to provide a homogeneous solution thereof is provided. According to such method, the

addition of an inert reagent to a liquid assay reagent has been found to achieve neutral density to thereby prolong the homogeneity thereof for from between about 1 hour and about 70 hours. In particular, such homogeneous liquid assay reagent is accomplished by providing an assay reagent comprising a liquid component and a particulate assay component, wherein the density of the liquid component and the density of the particulate assay component are substantially different. An inert reagent is added to the liquid assay reagent in an amount whereby the density of the liquid component and the density of the particulate assay component are substantially the same to thereby provide a substantially homogeneous liquid assay reagent having neutral density. According to a preferred embodiment, sucrose is added to a liquid assay reagent comprising microparticles for performing a heterogeneous immunoassay. Accordingly, the homogeneity of the microparticles in such liquid assay reagent can be maintained for a prolonged period of time between initial agitation of the assay reagent and subsequent agitation of the assay reagent. Means for automatically accomplishing such initial agitation and subsequent agitation with an automated, continuous and random access analytical instrument are also provided in order to maintain the homogeneity thereof, is also provided. In order to insure consistent, rapid resuspension and continued mixing of particulate assay reagents as well as non-particulate assay reagents with minimal operator involvement, the assay reagents can be mixed automatically each time a new assay reagent container or pack is introduced into the instrument, and periodically during instrument operation. The reagent container or pack is mounted onto a circular carousel which is capable of bidirectional rotation. The automatic mixing of assay reagents in the assay reagent container or pack can be accomplished by a back and forth motion of the carousel with asymmetric pauses between direction changes of such back and forth motion. The carousel acceleration, velocity, distance moved, and pause-asymmetry are optimized to

provide rapid assay reagent resuspension without foaming or bubble formation.

Automated assay reagent agitation eliminates the need of an operator to manually agitate the assay reagents which, for example, have been stored prior to their placement on the instrument. Such automated agitation allows the assay reagents to be loaded onto the instrument in less time and with less involvement of the operator. In addition, there is less tendency for reagents to foam or form bubbles, which is otherwise detrimental to instrument function and can negatively impact assay performance, with automatic mixing than with manual agitation such as inversion. Moreover, automated agitation insures that assay reagents are always mixed sufficiently and that they are mixed consistently. During operation of the instrument, not only can initial agitation be automatically accomplished to provide a homogeneous assay reagent as described herein prior to the use thereof, such as at the beginning of the day or after a period of time has elapsed during which the instrument has been non-operational or when a new reagent container or pack is introduced into the instrument, subsequent automatic mixing during instrument operation also maintains homogeneity of the assay reagents. Accordingly, periodic removal of assay reagent packs by an operator in order to mix the reagents is not necessary, and the period of time the reagent containers or packs can remain in the instrument without being removed for agitation is extended.

Brief Description of the Drawings

30

FIGURE 1 is an isometric view of an automated analytical system described herein illustrating the system cabinetry, exposed front end carousel, computer screen and keyboard.

35

FIGURE 2 is an isometric view of the automated analytical system apparatus frame and cabinet.

FIGURE 3 is a top plan view of the automated analytical system in section with component covers removed to show the automated analytical system apparatus in detail and relative position.

5 FIGURE 4 is a front elevational view of the automated analytical system in isolation and partial section of elements of the front end carousel.

10 FIGURES 4A and 4B represent a perspective side elevational view and partial end view of a reagent pack and reagent pack cover means for use with the automated analytical system.

FIGURE 5 is a top view in isolation and partial section of drive and guide elements of the front end carousel of the automated analytical system being removed.

15 FIGURE 6 is a cross-sectional side view of a process carousel of the automated analytical system in isolation with two reaction vessels in place, one of which is in position for an FPIA read.

20 FIGURE 7 is an isometric view of the probe, probe arm and pipettor of the automated analytical system in isolation.

FIGURE 8 is a schematic side view of the probe arm wiring and sensor means of the automated analytical system.

25 FIGURE 9 is a cross-sectional side elevational view of an automatic bubble flushing syringe apparatus of the automated analytical system.

FIGURE 9A is a sectional side view in isolation of the syringe bore end portion of the automatic bubble flushing syringe with the reciprocating piston near the end of travel toward the bore end portion.

30 FIGURE 9B is a sectional end view in isolation of the piston and bore of the automatic bubble flushing system syringe taken along line 9B-9D.

35 FIGURES 10 and 10A represent a top plan view of a reaction vessel and a side view of the reaction vessel for use with the automated analytical system, respectively, with reaction vessel compartments labeled where appropriate for FPIA processing.

FIGURES 10B and 10C present a top plan view and a side view of the reaction vessel, respectively, labeled and presented for MEIA processing.

FIGURE 11 is a sectional side view of the transfer element of the automated analytical system engaging a reaction vessel for transfer from the main carousel into the transfer station.

FIGURE 12 is a perspective side elevational view of the transfer station of the automated analytical system.

FIGURE 13 is a top plan view in section illustrating in isolation the controlled environment portion of the automated analytical system.

FIGURE 14 is a top plan view in section of the lower cabinet of FIGURES 1 and 2 illustrating water and/or buffer supply station as well as liquid and solid waste containers of the automated analytical system.

FIGURE 15 is a schematic view illustrating the system control environment airflow and temperature control of the automated analytical system.

FIGURE 16 is a side elevational view in partial section of a MEIA cartridge for use with the automated analytical system.

FIGURE 17 is a side elevational view in section of a MEIA cartridge feeder of the automated analytical system.

FIGURE 18 is a side sectional view in isolation of the MEIA cartridge feeder-cartridge orientation pin mechanism of the automated analytical system.

FIGURE 19 is a side sectional view in isolation of the MEIA cartridge ejector of the automated analytical system.

FIGURE 20 is a box diagram of the optics signal processor of the automated analytical system.

FIGURE 21 is a schematic of the FPIA optical system of the automated analytical system.

FIGURE 22 is a schematic of the FPIA read[er] sequence of the automated analytical system.

FIGURE 23 is a side sectional view in isolation of a MEIA cartridge carousel of the automated analytical system, MEIA cartridge and MEIA reader.

5 FIGURE 24 is a schematic of the MEIA system optical assembly of the automated analytical system.

FIGURE 25 is a schematic of the MEIA read sequence of the automated analytical system.

FIGURE 26 is a schematic reaction sequence of a FPIA for T4 performed on the automated analytical system.

10 FIGURE 27 is a schematic reaction sequence of a one-step sandwich MEIA performed on the automated analytical system.

FIGURE 28 is a schematic reaction sequence of a two-step sandwich MEIA performed on the automated
15 analytical system.

FIGURE 29 illustrates the comparison of centrifuged and non-non-centrifuged AFP microparticles.

FIGURE 30 illustrates an AFP microparticle accelerated settling study at various sucrose concentrations.

20 FIGURE 31 illustrates stability calibration curves for AFP microparticles.

FIGURE 32 illustrates the course of settling of AFP microparticles over a period of time.

25 FIGURE 33 illustrates the effect of reagent bottle size on the resuspension of a TSH microparticle reagent.

FIGURE 34 illustrates the effect of reagent bottle configuration on the resuspension of a TSH microparticle reagent.

30 FIGURE 35 illustrates the effect of reagent bottle fill volume on the resuspension of a TSH microparticle reagent.

FIGURE 36 illustrates the effect of sucrose concentration on the resuspension of a TSH microparticle reagent.

35 FIGURE 37 illustrates the effect of automatic agitation on the resuspension of a TSH microparticle reagent.

FIGURE 38 illustrates the effect of the range of automatic agitation on the resuspension of a TSH microparticle reagent.

Description of the Invention

Definitions

5 The following definitions are applicable to the present invention:

 The terms "homogeneity" and "homogeneous", as used herein to describe the physical characteristics of an assay reagent, refer to the uniformity of one or more assay
10 components substantially throughout the assay reagent (?). Such one or more components can be solid phase particulate materials, such as assay beads, assay particles, and assay microparticles for performing heterogeneous assays, which
15 tend to settle over a period of time, or liquid reagents having different densities which thereby result in separation thereof over a period of time.

 The term "spectrophotometric assay", as used herein, refers to the interaction in an assay solution between the analyte to be determined and a reagent system specific for the
20 analyte to produce a detectable change in the transmittive properties of the assay solution. The change in the transmittive properties refers to the amount of light absorbed or scattered by an assay solution within a particular wavelength band when a beam of light of known intensity is
25 passed through the assay solution. The change in the transmittive properties of an assay solution is measured by passing monochromatic light having a known intensity through the assay solution and determining the ratio of the intensity of the transmitted or scattered light to the intensity of the
30 incident light. Nearly all analytes either absorb energy of a specific wavelength or interact in an assay solution with a particular reagent system to produce a detectable change in the transmittive properties of the assay solution, characteristics which have resulted in the development of
35 numerous specific spectrophotometric assays. Spectrophotometric assays which rely upon the measurement of the change in the transmittive properties of an assay

solution as a measure of an analyte in the assay solution include, for example, assays wherein there is a change in the color of the assay when there is a change in the turbidity of the assay solution, that is, turbidimetric or nephelometric assays.

The term "colorimetric assay", as used herein, refers to the change in the transmittive properties of an assay solution which is generally referred to as the absorbance of the assay solution and is dependent upon the change in the color of the assay solution due to the interaction of the analyte to be determined and reagent system specific for the analyte. The absorbance of the assay solution is related to the concentration of the analyte in the assay solution. A colorimetric assay utilizes a chromogenic reagent system capable of interacting in an assay solution with the particular analyte of interest, to produce a detectable change in the transmittive properties, specifically the color, of the assay solution. Numerous chromogenic reagent systems useful in the determination of specific analytes have been developed and are commercially available.

The term "turbidimetric assay", as used herein, refers to the determination of the amount of light scattered or blocked by particulate matter as light passes through an assay solution. The analyte of interest interacts with a reagent system specific for the analyte to form a suspension of particles in the assay solution. As a beam of light having a known intensity is passed through an assay solution, the suspension of particles formed by the interaction of the analyte reagent system blocks or scatters the incident light, thereby reducing the intensity of the light transmitted through the assay solution. The change of the transmittive properties in a turbidimetric assay refers to the decrease in the intensity of the light transmitted through an assay solution and is related to the amount of incident light that is scattered or blocked by the suspension of particles, and depends upon the number of particles present and the cross-sectional area of such particles.

The term "nephelometric assay", as used herein, is similar to a turbidimetric assay in that the analyte of interest interacts with a reagent system specific for the ligand to form a suspension of particles in the assay solution.

5 The change in the transmittive properties of the assay solution is also related to the amount of incident light scattered or blocked by the suspension of particles. Unlike a turbidimetric assay wherein the intensity of the light transmitted through the assay solution is measured, the
10 scattered or blocked light is measured at an angle to the light incident to the assay solution. Therefore, in a nephelometric assay the change in the transmittive properties refers to the difference in intensities of light incident to the assay solution and light scattered at an angle to the incident light.

15 The term "fluorometric assay", as used herein, refers to the determination of an analyte in an assay solution which is chemically or immunologically transformed into a fluorescent complex or conjugate thereby producing a detectable change in the fluorescent properties of the assay solution. The change
20 in the fluorescent properties of the assay solution is measured by exciting the fluorescent complex or conjugate properties produced with monochromatic light of a wavelength within the excitation wavelength band of the fluorescer, and measuring the intensity of the emitted light at a wavelength
25 within the emission wavelength band of the fluorescer. The fluorescent intensity of the emitted light is related to the concentration of the analyte. However, the intensity of the fluorescence emitted by the assay solution may be inhibited when the analyte to be determined complexes with
30 nonfluorescent interferences such as protein or phosphates present in the sample, or when the sample containing the ligand to be determined has sufficient color so as to act as a filter and thereby reduce the intensity of the emitted fluorescence. It is well recognized that in order to maximize
35 the sensitivity and specificity of a fluorometric assay, these inhibiting factors, if present, must be overcome either by removal of the nonfluorescent interferences or color

producing material prior to the analysis, or by compensating for the presence of such factors using an internal standard added to a second aliquot of sample and carrying out the entire assay procedure using the aliquot containing the internal
5 standard.

The term "homogeneous immunoassays", as used herein, refers to a competitive immunoassay format involving a competition between an analyte from a test sample and a tracer for a limited number of receptor binding sites on an
10 antibody to the analyte. The tracer comprises the analyte or analog thereof labeled with a detectable moiety wherein the concentration of analyte in the test sample determines the amount of the tracer that will specifically bind to the antibody. The amount of the tracer-antibody conjugate
15 produced by such binding may be quantitatively measured and is inversely proportional to the amount of analyte present in the test sample. For example, fluorescent polarization techniques for making such determination, such as in fluorescent polarization immunoassays as described herein,
20 are based on the principle that a fluorescently labeled compound when excited by linearly polarized light will emit fluorescence having a degree of polarization inversely related to its rate of rotation. When a molecule such as a tracer-antibody conjugate having a fluorescent label is excited with
25 a linearly polarized fluorescent molecule it is constrained from rotating between the time light is absorbed and emitted. When a "free" tracer molecule (i.e., unbound to an antibody) is excited by linearly polarized light, its rotation is much faster than the corresponding tracer-antibody conjugate and the
30 molecules are more randomly orientated, therefore, the emitted light is polarized. Accordingly, when plane polarized light is passed through a solution containing the aforementioned reagents, a fluorescent polarization response is detected and correlated to the amount of analyte present in
35 the test sample. Various fluorescent compounds which can be employed for performing fluorescent polarization assays include, but are not intended to be limited to,

aminofluoresceins, such as described in U.S. Patent No. 4,510,251 and U.S. Patent No. 4,614,823, incorporated herein by reference; triazinylaminofluoresceins, such as described in U.S. Patent No. 4,420,568 and U.S. Patent No. 4,593,089,
5 incorporated herein by reference; carboxyfluoresceins, such as described in U.S. Patent No. 4,668,640, incorporated herein by reference; and the like.

The term "heterogenous immunoassays", as used herein, refers to immunoassay formats involving a labeled reagent or
10 tracer comprising an analyte, an analog of the analyte, or an antibody thereto, labeled with a detectable moiety, to form a free species and a bound species. In order to correlate the amount of tracer in one of such species to the amount of analyte present in the test sample, the free species must first
15 be separated from the bound species, which can be accomplished according to methods known in the art employing solid phase materials for the direct immobilization of one of the binding participants in the binding reaction, such as the antibody, analyte or analog of the analyte, wherein one
20 of the binding participants is immobilized on a solid phase material, such as a test tube, beads, particles, microparticles or the matrix of a fibrous material, and the like, according to methods known in the art. Heterogenous immunoassays can be performed in a competitive immunoassay format wherein, for
25 example, the antibody can be immobilized to a solid phase material whereby upon separation, the amount of the tracer which is bound to such solid phase material can be detected and correlated to the amount of analyte present in the test sample. Another form of a heterogeneous immunoassay
30 employing a solid phase material is referred to as a sandwich immunoassay, which involves contacting a test sample containing, for example, an antigen with a protein such as an antibody or another substance capable of binding the antigen, and which is immobilized on a solid phase material. The solid
35 phase material typically is treated with a second antigen or antibody which has been labeled with a detectable moiety. The second antigen or antibody then becomes bound to the

corresponding antigen or antibody on the solid phase material and, following one or more washing steps to remove any unbound material, an indicator material such as a chromogenic substance which reacts with the detectable moiety (e.g.,
5 where the detectable moiety is an enzyme, a substrate for such enzyme is added) to produce a color change. The color change is then detected and correlated to the amount of antigen or antibody present in the test sample. For example, the teachings of the present invention can be employed in a
10 heterogeneous immunoassay which can be performed by the automated analytical system described herein, in either a competitive or sandwich immunoassay format, or in a microparticle capture enzyme immunoassay, such as that described in Clinical Chemistry, Volume 34, No. 9, pages
15 1726-1732 (1988), employing microparticles as the solid phase material.

The term "test sample", as used herein, refers to a material suspected of containing the analyte. The test sample can be used directly as obtained from the source or following
20 a pretreatment to modify the character of the sample. The test sample can be derived from any biological source, such as a physiological fluid, including, blood, saliva, ocular lens fluid, cerebral spinal fluid, sweat, urine, milk, ascites fluid, raucous, synovial fluid, peritoneal fluid, amniotic fluid or the
25 like. The test sample can be pretreated prior to use, such as preparing plasma from blood, diluting viscous fluids, or the like; methods of treatment can involve filtration, distillation, concentration, inactivation of interfering components, and the addition of reagents. Besides physiological fluids, other liquid
30 samples can be used such as water, food products and the like for the performance of environmental or food production assays. In addition, a solid material suspected of containing the analyte can be used as the test sample. In some instances it may be beneficial to modify a solid test sample to form a
35 liquid medium or to release the analyte.

The term "analyte" or "analyte of interest", as used herein, refers to the compound or composition to be detected

or measured and which has at least one epitope or binding site. The analyte can be any substance for which there exists a naturally occurring binding member or for which a binding member can be prepared. Analytes include, but are not limited to, toxins, organic compounds, proteins, peptides, microorganisms, amino acids, nucleic acids, hormones, steroids, vitamins, drugs (including those administered for therapeutic purposes as well as those administered for illicit purposes), virus particles and metabolites of or antibodies to any of the above substances. In particular, such analytes include, but are not intended to be limited to, ferritin; creatinine kinase MIB (CK-MB); digoxin; phenytoin; phenobarbital; carbamazepine; vancomycin; gentamycin; theophylline; valproic acid; quinidine; leutinizing hormone (LH); follicle stimulating hormone (FSH); estradiol, progesterone; IgE antibodies; vitamin B2 micro-globulin; glycated hemoglobin (Gly. Hb); cortisol; digitoxin; N-acetylprocainamide (NAPA); procainamide; antibodies to rubella, such as rubella-IgG and rubella-IgM; antibodies to toxoplasmosis, such as toxoplasmosis IgG (Toxo-IgG) and toxoplasmosis IgM (Toxo-IgM); testosterone; salicylates; acetaminophen; hepatitis B virus surface antigen (HBsAg); antibodies to hepatitis B core antigen, such as anti-hepatitis B core antigen IgG and IgM (Anti-HBC); human immune deficiency virus 1 and 2 (HIV 1 and 2); human T-cell leukemia virus 1 and 2 (HTLV); hepatitis B e antigen (HB_eAg); antibodies to hepatitis B e antigen (Anti-HB_e); thyroid stimulating hormone (TSH); thyroxine (T₄); total triiodothyronine (Total T₃); free triiodothyronine (Free T₃); carcinoembryonic antigen (CEA); and alpha fetal protein (AFP). Drugs of abuse and controlled substances include, but are not intended to be limited to, amphetamine; methamphetamine; barbiturates such as amobarbital, secobarbital, pentobarbital, phenobarbital, and barbital; benzodiazepines such as librium and valium; cannabinoids such as hashish and marijuana; cocaine; fentanyl; LSD; methaqualone; opiates such as heroin, morphine, codeine, hydromorphone, hydrocodone, methadone, oxycodone,

oxymorphone and opium; phencyclidine; and propoxyhene. The term "analyte" also includes any antigenic substances, haptens, antibodies, macromolecules and combinations thereof.

5 The term "analyte-analog", as used herein, refers to a substance which cross-reacts with an analyte-specific binding member, although it may do so to a greater or lesser extent than does the analyte itself. The analyte-analog can include a modified analyte as well as a fragmented or
10 synthetic portion of the analyte molecule, so long as the analyte-analog has at least one epitopic site in common with the analyte of interest. An example of an analyte-analog is a synthetic peptide sequence which duplicates at least one epitope of the whole-molecule analyte so that the analyte-
15 analog can bind to an analyte-specific binding member.

 The term "binding member", as used herein, refers to a member of a binding pair, i.e., two different molecules wherein one of the molecules specifically binds to the second molecule through chemical or physical means. In addition to
20 antigen and antibody binding pair members, other binding pairs include, as examples without limitation, biotin and avidin, carbohydrates and lectins, complementary nucleotide sequences, complementary peptide sequences, effector and receptor molecules, enzyme cofactors and enzymes, enzyme
25 inhibitors and enzymes, a peptide sequence and an antibody specific for the sequence or the entire protein, polymeric acids and bases, dyes and protein binders, peptides and specific protein binders (e.g., ribonuclease, S-peptide and ribonuclease S-protein), and the like. Furthermore, binding
30 pairs can include members that are analogs of the original binding member, for example, an analyte-analog or a binding member made by recombinant techniques or molecular engineering. If the binding member is an immunoreactant it can be, for example, a monoclonal or polyclonal antibody, a
35 recombinant protein or recombinant antibody, a chimeric antibody, a mixture(s) or fragment(s) of the foregoing, as well as a preparation of such antibodies, peptides and nucleotides

for which suitability for use as binding members is well known to those skilled in the art.

The term "detectable moiety", as used herein, refers to any compound or conventional detectable chemical group having a detectable physical or chemical property and which can be used to label a binding member to form a conjugate therewith. Such detectable chemical group can be, but is not intended to be limited to, enzymatically active groups such as enzymes, enzyme substrates, prosthetic groups or coenzymes; spin labels; fluorescers and fluorogens; chromophores and chromogens; luminescers such as chemiluminescers and bioluminescers; specifically bindable ligands such as biotin and avidin; electroactive species; radioisotopes; toxins; drugs; haptens; DNA; RNA; polysaccharides; polypeptides; liposomes; colored particles and colored microparticles; and the like.

The term "continuous access", as used herein, refers to the ability to add additional test samples or reagents to the automated analytical system described herein without the interruption of assays which are being performed by the automated analytical system at the time of such addition.

The term "random access", as used herein, refers to the ability of the automated analytical system described herein to simultaneously perform more than one scheduled assay in any order in which such plurality of scheduled assays are presented into the automated analytical system described herein.

The term "simultaneous", as used herein, refers to the ability of the automated analytical system described herein to independently perform two or more scheduled assays at the same time.

The term "kitting", as used herein, refers to the ability of the automated analytical system described herein to create a unit dose disposable by separately transferring test samples and reagents to a reaction vessel described herein without initiation of an assay reaction sequence.

The term "quat" refers to a polycationic material solution for assays.

The term "flexible protocols" refers to the variety of different assay protocols capable of being processed in accordance with the inventive system. Examples include MEIA formats configured in 1-step and 2-step sandwich and competitive assay formats; order of activity processing, including the ability to initiate sample processing for both MEIA formats and FPIA formats on the front-end carousel prior to transfer onto the process carousel; variable incubation periods; optical read formats and wash sequences. This contrasts to some prior art, known random access systems which force all assay protocols to adhere to a strict "lock-step" format, in which assay configuration (i.e. 1-step versus 2-step formats), activity order, incubation timing, and other similar protocols are fixed by the instrument.

Methods

According to the present invention, methods for maintaining the homogeneity of a liquid assay reagent are provided. According to one embodiment of the present invention, the addition of an inert reagent to a liquid assay reagent comprising one or more particulate assay components such as beads, particles, microparticles, and the like, has been found to achieve neutral density of such particulate assay components. Such inert reagents contemplated by the present invention include, but are not intended to be limited to, sucrose, metrizamide, metrioic acid, and the like.

Neutral density of the assay reagent as described herein is achieved by determining the optimum concentration of the inert reagent which will eliminate the settling of the particulate materials of the liquid assay reagent. When the density of the liquid portion of the liquid assay reagent and the particulate materials are equivalent, the particulate materials will be in a suspended state with substantially little or no settling thereof wherein the homogeneity of such liquid assay reagent can be maintained for an extended period of time between, for example, initial agitation of the liquid

assay reagent and subsequent agitation of the liquid assay reagent.

The concentration of the inert material required to achieve neutral density will depend upon the particular assay components in the liquid assay reagent. In addition, the concentration of the inert material will also depend upon the assay kinetics of a particular assay being performed employing such particulate assay reagents. For example, in a heterogeneous immunoassay employing a solid phase assay reagent comprising a protein, such as an analyte or antibody thereto, immobilized to a particle or microparticle, the concentration of the inert material will depend upon the protein to microparticle ratio. Where such heterogeneous immunoassay employing such solid phase assay reagent is performed in either a sandwich format or a competitive format, the concentration of the inert material in the solid phase assay reagent in each of such formats will vary because the protein:particle ratio in the sandwich format will be higher than the protein:particle ratio in the competitive format.

According to a preferred embodiment of the present invention, sucrose is employed to provide a homogeneous microparticle assay reagent according to the present invention for use in a microparticle enzyme immunoassay (MEIA) as described above for the determination of various analytes.

Analytical Systems

According to the present invention, various known assay techniques and formats described herein, employing various analytical apparatus described herein, as well as other analytical apparatus known in the art, can be performed. Such analytical apparatus employ various detection systems including, but not limited to, spectrophotometric absorbance assays such as end-point reaction analysis and rate of reaction analysis, turbidimetric assays, nephelometric

assays, radiative energy attenuation assays (such as those described in U.S. Patent No. 4,496,293 and U.S. Patent No. 4,743,561 and incorporated herein by reference), ion capture assays, colorimetric assays, fluorometric assays, electrochemical detection systems, potentiometric detection systems, amperometric detection systems, and immunoassays. Immunoassays include, but are not intended to be limited to, heterogeneous immunoassays such as competitive immunoassays, sandwich immunoassays, immunometric immunoassays, and the like, where the amount of a detectable moiety employed therein can be measured and correlated to the amount of analyte present in a test sample.

When performing a heterogeneous immunoassay as described herein, separation of the bound and free species is preferably accomplished by capture of the microparticles on a glass fiber matrix of an MEIA cartridge, a process that relies on the high affinity of glass fibers for the microparticles, wherein the microparticles adhere to the surface of the fibers irreversibly, and nonspecifically bound material can be effectively removed by washing the matrix. The matrix also provides a precisely located mechanical support for the microparticles during the optical quantification phase of the assay protocol as described herein. In particular, microparticles coated with antibody to the analyte in the test sample are incubated with the test sample containing the analyte of interest to form a capture complex with the analyte from the test sample. A conjugate comprising antibody to the analyte labeled with a detectable moiety, preferably an enzyme, is then incubated with the capture complex to form the second of a sandwich complex. When performing a competitive immunoassay, microparticles coated with antibody to the analyte in the test sample are incubated with the test sample containing the analyte of interest and a conjugate comprising the analyte or analog thereof labeled with a detectable moiety, preferably an enzyme. Removal of unbound conjugate is accomplished with the glass fiber matrix of an MEIA cartridge and, where the detectable moiety is an

enzyme, a substrate for the enzyme capable of providing a detectable signal is added and the signal provided thereby is measured and correlated to the amount of analyte present in the test sample. Preferably, the enzyme-substrate system employed by the competitive and sandwich MEIA formats is alkaline phosphatase and 4-methylumbelliferyl phosphate (MUP), although other enzyme-substrate systems known in the art can be employed as well.

The MEIA cartridge comprises a reaction well for retaining and immobilizing microparticle-analyte complexes. The reaction well has an entrance port and means for holding a quantity of sample and assay reaction mixtures positioned over a fibrous matrix which retains and immobilizes microparticle-analyte complexes as described above. The fibrous matrix is composed of fibers having an average spatial separation greater than the average diameter of the microparticles. Preferably, the average fiber spatial separation is greater than 10 microns. The reaction well further comprises an absorbant material positioned below the fibrous matrix to enhance the flow of sample and assay reaction mixtures through the fibrous matrix. Preferably, the absorbant material is a fibrous material whose fibers predominantly lie in a plane perpendicular to the lower surface of the fibrous matrix. The absorbant material is in fluid communication with the fibrous matrix. Generally, the absorbant material is in physical contact with the lower surface of the fibrous matrix. The interior of the reaction well, therefore, is generally sized or contains positioning means to maintain the fluid communication between the absorbant material and the fibrous matrix. Preferably, a spike located at the bottom of the reaction well can be used to force the absorbant material into contact with the lower surface of the fibrous matrix. Additionally, it is preferable to vent to the atmosphere the gases displaced in the absorbant material by the liquids absorbed therein during the performance of an immunoassay.

The teachings of the present invention are particularly useful when performing the various assay techniques described herein on a continuous and random access analytical system apparatus as described below and as shown in the
5 Figures hereof.

Such instrument comprises a front end carousel assembly inclusive of a sample cup carousel, an assay reagent pack carousel and a reaction vessel carousel mounted
10 concentrically and serviced by a transfer pipetting means suitable for kitting and/or mixing reagents with a test sample. The reagent pack carousel includes a plurality of reagent packs containing one or more containers which are capable of holding assay reagents for performing various assays described herein and as known in the art. The kitted
15 and pipetted reaction vessels are transferred through a transfer station which provides means for transferring the kitted and pipetted reaction vessels to a processing work station which includes a controlled environment for maintaining temperature and provides timing for mixing of
20 reagents and incubation. At least two assay procedural apparatus are provided which are scheduled for the various samples and kitted reagents in a unit dose disposable means for analyzing the incubated reaction mixtures. The unit dose disposable reaction vessels are removed from the process
25 carousel by operation of the transfer station, which includes means for removing the disposable reaction vessel from the system. According to such analytical system apparatus, a system scheduler generates and optimizes the workload for the system's mechanical resources from all the tests ordered
30 to run on the system. The main goal of the scheduler is to keep the system's resources from sitting idle while there are tests remaining to be processed by the system. Keeping each of the resources busy also minimizes the time required by the instrument to perform the tests.

35 In particular, initial agitation and subsequent agitation of assay reagents contained in the reagent pack is accomplished by a back and forth motion of the assay reagent

pack carousel with asymmetric pauses which can be completed within a short period of time. The carousel acceleration, velocity, distance moved, and pause-asymmetry are optimized to provide rapid assay reagent resuspension without foaming or bubble formation. In order to insure consistent, rapid resuspension and continued mixing of particulate assay reagents as well as non-particulate assay reagents with minimal operator involvement, the assay reagents can be mixed automatically each time a new reagent pack is added to the reagent carousel, and periodically during instrument operation. It is to be understood that such automatic agitation is particularly useful where a particulate assay reagent which has been modified according to the teachings of the present invention, such as the addition of sucrose to a microparticle assay reagent for performing an MEIA as described above, is employed.

As would be understood by one skilled in the art, the speed of rotation of the carousel will be dependent upon a number of considerations such as, for example, the shape of the assay reagent container; the fill volume of the assay reagent container, preferably at fill volumes of about one-half fill maximum or less; the distance, assay reagent density, acceleration and final velocity of movement; and the duration of pauses between successive movements, which can be determined visually or with an optical instrument. For example, when employing the continuous and random access analytical system apparatus described herein, the maximum acceleration of the carousel is from between about 7,500 steps per second and about 8,500 steps per second, preferably about 7,850 steps per second, and the velocity of the carousel is from between about 3,500 steps per second per second and about 4,500 steps per second per second, preferably about 4,000 steps per second per second. Although agitation of assay reagents can be accomplished with slower accelerations and velocities, such agitation will of course occur at a slower rate. The carousel movement is from between about 25 steps of carousel movement and about 150 or more steps of carousel

movement, with from between about 25 steps and about 75 steps in alternating directions, preferably about 50 steps in alternating directions. It is to be understood that the distance traveled in opposition directions does not have to be equal, and the pauses between successive movements can be of varying lengths of time which can be controlled by, for example, either timers or loops (x = 1 to y , next x) in the software. Preferably, such pauses are not equal, known as assymetric pauses, between direction changes to provide efficacious mixing. Where the alternating directions comprise about 50 steps in each direction, assymetric pauses are preferably from between about 75 steps and about 100 steps.

A high-level view of the scheduling process, including the scheduling of such back and forth motion of the assay reagent pack carousel, can be broken into two steps: (1) proper scheduling of each of the activities in a test is ensured before the test is kitted, and (2) an attempt to perform each test activity prior to its original scheduled execution time, to minimize resource idle time and increase test throughput in the system. To enable scheduling a test in advance of its performance in the system, each test's assay protocol contains several timing parameters used in the scheduling process. Each activity of the test contains time values which the scheduler uses to determine which resources the activity requires and the time period that these resources are needed. Also, each activity in the test can be tied to other activities by incubation periods. These incubation periods, which are dictated by the chemistry of the assay, help the scheduler determine the amount of time that must elapse between the execution of two activities. Each incubation period in the assay protocol provides for the minimum and maximum time that may elapse between the execution of each activity. These limits are referred to in the scheduling process as the incubation window for the activities.

When operating such analytical apparatus system, the operator chooses the order that tests are prepared to run on the instrument by selecting the placement of samples on the

instrument. The sample placed closest to the pipette station is the first sample prepared to run on the instrument. To guard against evaporation, a test will not be prepared until the scheduler ensures that all resources used by the test's activities will be available at the required times set forth in the test's assay protocol. Preparation of a particular test will be postponed whenever an activity of another test already in the instrument has a resource scheduled at the time it is needed by an activity on that test. The sample preparation area of the instrument will remain idle until the test can be scheduled without conflicting with tests already in the instrument. When proper scheduling of the test can be achieved, the test will be prepared and transferred into the process area.

The second step in the scheduling process is to optimize the workload for each system resource to minimize both the resource's idle time and the time required to perform the resource's workload. Once tests are transferred into the process area, the scheduler optimizes the existing schedule for each resource. At predetermined intervals, the scheduler examines the next interval of work for each resource. If there is any idle time in this interval, the scheduler attempts to minimize the idle time by rearranging the resource's workload to eliminate idle time, providing the activities remain within their allowed incubation windows. When optimization of this interval is complete, this section of the workload is performed by the resource at the designated times. The scheduler continues to prepare samples as long as there are samples on the instrument that have tests ordered to be run. optimization of the resources' workloads will continue until all tests transferred into the system have finished processing. The analytical apparatus system described herein allows special priority handling of specific samples identified by the user as being stat samples. A stat sample is a sample that must be processed by the instrument in the shortest amount of time possible. Special handling of stat samples occurs both in

the front sample entry area and in the processing area of the instrument.

When performing a stat procedure, the operator chooses the order that tests are prepared to run on the instrument by selecting the placement of samples on the instrument. The sample placed closest to the pipette station is the first sample prepared to run on the instrument. This pattern of sample preparation is interrupted whenever the user places a stat test on the instrument. Whenever a stat test is ordered, the system will finish preparing the test on the current sample, and then move directly to the stat sample to prepare all its tests. To guard against evaporation, sample preparation will not begin for a test before proper scheduling of the test's activities in the processing area is ensured. The system scheduling algorithm is also modified for stat processing. The scheduling algorithm used for normal tests attempts to maximize the number of tests processed in the instrument each hour. This occurs by allowing sufficient time between test activities to enable other tests' activities to be performed in these gaps. The scheduling approach used for stat tests attempts to process this one test in the shortest amount of time possible. Each activity of a stat test is scheduled at the earliest possible time of execution as defined in the test's assay definition. When all activities of a test are guaranteed proper scheduling in the instrument, sample preparation of the test will begin. After all tests on the stat sample are prepared, the system will return to the sample it was working on before it serviced the stat.

Stat tests receive special consideration in the processing area when there is idle time in a resource's workload. At predetermined intervals, the scheduler examines the next interval of work allocated to each resource in the processing area of the system. If there is any idle time during this interval, the scheduler attempts to minimize it by rearranging the resource's workload. Test activities scheduled for this resource that can be performed earlier than they are currently scheduled, as defined by their assay protocols, are

moved forward to fill the idle time. Stat test activities are the first candidates to be pulled forward in the workload, thus further decreasing the amount of time needed to process the stat test in the instrument. The system stat test handling algorithms have been shown to allow stat tests to be processed in the minimum amounts of time possible, without having a negative effect on the instrument's overall throughput of tests per hour.

According to the immunoassay methodologies described above, standard solutions of the analyte of known concentrations covering the clinical concentration range are typically prepared and assayed as is the test sample to be assayed. This blank assay provides a series of signal measurements corresponding to the known concentrations from which a standard curve is drawn. The optical signal corresponding to the unknown sample is correlated in a concentration value through interpretation from the blank or standard curve.

Automated analytical methodology for effecting analysis of a plurality of test samples according to analytical apparatus system described herein is achieved by introducing reagent packs, test sample container and reaction vessels onto concentric carousels of a main carousel. The test sample container can be a test tube, cuvette, vacutainer tube, and the like, for holding a test sample. The reagent packs and test sample containers are identified and aligned respectively with a reaction vessel for transfer and kitting of the reaction vessel by transfer of test sample and specific reagents from the reagent pack for preparation of a predetermined test. The reaction vessel containing the test sample and one or more reagents is transferred to a process carousel wherein controlled environment conditions exist for incubation once the sample has been appropriately mixed with various reagents to form a reaction mixture. When all assay processing steps have been completed, the reaction mixture is identified and transferred to at least, for example, one of a fluorescent polarization immunoassay reader or a

microparticle enzyme immunoassay cartridge positioned on a separate cartridge wheel or carousel for further preparation before reading. The processed test samples are read and the readings are calculated with the resulting data being recorded and/or printed.

The methodology of the automated immunoassay analytical system is achieved through the use of a self-contained, fully automated, continuous and random access instrument comprising a main carousel assembly consisting of the reagent pack carousel, a reaction vessel carousel and a test sample container carousel concentrically and independently rotatable. The main carousel assembly is provided with a transfer pipette operated by a boom arm for transferring and kitting test sample and reagents into the reaction vessel automatically following a predetermined test schedule. The main carousel assembly is provided with bar code readers for reagent packs and test sample containers and has the capability of aligning the reagent pack carousel and test sample container carousel and a reaction vessel for pipette transfer operations. Once the assay to be performed is scheduled, the reaction vessel carousel, the reagent pack carousel and the test sample container carousel are rotated until the reaction vessel, a reagent pack and a test sample container, respectively, are determined to be in the transfer pipette access position. The transfer pipette then transfers the test sample from the test sample container and, depending upon the assay to be performed, the reagents from the reagent pack are transferred to the reaction vessel. The reaction vessel carousel is then rotated to a transfer station position which contacts the reaction vessel with a transfer mechanism and pulls the reaction vessel into the transfer station. The reaction vessel is then loaded onto the process carousel by the transfer mechanism.

When performing a fluorescent polarization immunoassay (FPIA) with the automated analytical system, various pipetting activities are performed by a second transfer pipette apparatus which is in service for the process

carousel, and the process carousel is rotated so that the reaction vessel, when properly pipetted with, for example, FPIA reagents, is at the read station of the FPIA processing stations and the FPIA determination on reading, is made on the reaction vessel. The process carousel is then rotated so that the read reaction vessel is at the transfer station. The reaction vessel is again contacted and transferred by the transfer station. The transfer station is rotated and pushes the reaction vessel into a release container opening.

For a microparticle enzyme immunoassay (MEIA) performed with the automated analytical system described herein, after the various pipetting activities for the MEIA, which can be completed at the main carousel assembly, the reaction vessel is transferred to the process carousel as described in the FPIA process. Pipetting can also be accomplished in the process carousel or jointly between the two carousels. To complete the MEIA, the reaction mixture is transferred from the reaction vessel to a matrix of an MEIA cartridge on a cartridge carousel with the second transfer pipette. The matrix is washed with a buffer and a substrate, such as MUP (defined earlier), or other suitable substrate known in the art. The cartridge carousel is then rotated so that the MEIA cartridge is positioned at an MEIA processing assembly and the MEIA determination is made. The MEIA reaction vessel is ejected into the waste container as described for the FPIA reaction vessel. The MEIA cartridge is independently ejected from the cartridge wheel by an ejector at an appropriate ejector station into a waste container.

Preferably, two distinct analytical technologies as described above, FPIA and MEIA, are incorporated into the automated analytical system described herein and with which the method of the present invention can be employed. However, more than two distinct analytical technologies can be incorporated into the analytical system. These methods are complimentary and share a commonality of apparatus and procedural steps, with the FPIA generally being the method of choice for analytes of low molecular weight and MEIA for

molecules such as protein hormones, antibodies or analytes of low molecular weight requiring higher sensitivity. The two technologies share system components including the operator control panel, pipetting boom assemblies, fluidic systems, air and liquid reagent heaters, printers, bar code reader and stepper motors. Such commonality of use of system components allows for a compact instrument despite the dual FPIA and MEIA capability.

The FPIA optic systems (such as described in U.S. Patent No. 4,269,511 and incorporated herein by reference) can utilize a polarizing filter which is an electrically switched liquid crystal, maintaining a compact size and avoiding complex and potentially unreliable moving parts. When performing FPIA assays utilizing the automated analytical system described herein, the FPIA reagent packs will typically include a tracer comprising the analyte or analog thereof, coupled to a detectable moiety, an antibody specific to that analyte, and a specimen pretreatment reagent. In a preferred FPIA format, the analyte being determined competes with the tracer for a limited number of binding sites on the antibodies specific to the portion or portions of the analyte and tracer. The detectable moiety component of the tracer is preferably a fluorescent moiety selected from the group consisting of fluoresceins, aminofluoresceins, carboxyfluoresceins, fluoresceinamines, and the like, more preferably carboxymethyl-aminomethyl-fluorescein, carboxyethylaminomethyl-carboxyfluorescein, 6-carboxyfluorescein, 5-carboxyfluorescein, succinylanimomethyl-fluorescein, thiourea-aminofluorescein, methoxytrianolylaminofluorescein, aminofluorescein, and the like.

MEIA results can be determined by quantifying the rate of fluorescence developed when fluorogenic substrate is converted by the action of an enzyme labeled conjugate. For example, when performing either a competitive MEIA or sandwich MEIA, the specifically bound alkaline phosphatase on the microparticles is detected by addition of the fluorogenic

substrate MUP to the matrix. The alkaline phosphatase catalyzes hydrolysis of the MUP to inorganic phosphate and fluorescent 4-methylumbelliferone (4-MU). The liberated 4-mu is detected by the MEIA optics assembly front surface
5 fluorometer which is designed to detect fluorescence of low concentrations of 4-MU without interference by fluorescence of 4-MUP at a wavelength of 367. A system of lenses and optical filters focus filtered light (wavelength = 365) from a mercury arc lamp on to the surface of the matrix and focus
10 emitted fluorescence from 4-MU (wavelength = 448) on to a photo multiplier tube. Like the FPIA optics assembly, the MEIA optics system is compact and has no moving parts. About five percent of the excitation light is detected by a photodiode, allowing normalization of the fluorescence data and
15 generation of a control signal used by the lamp power supply to maintain the intensity of the excitation light within five percent over the useful life of the lamp. The MEIA post-processor uses linear regression analysis to convert the data from multiple successive determinations of 4-MU
20 fluorescence to a rate which is proportional to the concentration of alkaline phosphatase conjugate specifically bound to the microparticles.

MEIA formats can be run with a multi-position MEIA auxiliary carousel and process carousel as well as a MEIA
25 reagent pack containing a homogeneous microparticle assay reagent according to the present invention, an alkaline phosphatase conjugate and, in some cases, a dilute buffer specific for the assay being performed. Since the microparticles of the homogeneous microparticle assay
30 reagent remain suspended during the course of the assay to provide a homogeneous solution thereof, they can accurately and reliably pipetted. The effective surface area of polystyrene latex microparticles is several fold greater than that of a large diameter polystyrene bead (e.g., one quarter
35 inch beads) commonly used in commercial immunoassays. Because of this large surface area and the very small diffusion distance between analyte and the capture molecules

on the surface of the microparticles, the capture phase employed in many of the MEIA methods being performed reaches equilibrium within several minutes, allowing for a full carousel of test samples to be completed in a very short time frame.

Unlike an FPIA, the heterogeneous immunoassays, such as a MEIA, require a separation step as described above. In particular, after incubation of the microparticles with a test sample, the microparticles are separated from the reaction mixture by transfer to the matrix contained in the MEIA cartridge as described above. The matrix provides a precisely located mechanical support for the microparticles during the subsequent optical read phase of the assay. This precisely located mechanical support, i.e. the cartridge, is fit into the auxiliary carousel at a predetermined spacing from the reader apparatus by camming means.

Referring to the drawings, FIGURES 1 and 2 present isometric views of an automatic immunoassay analytical system apparatus with which the teachings of the present invention are particularly useful. It is to be understood that the automated immunoassay analytical system described herein is presented only with those components of primary interest with respect to the assay cuvette of the present invention. The drawings do not illustrate all of the mechanical and electrical elements for driving and controlling the various components of the system, wherein an of such omitted elements may have various known forms which can be readily realized by one of ordinary skill in the art having knowledge of the information provided herein with regard to the mode of operation of the system and the various components and related processes utilized for treating samples and determining analytical results.

The system apparatus as it appears in FIGURE 1 presents the system apparatus as used by the technician, with FIGURE 2 illustrating an isometric view of the frame and cabinetry with component parts removed. The system apparatus described herein is identified generally by the reference numeral 2 in

FIGURE 1. The system apparatus 2 has an exposed front end carousel 4 which is serviced by a first transfer pipette mechanism 6 for kitting scheduled tests along with samples into a reaction vessel. The system provides a computer screen 8 and computer keyboard 10 along with access panels 12 for accessing storage and waste compartments. The system apparatus 2 is provided with rollers 14 for movement of the system apparatus within a laboratory complex as required. The freedom of movement of the system apparatus through rollers 14 is allowed since the system is fully self-contained but for power requirements.

In FIGURE 2, the system apparatus 2 cabinet frame 16 is illustrated with substantially all functioning components of the system apparatus removed. A controlled environment zone 18 is a closed unit during operation with light shielding and rigid control of airflow as well as temperature as opposed to the open front end carousel 4. The front end carousel 4 communicates with the controlled environment zone 18 through a transfer port 20. The front end carousel 4 is mounted to an aluminum base plate which rests on a support platform 22 and the first transfer pipette mechanism is mounted on means 24.

The top plan view in section of FIGURE 3 presents the functioning component system apparatus in some detail with relative positioning of the system apparatus to further illustrate the process flow of the system apparatus. For example, sample cups 26 are mounted on a sample cup carousel 28 which is concentrically fitted within the front end carousel 4 along with reagent pack carousel 32 and reaction vessel carousel 36. The reagent pack carousel 32 is concentrically fitted between the sample cup carousel 28 and the reaction vessel carousel 36. The reagent pack carousel carries reagent packs 30 and the reaction vessel carousel 36 carries reaction vessels 34. The front end carousel 4 has an operable bar code reader 38 for automatically identifying reagent pack carousel 32 and sample carousel 28. A wash cup 40 is provided for the first transfer pipette mechanism 6 for

washing as required between transfer of various sample and reagents. The first transfer pipette mechanism 6 is utilized in kitting the various reagent pack liquid materials and sample into a reaction vessel 34. The reagents and the sample are properly kitted through means of the first transfer pipette mechanism 6 inclusive of pump means. The various carousels are rotated and aligned for kitting at the pipetting station. The kitted reaction vessel 34 is positioned by reaction vessel carousel 36 into the proper position for transfer to the transfer station 42. The reaction vessel 34 is transferred to the transfer station 42 through transfer means wherein the transfer station 42 is then rotated to move the reaction vessel onto process carousel 46. As shown, the process carousel is driven by a stepper motor 48 and is serviced by a second transfer pipette mechanism 50. Both the FPIA and MEIA procedures utilize the system apparatus commonly up through and including the process carousel 46. The process carousel 46 includes FPIA processing 52 and FPIA processing lamp 54 for direct reading of FPIA analysis of kitted, pipetted and properly reacted reagents sample from the reaction vessel 34. The controlled environmental zone 18, which includes the transfer station 42 and process carousel 46, provides FPIA processing with air circulation under temperature control by cabinet air circulation fan 56. A wash cup 58 for the second transfer pipette mechanism 50 is provided. The second transfer pipette 50 is utilized for adding reagents (pipetting) under conditions of incubation and timing to the sample in the FPIA test schedule reaction vessel 34 for FPIA processing. MEIA processing can also utilize the second transfer pipette 50 for adding reagents to the sample before the reaction mix is added to MEIA cartridges 68 which are mounted on the cartridge wheel carousel 64. The transfer of the MEIA reagent mixed sample to the MEIA cartridge 68 is by the function of the second transfer pipette 50. A motor 60 drives the cartridge wheel 64. The cartridge wheel 64 is provided with MEIA cartridges 68 through the operation of a cartridge hopper 66 which automatically feeds and positions the MEIA

cartridges 68 onto the cartridge wheel 64. The process area includes the second transfer pipette mechanism 50 and heater/pump 44. The cartridge wheel carousel 64 is further serviced by a MEIA buffer heater and dispenser 70, MUP heater and dispenser probe 72, and MEIA reader 74. The MEIA
5 cartridges 68 are removed from the cartridge wheel 64 by a cartridge ejector 62 after the MEIA read has been completed.

It is to be understood that the utilization of the first transfer pipette mechanism 6 and the second transfer pipette
10 mechanism 50 as described herein provide a safety mechanism to ensure that test samples and reagents are pipetted to thereby prevent false negative results in the event there are incorrect amounts of the respective sample and reagents for a particular assay.

15 Approaching the operable elements of the system apparatus in greater detail, FIGURE 4 provides a front elevational view in isolation and partial section of elements of the front end carousel 4. FIGURES 4A and 4B illustrate a reagent pack with a cover means 31 which is opened and
20 closed pivoting along axis 37. A return notched drive arm 35 is utilized to open and close the cover means 31 by contact with the cover contact surface 33.

FIGURE 5 provides a top view in isolation and partial section of elements of the drive and guide systems of the
25 main carousel 4 with the various carousels removed. In FIGURE 5 a sample cup carousel stepper motor 76 is shown mounted with mounting spring 78. The reagent pack carousel motor 80 is also shown with a mounting spring 82. The reaction vessel carousel motor 84 and mounting spring 86 are positioned to
30 the exterior of the two inner carousels, i.e. the sample cups carousel 28 and the reagent pack carousel 32. Roller guides 88 are provided for the sample cup carousel 28 and a tensioning spring 90. The reagent pack carousel is provided with roller guides 92 and tensioning means 94. The reaction vessel roller
35 guides 96 are also provided with spring elements 98, the purposes of the guide and these various spring elements being

to maintain very finite tracking of the concentric carousels when motivated by the individual stepper motors.

The front end carousel 4 inclusive of the three front end carousels, the sample cup carousel 28, reagent pack carousel 5 32 and reaction vessel carousel 36 can by example contain the following capacities. The sample cup carousel 28 can hold 60 blood collection tubes, such as Vacutainer® blood collection tubes, or 90 sample cups which are injection molded as one piece and can be provided with standalone base mounts. 10 Standalone base mounts are suitable for technician handling and pipetting of samples into the sample cups. The reagent pack carousel 32 provides for 20 different reagent packs 30. The reaction vessel carousel 36 provides 90 reaction vessels 34.

15 The process carousel 46 as shown in FIGURE 6 is an isolational cross-sectional side view. One reaction vessel 34 is at rest or nonoperative position and a second reaction vessel 34 is in position for FPIA read. The process carousel 46 is capable of bidirectional motion for timely movement of the 20 various reaction vessels 34 to pipettor action, read, or transfer to and from the carousel. Up to about 36 or more reaction vessels 34 can be processed at one time on the process carousel 46 depending on diameter and sizing of the reaction vessels 34.

25 The first transfer pipette mechanism 6 of FIGURE 7 includes a transfer pipette Z axis motor 102 which moves the probe arm 104, probe 106 and probe tip 108 in a vertical direction while transfer pipette R axis motor 100 drives the probe arm 104, probe adjustment means 106 and probe tip 108 30 in a horizontal motion. The first transfer pipette mechanism 6, sometimes labeled "Sample Probe Arm Mechanism", moves the probe between the sample cup 26, the reagent pack 30, the reaction vessel 34 and the wash cup 40. The wash cup 40 is used to wash the interior and exterior surfaces of the first 35 transfer pipettor mechanism 6 probe. The drive of the first transfer pipette mechanism is a rack-and-pinion drive means along the Z and R axis by two- stepper motor drivers. A brake

is provided to hold the Z axis position when power is lost, thus avoiding damage to the system apparatus. For example, the first transfer pipette mechanism can be designed to have a Z axis travel of about 3 inches and an R axis travel of about 11-1/2 inches.

The first transfer pipette mechanism 6 and the second transfer pipette mechanism 50 are closely related in general system apparatus function and design, with variation on travel and size being the only substantial differences. Both units have a probe arm circuit 110 as illustrated by the schematic side view of FIGURE 8. The schematic illustrates the R axis motor 100 and the Z axis motor 102 in relationship to an upper PCB 112 and a R axis home sensor 114. A lower PCB 116 is illustrated in relationship to the Z axis home sensor 118 with a coil cable 120 connecting the various elements.

Various elements of syringe 122 which provides automatic bubble flushing and fluids to the various pipetting mechanisms is provided in various views in FIGURES 9, 9A and 9B. The ability of diagnostic instrumentation to accurately perform an assay is critically dependent on the precision and accuracy with which syringes, i.e. pipetting, can aspirate and dispense reagents and samples. The precision and accuracy of a syringe is severely degraded by the presence of small air bubbles inside a syringe. Bubbles, unfortunately, are all too common and are difficult to remove or avoid. Syringe 122 avoids these problems by automatically flushing bubbles completely out of the fluidics system. The syringe 122 is configured such that a piston 124 reciprocates through a seal 126 and into a close-fitting bore 128. The end of the bore 130 is closed. The piston 124 has a piston end 132 which approximates the geometry of the closed bore end 130. Two ports to the bore are 180° apart and are located near the seal and are comprised of a fluid entry port 134 and a fluid exit port 136. An annulus 138 exists between the piston 124 and bore 128. Pressurized line diluent is introduced to the fluid entry port 134. The fluid flows out into the annulus 138 around both sides of the piston 124 and then into the fluid exit

port 136. This crossflow flushes bubbles from the area near the seal. While the crossflow is occurring, the piston 124 is reciprocated inside the bore 128. This reciprocation causes high fluid flow velocities in the annulus 138 between the piston 124 and the bore 128. The high flow velocity dislodges any bubbles that may be adhering to the piston 124 or bore wall. The inward stroke of the piston 124 pushes these dislodged bubbles across the crossflow area where they are swept out of the syringe. The piston end 132 and the bore end 130 have similar spherical shapes. When the piston 124 strokes to its full inward extension, it comes very close to the bore end 130. Any bubble that may be stuck on the bore end 130 is disrupted and dislodged. Likewise, when the piston strokes to its full outward extension, its end is flush with the seal 126. The sequence of reciprocating the piston while crossflowing can be automatically executed any time by the system apparatus.

Once the fluid leaves the fluid exit port 136 of the syringe 122, it must travel through a tube fitting, through a length of tubing, through another tube fitting, into a probe 106 and out the probe tip 108. It is at the probe tip 108 that the aspirating and dispensing of reagents actually occurs. Any bubbles trapped between the syringe and the probe tip will also degrade performance, so there must be no place for the bubbles flushed out of the syringe to lodge. It is therefore necessary to use zero dead volume tubing fittings on the tubing between the syringe and the probe.

The reaction vessel 34 is discussed in detail relative to either the MEIA scheduling or the FPIA scheduling in FIGURES 10, 10A, 10B and 10C. FIGURES 10 and 10A present the FPIA kitting utilization. The reaction vessel is illustrated in both the top plan view (FIGURE 10) and the side view (FIGURE 10A). S reagent antiserum is deposited in well 142 while T reagent tracer is deposited in well 144 with P reagent popper being deposited in well 146. Wells 150 and 152 can serve for providing a variety of reagents, buffers and/or dilution liquids to the apparatus. The sample is deposited in well 148 and

predilution liquid in well 154. The utilization of the transfer pipettor in depositing the required reagents into a reaction vessel along with the sample is called kitting. The depositing of the various required reagents and the like into a single
5 reaction vessel along with a sample to be analyzed is called pipetting.

The MEIA reaction vessel as shown in top and side views of FIGURES 10B and 10C, respectively, contains prediluent in well 156; microparticle materials being deposited in well
10 158; conjugate directly in the reaction well 166; assay diluent in well 162; and the sample in well 164. The buffer well is 168 and predilution well is 170. Once kitting is complete, many of the subsequent FPIA and MEIA pipetting steps can be performed either in the main carousel or in the
15 process carousel utilizing the pipetting mechanisms of both carousels. This is possible because the kitted reaction vessel, once kitted, is transferred immediately into the transfer station and thus into the process carousel which exists in a controlled temperature environment .

The transfer station 42 plays a key role in apparatus and process function. In FIGURE 11, a sectional side view of the transfer element of the transfer station 42 is shown engaging reaction vessel 34 by means of a reaction vessel transfer projection 172. The transfer arm 173 is projected out
25 between reaction vessel elements of the reaction vessel carousel 36 and, by rotation of the transfer station 42, engages the reaction vessel transfer projection 172. By means of a transfer arm drive gear 174, the transfer arm 173 rack gear 176 moves the transfer arm 173 out and in relationship
30 to the transfer station 42. The transfer station 42 has a rotation axis 178. In FIGURE 11A, a reaction vessel is shown in phantom as would be mounted on the front end carousel 4, reaction vessel carousel 36 engaged by the transfer arm 173 by means of reaction vessel transfer projection 172. The
35 reaction vessel 34 in FIGURE 11 is illustrated onboard the transfer station by reaction transfer station 42 moves the reaction vessel 34 between the front end carousel 4 and the

process carousel 46. The transfer station 42 moves the discarded reaction vessel 34 from the process carousel 46 to the waste ejection station (not shown). The transfer station 42 is driven by a stepper motor drive and is supported by precision linear ball bearings and axis of rotation ball bearings.

The process carousel 46 holds, for example, 36 reaction vessels 34 and has a carousel diameter of about 12.5 inches. The process carousel 46 moves the reaction vessel 34 between the transfer station 42, the second transfer pipettor mechanism 50, the point of pipetting, and the FPIA reader processing 52. The process carousel 46 is driven by a stepper motor and supported by three wheels for height control and control of any radial movement caused by irregularly shaped carousel elements.

The second transfer pipette mechanism 50 moves the pipette probe between the wells in the reaction vessel 34 on the process carousel 46 to the MEIA cartridge 68 on the auxiliary carousel 64 and to the wash cup 58. A rack-and-pinion drive through two axis stepper motor drives achieves precision drive on both the R and Z axis. Travel, for example, on the Z axis can be about 3 inches and on the R axis about 4.5 to 5.0 inches.

The auxiliary carousel 64 holds, for example, 32 MEIA cartridges 68 and has a diameter of about 9.5 inches. The auxiliary carouser 64 moves the MEIA cartridges 68 between various stations including the second transfer pipettor mechanism pipette point, the MUP dispense station 72, the MEIA washstation 70 and the MEIA reader 74 and the MEIA cartridge ejection point 62. The auxiliary carousel 64 is stepper motor driven and is carried by three wheels with one wheel located at the Z axis height control at the cartridge insertion point, the second wheel at the pipette point, and the third wheel at the MEIA reader in order to maintain the auxiliary carousel 64 within desired geometric relationships to these various functions.

MEIA cartridges 68 are loaded into a cartridge hopper 66 which feeds the MEIA cartridges 68 into the auxiliary carousel 64. The automatic feeding of the MEIA cartridges 68 is provided with a proper height adjustment of the cartridge 68 into the auxiliary carousel 64 as required by MEIA reading. The cartridge hopper 66 feeds individual cartridges 68 to the auxiliary carousel 64 and changes the axis of orientation of the cartridge 68 from horizontal to vertical by automatic means. Removal of the MEIA cartridges 68 is achieved through the use of an ejector 62 which operates through an ejection rod and forces the MEIA cartridge 68 from the auxiliary carousel 64 which is dropped into a solid waste container.

Buffer supply stations are presented in FIGURE 14 which is a top plan view in section of the apparatus showing the cabinet frame 16, front end carousel 4 in partial phantom and a power supply element 192 along with diluent system or buffer pressurization means 194. A supply bottle 196 is also mounted in the lower cabinet of frame 16 as well as solid waste 198 and liquid waste 200 containers for receiving processed liquids and solid waste.

A schematic view illustrating the environmental airflow and temperature control system is shown in FIGURE 15 wherein make up air 204 enters and hot air exits at exhaust 206. Airflow 202 is indicated by arrows and the controlled environmental airflow schematic 214 is provided with at least one heater element 208 and fan element 210. At least one temperature sensor 212 is provided for control of the air temperature and can be correlated with the airflow 202 control.

The MEIA cartridge 68 is shown in a side elevational view in FIGURE 16. The MEIA cartridge 68 has a funnel throat 216 and a cartridge opening 218. The MEIA cartridge 68 contains support matrix material 222.

A MEIA cartridge 68 and cartridge hopper 66 are shown in a side elevational view in FIGURE 17. The MEIA cartridges are positioned horizontally in the cartridge hopper 66 and are manipulated from the bottom of the V-shaped cartridge hopper

66 one-by-one through a cartridge shuttle 222. The cartridge feeder has a cartridge cam block 224 and a cartridge orientation shoot 226 which functions through cartridge orientation pin 228 and cartridge orientation pin 230 for providing the MEIA cartridge 68 in vertical alignment= for insertion into the auxiliary carousel 64. The orientation pins 228 and 230 are illustrated in FIGURE 18 which is a side sectional view in isolation of the MEIA cartridge feeder cartridge orientation mechanism. The MEIA cartridge 68 is shown in an enlarged view in FIGURE 18 as being engaged and disengaged by cartridge orientation pin 228 and cartridge orientation pin 230. The cartridge orientation pin 230 is shown in engagement position at position 232 against the base 236 of the MEIA cartridge 68 while cartridge orientation pin 228 is shown in engagement position 234 of the cartridge funnel throat portion 216. Upon withdrawal of these pins from the engaging positions, the MEIA cartridge 68 is released from the bottom portion first, i.e. the withdrawal of cartridge orientation pin 230, thus allowing the bottom of a cartridge 68 to drop by gravity before the top of the cartridge is released which is engaged by cartridge orientation pin 228 in the cartridge funnel throat 216. The rounded or semicircular holding surfaces of the orientation pin allow the release of the bottom of the MEIA cartridge and the rolloff of the funnel throat portion 216 from the cartridge orientation pin 228. The vertically aligned MEIA cartridge 68 is then inserted into the auxiliary carousel 64 to a controlled height by the action of an insertion cam means 227 as shown in FIGURE 17.

A side view of a MEIA cartridge ejector 62 is illustrated in FIGURE 19. The cartridge ejector 62 functions through an ejector rod 240 and can be driven by manual or automatic drive means 242. The ejected MEIA cartridge is ejected through an ejection passage to the solid waste 198 container.

A box diagram of the optics signal processor of the apparatus is provided in FIGURE 20 wherein the signal from the FPIA optics 248 is fed to a DSP A/D 250 which also sends serial bus signal 252 from an optic signal processor 8-bit

microcontroller 254. The controller 254 is connected to computer elements through 256. Signal from the MEIA optics 258 are fed into a DSP A/D element 260 which also sends serial bus 262 from the controller 254. Signal is fed to the FPIA optics through 264 from high voltage power supply 266 and serial bus 268 which is in communication between the microcontroller 254 and the optics power supply board 270A. The FPIA tungsten lamp power supply FPIA 270 is in electronic communication with the FPIA optics 272. Signal is sent to the MEIA optics through 274 from high voltage power supply 276 which is in communication through serial bus 268 to the microcontroller 254 and mercury lamp power supply MEIA 280. The MEIA mercury lamp power supply 280 is also in electronic communication with MEIA optics through 282.

A schematic view of the FPIA optical system 284 is shown in FIGURE 21. The FPIA optical system 284 has a tungsten halogen source lamp 286 which focuses light through a heat reflector 288, an aperture 290 and heat absorber 292 to a lens 293 for introduction into an excitation filter 294. The light energy is then contacted with a beam splitter 296 which presents part of the beam to a polarizer 298 and liquid crystal 300. The light continues into another lens 301 before being focused on the cuvette 140 containing the FPIA reaction mixture. Light is emitted from the cuvette through lens means 303 before entering an emission filter 302. The reflected light from the emission filter 302 passes through a polarizer 304 before going to a focusing lens 306 and being focused for feed into photo multiplier tube 308. The beam splitter 296 splits out part of the light from the original source through lens 310 into a reference detector 312 which, in turn, controls the tungsten halogen source lamp.

A schematic view of the FPIA read sequence 314 is presented in FIGURE 22. The FPIA read sequence 314 has a pre-read time 316 divided into carousel move time 318 and carousel settle time 320. Sub-read interval 340 is divided into a horizontal sub-read 342, A/D converter settle time 344, and a liquid crystal activation time 346. A vertical sub-

read interval is identified by 348 which is inclusive of A/D converter settle time 350. Liquid crystal relaxation time is indicated by 352. The liquid crystal relaxation time 352 is illustrated in a preread time sequence. High voltage settle
5 time 324 is further illustrated by lamp settle time 326 that shows the lamps in a sinner 328 and full burn 330 activation. Activities of the FPIA read sequence 314 provide for activities where scheduling windows 332 as exemplified by read prep 334, read parameter 336 during which the lamps are
10 at full burn, and collection results 338 during the lamp settlement time and liquid crystal relaxation time 352.

FIGURE 24 is a schematic view of the MEIA system optical assembly 364. An MEIA light source is provided by mercury source lamp 364 which passes light through an
15 excitation filter 362 to a filter reflector 360 before being fed through lens 358 into MEIA cartridge 68. Reflected fluorescent light is fed back through the filter 360 to a photomultiplier tube 374 after passing through a wide band-pass emission filter 370 and narrow band-pass emission filter 372. Part of
20 the light energy from the mercury source lamp 364 passes directly through filter 360 to a bandpass filter 368 before influencing the photo diode 366.

An MEIA read sequence schematic is presented in FIGURE 25 wherein the MEIA read sequence 376 has a preread time
25 378 inclusive of carousel move time 380 and carousel settle time 382. High voltage settle time is indicated by graph 384 which is coincident with the lamp settlement time 386 showing lamp simmer 388 and lamp full burn 390. MEIA read sequence 376 has activities with scheduling windows 392
30 inclusive of read prep 394, read parameter 396 and collection results 398. The actual MEIA read sequence 376 is inclusive of sub-read interval 400 having a sub-read 402 and a dwell time 404. Another segment of the MEIA read sequence 376 is indicated by sub-read interval 406 inclusive of sub-read
35 number to 408 and dwell time 410 with additional sub-reads 412 as indicated by number 3 through (N-1) and partial sub-

read interval 414 inclusive of sub-read number N-416. The next possible preread time is indicated by 418.

In addition to the continuous and random access analytical system apparatus described herein, immunoassay formats described above can be performed on the Abbott IMx® analyzer and the Abbott TDx® analyzer (Abbott Laboratories, Abbott Park, Illinois, USA) according to the teachings of the present invention. The Abbott IMx® analyzer utilizes MEIA technology for high and low molecular weight analytes requiring greater sensitivity, and FPIA technology, such as that used on the Abbott TDx® analyzer, is used primarily for lower molecular weight analytes. A front surface fluorometer is used to quantify a fluorescent product generated in the MEIA assays, while a fluorescence polarization optical system is used to quantify the degree of tracer binding to antibody in the FPIA assays. The test samples are automatically processed by a robotic arm with a pipetting probe and a rotating carousel which positions the samples for processing and permit the analysis of multiple samples and provides for access to the test samples for the formation of subsequent reaction mixtures. In particular, test samples are automatically processed with a robotic arm comprising a pipetting probe and a rotating carousel which positions the samples for processing. The assay reagents for performing MEIA and FPIA procedures are stored in a stationary reagent pack from which the pipetting probe removes the appropriate assay reagents for performance of a particular assay. The reagent pack typically includes a plurality of containers which separately contain various assay reagents such as, for example, an antibody reagent, a labeled reagent, buffer, diluent, and the like, for performing MEIA and FPIA methodologies.

Other assay formats, such as homogeneous assays, the detection of precipitate formed by reaction between antigens and antibodies in a test sample-cell to form light scattering centers, and methods and apparatus for detecting immunological agglutination reactions known in the art can

also be performed according to the teachings of the present invention. Such apparatus and methods include, for example, the steps of measuring light absorption of the liquid medium with antibody before and after the antigen-antibody reaction
5 by using light which is absorbable by the antibody, and calculating the difference of the absorptions. In this way, the presence or absence of agglutination can be detected based on the fact that the agglutination reaction reduces the concentration of antibody, which affects the light absorption
10 of the liquid medium. As is typical of methods and apparatus for performing homogeneous assays, these procedures do not require separation of a solid phase from the reaction mixture for further analysis. Spectrophotometric assays can also be performed according to the teachings of the present invention
15 on the Abbott Spectrum clinical analyzer and the Abbott Spectrum Series II clinical analyzer (Abbott Laboratories, Abbott Park, IL, USA). In addition, turbidimetric and nephelometric assays according to the teachings of the present invention can be utilized in the analysis of blood,
20 urine, spinal fluid, and the like, for the determination of analytes such as proteins wherein there is no comparable colorimetric assay due to the lack of an effective chromogenic reagent system. Yoe and Klimman, Photoelectric Chemical Analysis, Vol. II: Nephelometry, Wiley & Sons, Inc., New York,
25 1929, describe various nephelometric assays. various reagents and reagent systems which can be employed for performing spectrophotometric assays on the automated analytical systems described herein include, but are not intended to be limited to, those for the simultaneous determination of
30 glucose and urea, such as described in U.S. Patent No. 5,037,738 and incorporated herein by reference. The simultaneous determination of calcium and phosphorous; the simultaneous determination of cholesterol and triglycerides; determining isoenzymes; determining blood ammonia levels,
35 and the like, can be also be performed according to the teachings of the present invention.

It is to be understood that the method of automatically agitating an assay reagent container with the continuous and random access analytical system apparatus described herein is not limited thereto, and that such automatic agitation can be accomplished with other automated instruments known in the art having a carousel or the like on which assay reagent containers or packs are mounted or having automated mechanisms for the agitation of assay reagents thereon. It is to be also understood that although automated instruments for performing various assay formats and techniques have been described herein, such assay techniques and formats or portions thereof can be performed manually in accordance with the teachings of the present invention. When such assays are performed manually, various means for agitation of a homogeneous liquid assay reagent according to the present invention can be agitated with various automated agitation means known in the art such as by rotation, rocking, contacting with a Vortex® mixing apparatus, and the like.

The present invention will now be illustrated, but is not intended to be limited by, the following examples.

Example 1

Preparation of Homogeneous Microparticle Assay Reagent

The optimum sucrose concentration of a microparticle assay reagent necessary to maintain an even suspension of microparticles was determined employing an Abbott IMx analyzer as follows:

(1) The effect of centrifugation on the original microparticle solution:

- a. Pipet 1.5 mL original microparticle solution into an eppendorf centrifuge tube.
- b. Centrifuge in TDx centrifuge until a pellet of microparticles and a clear supernatant are obtained (usually 2-3 minutes, time varies for each assay). If a clear supernatant is

not obtained after 5 minutes (TSH for example), it may be necessary to dilute the original microparticles 1:1 before centrifugation:

- i) pipet 0.75 mL original microparticles +
5 0.75 mL IMx MEIA buffer into the eppendorf centrifuge tube
- ii) centrifuge until a clear supernatant is obtained (usually 2-3 minutes)
- iii) carefully pipet off supernatant, repeat step i) and ii)
- 10 c. Resuspend microparticles.
- d. Perform an assay calibration with the original microparticles and the centrifuged microparticles. (Place the eppendorf tubes in reagent bottles containing water. Secure the tubes with parafilm. HVR tables
15 may have to be altered if level sense errors are encountered.)
- e. Compare calibrator curves (i.e., curve shapes, avg. rates, and A-F span;. Table 1, Table 2 and Figure 29.

(2) Prepare 1 liter of assay microparticle buffer solution
20 (excluding the sucrose).

(3) Prepare the test microparticle buffer solutions with varying concentrations of sucrose. The sucrose concentrations to test will vary for each assay. A range of sucrose concentration
25 should be checked. (for example, where AFP sucrose concentration is 13.6%, test solutions are prepared at 13%, 14%, 15%, 16%, and 17% sucrose).

(4) Determine the optimal sucrose concentration from an
30 accelerated settling experiment:

- a. Pipet 1.5 mL of the original microparticles into an eppendorf centrifuge tube.
- b. Centrifuge and pellet out the microparticles
35 with an Abbott TDx analyzer centrifuge.
- c. Remove the original supernatant.

- d. Replace the original supernatant with the test microparticle buffer solutions containing various concentrations of sucrose.
- e. Resuspend.
- 5 f. Centrifuge for 31 minutes /1400g at 2-8 C (equivalent to approximately 31 days real time).
- g. By visual observation select the tube with the lowest concentration of sucrose where the microparticles have remained evenly suspended, (not floating) which is the optimal sucrose concentration (Table 3).
- 10 h. If no tube has achieved neutral density, steps 3 and 4 are repeated using the process of elimination to find the optimal sucrose concentration.
- 15 (5) Determine the optimal sucrose concentration with several different lots of microparticles. Because the density of the microparticle reagent is dependant on the amount of protein bound per microparticle, some lot to lot variation in the optimum sucrose concentration may be seen.
- 20

Example 2

Assay Performance Of Homogeneous Microparticle Assay Reagent

25 Assay performance with (1) the optimal sucrose concentration and (2) current sucrose concentration was evaluated as follows:

30 (1) Prepare microparticle diluent which contains original sucrose concentration and optimal sucrose concentration (as determined in Example 1). Prepare according to original diluent MF. Adjustment to sucrose will have to be made to achieve new concentration, optimal sucrose concentration.

35 (2) Obtain concentrated microparticle solution and also dilution required to prepare working microparticle reagent. Obtain master-lotted conjugate for that lot of microparticles and use throughout study.

(3) Prepare x mL of optimum sucrose microparticles, and current sucrose microparticles with new supernatant (where x = normal # of mL in a full reagent pack.

5 (4) Perform a calibration run and evaluate precision. Execute 3 runs of a full carousel of low, medium, and high controls on:

- 10 a. microparticles with the current sucrose concentration;
b. microparticles with the new optimal sucrose concentration;

15 (5) Repeat steps 2 - 4 with the second lot of microparticles.

(6) For each run calculate the mean, standard deviation, and %CV (Table 4, Table 5 and Figure 30).

20

Example 3

Settling Time Of Homogeneous Microparticle Assay Reagent

An 18 day real time stability study was performed as follows:

25

	<u># Kits</u>	<u>Sucrose Concentration</u>	<u>Storage Temp.</u>	<u>Inversion Before Run</u>
	1	Original	2-8°C	yes
30	1	Original	RT	no
	1	Optimal	2-8°C	yes
	1	Optimal	2-8°C	no
	1	Optimal	RT	yes
35	1	Optimal	RT	no

[Store all reagent bottles capped; recommended time points are at 0, 1, 2, 4, 8, 10, 15, 17, and 21 days]

- (1) Prepare sufficient microparticle reagent, original sucrose concentration and optimal sucrose concentration, such that a total of 3 kits can be prepared with original sucrose and 5 kits can be prepared with optimal sucrose as described above.
- 5 (2) At time 0, perform a calibration curve for the original and optimal sucrose reagent condition (curve #1 & curve #2).
- 10 (3) At each time point run the A&D calibrator in duplicate and control in repetitions of 2. Run the above for each test condition.
- 15 (4) All runs except the calibration runs, should be in mode 2 (If the assay is not normal run as mode 2, edit assay parameter #102 RUN MODE to 1 and select mode 2 on the front panel before pressing the RUN key.)
- 20 (5) All runs (excluding the calibration run) for each condition should be from only one reagent pack (total # of test per reagent pack = 90 after the last testpoint). The calibration runs are done on a separate reagent pack to provide:
- a. 1 reagent pack for calibration of original sucrose concentration;
 - 25 b. 1 reagent pack for calibration of optimal sucrose concentration;
 - c. 1 reagent pack for each storage test condition (total of 6, 2 original sucrose and 4 optimal sucrose)
- 30 (6) Calculate the means, std. dev., within run CVs and between run CVs for each test condition.

35 The results listed below were obtained with a different experimental design and are only used as an example of data already obtained.

(1) Centrifugation and replacement of the original supernatant should not significantly affect assay performance. If for some reason, centrifugation and the new supernatant do affect assay performance, (HCG for example) continue on with the accelerated settling experiment. (step D) to find the optimum sucrose concentration. In this case, for the evaluation of assay performance, sucrose may be added directly to the original microparticles and brought up to the optimum sucrose concentration found in step D (it is to be understood that this may result in a slight decrease in rates since the addition of sucrose directly to the original microparticles increases the total reagent volume, thus slightly diluting the microparticle concentration).

(2) The optimal sucrose concentration determined during the accelerated settling experiment will be:

- a. Be visually evenly suspended after the 31 min centrifugation at 1400 g.
- b. Have comparable calibration curves as the control runs.
- c. Have comparable assay performance (precision) as the control runs (compare mean, std, dev. and within run CVs)

(3) For the Real Time stability study:

a. Rates for the Non-Inverted original microparticles may decrease as the microparticles settle over time and between run CVs will be higher than the normally inverted microparticles. It is to be understood that not all assays will be significantly affected by microparticle settling throughout this stability study. In the examples shown, AFP rates were not significantly affected; Figure 31).

b. The rates for the original inverted microparticles are not expected to change significantly with respect to time, and should have acceptable within run and between run CVs.

c. The rates for the optimal sucrose microparticles are expected to behave similarly as the original inverted microparticles.

(4) Various sucrose concentrations in TSH, HCG, T3, and Toxo-IgG are illustrated in Table 6 and a summary of precision for TSH, HCG, T3 and ToxoligG is shown in Table 7 (Figures 32 and 33).

Table 1
A F P - CENTRIFUGATION VS. NON-CENTRIFUGATION

IMx ASSAY PRINTOUTS: (Example-AFP CAL, NO CENTRIFUGATION)

DATE: 0/00/00

TIME: 0:00:00

TECH ID : _____

RGNT LOT : _____

SERIAL # : 000

ASSAY 65 AFP

LOC	CALIB	RATE
1	A	14.2
1	A	13.7
2	B	125.8
2	B	128.1
3	C	337.9
3	C	337.9
4	D	526.8
4	D	550.2
5	E	841.7
5	E	850.4
6	F	1123.3
6	F	1142.8

CALIB	ng/mL	AVGR
A	0	13.9
B	15	127.0
C	50	337.9
D	100	538.5
E	200	846.1
F	350	1133.1

CALIBRATION ACCEPTED

LOC	ID	ng/mL	RATE NOTE
13	_____	18.71	149.3
14	_____	19.04	151.3
15	_____	83.15	470.9
16	_____	80.26	459.3
17	_____	194.59	829.4
18	_____	196.05	833.9

54

A F P - CENTRIFUGATION VS. NON-CENTRIFUGATION

IMx ASSAY PRINTOUTS: (Example-AFP CAL, +1 MIN CENTRIFUGATION)

DATE: 0/00/00

TIME: 0:00:00

TECH ID : _____
RGNT LOT : _____
SERIAL # : 000

ASSAY 65 AFP

LOC	CALIB	RATE
1	A	14.3
1	A	16.0
2	B	123.9
2	B	133.6
3	C	327.9
3	C	334.1
4	D	541.1
4	D	523.9
5	E	857.2
5	E	828.1
6	F	1095.3
6	F	1093.3

CALIB	ng/ml	AVGR
A	0	15.2
B	15	128.8
C	50	331.0
D	100	532.5
E	200	842.7
F	350	1094.3

CALIBRATION ACCEPTED

LOC	ID	ng/mL	RATE NOTE
13	_____	18.28	147.7
14	_____	17.85	145.2
15	_____	81.04	456.1
16	_____	80.92	455.6
17	_____	189.18	809.1
18	_____	186.25	800.0

Table 2
**ACCELERATED MICROPARTICLE
SETTLING**

5	Visual Observations after 31 min. centrifugation (1400g) at 2-8°C	
10	Sucrose	AFP
15	13%	pellet*
	14%	pellet
	15%	pellet
	16%	pellet
	17%	even
20	suspension*	
	control	pellet

[* tubes chosen for assay runs]

Table 3
AFP

Original	Microparticles			17% SUCROSE			17% SUCROSE		
	low	med	high	low	med	high	low	med	high
Run 1									
	19.69	86.73	193.11	19.51	79.61	174.71	20.13	78.00	157.85
	19.81	82.45	188.30	18.76	82.78	181.99	18.17	77.52	168.35
	19.32	83.52	184.87	20.88	83.55	180.27	20.53	77.18	160.69
	18.84	87.99	195.96	18.91	88.22	187.89	20.47	80.98	168.38
	20.17	86.68	198.04	19.20	88.20	184.32	19.25	80.08	164.68
	20.35	84.81	194.01	22.47	88.69	192.76	19.84	79.79	169.37
	23.13	89.83	206.79	19.92	87.44	198.00	22.02	81.98	182.90
Mean	20.19	86.00	195.58	19.95	85.50	185.71	20.07	79.36	167.46
SD	1.39	2.58	5.60	1.32	3.52	7.87	1.19	1.83	8.06
%CV	6.90	2.99	2.86	6.62	4.12	4.24	5.92	2.31	4.81
Run 2									
	18.93	76.10	168.66	17.89	73.62	163.49	19.29	74.14	153.28
		71.76	170.85	15.62	73.46	159.01	18.62	73.19	151.02
	17.37	72.28	169.56	17.78	74.32	178.01	18.71	76.92	168.01
	16.58	83.62	176.74	16.94	79.45	171.88	19.52	75.24	168.01
	16.76	79.44	171.02	17.40	76.81	174.14	20.36	77.89	167.17
	19.55	79.44	181.98	18.67	81.48	178.75	20.58	80.41	156.35
	18.80	82.67	198.14	20.13	82.22	190.51	21.15	78.00	177.57
Mean	18.00	77.90	176.71	17.78	77.34	173.68	19.75	76.54	162.77
SD	1.25	4.70	10.56	1.40	3.73	10.41	0.97	2.51	9.54
%CV	6.97	6.04	5.98	7.90	4.83	5.99	4.91	3.28	5.86
Run 3									
	19.40	78.98	179.20	17.62	66.34	152.39	21.20	78.34	166.46
	18.45	75.88	171.08	15.87	72.05	156.68	22.34	78.48	163.09
	17.07	76.75	179.55	13.95	68.92	152.42	20.57	80.57	172.79
	17.14	75.31	179.24	17.16	73.80	163.74	21.24	84.83	177.81
	16.55	80.23	170.24	15.95	70.78	153.87	21.28	85.06	172.10
	18.60	74.55	179.10	16.48	73.69	155.54	22.49	81.51	169.94
	19.14	84.14	188.27	17.94	78.75	180.21	23.10	82.92	189.02
Mean	18.05	77.98	178.10	16.42	72.05	159.26	21.75	81.67	173.03
SD	1.12	3.39	6.06	1.35	3.97	10.02	0.90	2.76	8.48
%CV	6.20	4.34	3.40	8.22	5.51	6.29	4.15	3.37	4.90

Reagent Lot#
32072M100/ 11Jan90

Reagent Lot#
32072M100/ 11Jan90

Reagent Lot#
33172M100/ 20Jun90

SUBSTITUTE SHEET

3/21/2005, EAST Version: 2.0.1.4

Table 4
ORIGINAL MICROPARTICLES - A F P

		<u>LOW</u>	<u>MEDIUM</u>	<u>HIGH</u>
Run 1	MEAN	20.19	86.00	195.58
	STD. DEV.	1.39	2.58	5.60
	%CV	6.90	2.99	2.86
Run 2	MEAN	18.00	77.90	176.71
	STD. DEV.	1.25	4.70	10.56
	%CV	6.97	6.04	5.98
Run 3	MEAN	18.05	77.98	178.10
	STD. DEV.	1.12	3.39	6.06
	%CV	6.20	4.34	3.40
Avg between run CV		6.67	5.77	5.73
		LOT# 32072M100		

MICROPARTICLES - OPTIMUM CONC. AT 17% SUCROSE

		<u>LOW</u>	<u>MEDIUM</u>	<u>HIGH</u>
Run 1	MEAN	19.95	85.50	185.71
	STD. DEV.	1.32	3.52	7.87
	%CV	6.62	4.12	4.24
Run 2	MEAN	17.78	77.34	173.68
	STD. DEV.	1.40	3.73	10.41
	%CV	7.90	4.83	5.99
Run 3	MEAN	16.42	72.05	159.26
	STD. DEV.	1.35	3.97	10.02
	%CV	8.22	5.51	6.29
AVG. BETWEEN RUN CV:		9.86	8.65	7.66
		LOT# 32072M100		

MICROPARTICLES - OPTIMUM CONC. AT 17% SUCROSE

		<u>LOW</u>	<u>MEDIUM</u>	<u>HIGH</u>
Run 1	MEAN	20.07	79.36	167.46
	STD. DEV.	1.19	1.83	8.06
	%CV	5.92	2.31	4.81
Run 2	MEAN	19.75	76.54	162.77
	STD. DEV.	0.97	2.51	9.54
	%CV	4.91	3.28	5.86
Run 3	MEAN	21.75	81.67	173.03
	STD. DEV.	0.90	2.76	8.48
	%CV	4.15	3.37	4.90
AVG. WITHIN RUN CV:		5.23	3.24	3.06
		LOT# 32172M100		

REAL TIME STABILITY/ A F P

Control
Original uParticles
NON-Inverted

Control	Original uParticles	Inverted
---------	---------------------	----------

Optimum Conc.
17% Sucrose

		<u>LOW</u>	<u>MED</u>	<u>HIGH</u>	<u>LOW</u>	<u>MED</u>	<u>HIGH</u>	<u>LOW</u>	<u>MED</u>	<u>HIGH</u>
avg. rate	0	19.84	78.45	183.55	19.15	77.73	171.68	17.97	74.44	170.3
Std. Dev.		0.72	3.06	5.40	0.34	2.35	5.92	0.85	3.76	6.45
%CV		3.61	3.90	2.94	1.80	3.03	3.45	4.73	5.05	3.78
avg rate	5	18.49	78.33	174.89	18.51	73.10	179.81	18.25	78.45	176.5
Std. Dev.		0.82	3.27	9.12	0.46	3.75	4.90	0.21	1.53	7.23
%CV		4.42	4.18	5.22	2.47	5.13	2.73	1.18	1.95	4.10
avg rate	23	18.33	72.92	172.73	18.31	76.76	174.44	18.21	75.35	171.3
Std. Dev.		0.51	2.46	3.63	0.66	2.38	5.52	0.44	2.09	8.68
%CV		2.80	3.37	2.10	3.63	3.10	3.16	2.44	2.78	5.06
avg rate	28	19.33	75.82	180.39	18.75	75.77	174.31	19.55	78.44	173.8
Std. Dev.		0.50	2.13	6.76	0.53	2.79	1.81	0.97	1.75	7.39
%CV		2.61	2.81	3.75	2.84	3.68	1.04	4.95	2.23	4.25
avg rate	46	18.92	75.66	169.90	18.51	76.21	172.70	17.62	78.58	171.5
Std. Dev.		0.74	2.42	9.47	0.27	3.00	8.31	0.58	4.82	8.65
%CV		3.92	3.19	5.58	1.48	3.94	4.81	3.30	6.14	5.04
avg rate	68	17.95	72.35	170.34	17.95	75.98	171.28	17.95	72.35	170.3
Std. Dev.		0.58	1.83	4.05	0.68	3.79	5.22	0.58	1.83	4.05
%CV		3.23	2.53	2.38	3.79	4.98	3.05	3.23	2.53	2.38
tot mean		18.81	75.59	175.30	18.53	75.93	174.04	18.26	76.27	172.3
SD		0.70	2.58	5.56	0.40	1.55	3.11	0.67	2.62	2.42
Bet Run %CV		3.70	3.42	3.17	2.18	2.04	1.79	3.68	3.44	1.40

Table 6

CURRENT AND RECOMMENDED SUCROSE CONCENTRATIONS

5			
	<u>ASSAY</u>	<u>CURRENT SUCROSE CONC.</u>	<u>RECOMMENDED SUCROSE CONC.</u>
10	AFP	13.6%	17%
	TSH	13.6%	17%
	HCG	13.6%	20%
15	T3	13.6%	5%
	TOXO-IgG	0%	17%
20			

Table 7

**SUCROSE OPTIMIZATION
EFFECT ON ASSAY PERFORMANCE**

5 10 15 20 25	ASSAY SUCROSE	OPT SUCROSE CONC. (%)	ORIGINAL UPART			ORIGINAL UPART			OPTIMUM		
			L	M	H	L	M	H	L	M	H
	AFP	17 (13.6)	3.7	3.5	3.2	2.2	2.0	1.8	3.7	3.4	1.4
	TSH	17 (13.6)	11.0	3.9	3.5	8.1	1.3	2.2	3.3	1.4	2.1
	hCG	20 (13.6)	6.8	6.6	6.1	1.8	2.2	2.2	2.1	1.1	0.0
	T3	5 (13.6)	9.3	3.4	2.8	11.2	3.3	3.1	9.9	6.0	2.9
	TOXO IgG*		17 (0)		26.4	8.0		11.0	3.5		
	5.9 4.8										

* Negative and positive control

- No effect on sucrose on cal curve rates except for T3
(ACal decreased 21% with 5% sucrose)

Example 4
Settling and Resuspension of Microparticles

The effects of assay reagent container size, assay
5 reagent container fill volume, and the effect of various
automated agitation steps of an assay reagent container; were
demonstrated employing the following materials and methods:

(i) Microparticles were allowed to settle at room
10 temperature or from between 4°C - 8°C for varying lengths of
time

(ii) Mechanical resuspension of TSH microparticle
reagent in assay reagent containers for resuspension of
microparticles was performed on a reagent pack carousel as
described herein.

15 (iii) Microparticle reagent comprising TSH antibodies
immobilized to microparticles.

(iv) Sucrose reagent concentrations of 0%, 6.5%, and
13%.

(v) Fill volumes of assay reagent containers were 5
20 mL and 10 mL.

(vi) Mixing time course was determined by combining
duplicate 100 uL aliquots of assay reagent solution with 900
uL distilled water at various times, sonicated for 5 minutes,
and spectrophotometrically read at an O.D. of 700.

25

(a) Effect Of Assay Reagent Container Size And
Configuration On Microparticle Settling

(i) The effect of the size of a reagent container was
30 demonstrated employing a large Abbott IMx reagent bottle (30
mL capacity) and a small Abbott IMx reagent bottle (12 mL
capacity) . A TSH microparticle reagent comprising 0%
sucrose was allowed to settle for a 6 week period at fill
volumes of 10 mL in each of the reagent bottles, followed by
35 mechanical resuspension (50/50 step protocol) with
asymmetric pauses. As shown in Figure 34, resuspension was

completed at 2.5 minutes with the large reagent bottle and at 20 minutes with the small reagent bottle.

(ii) The effect of reagent container configuration on the resuspension of a microparticle reagent was demonstrated employing a rectangular 20 mL reagent bottle and a large Abbott IMx reagent bottle. A TSH microparticle reagent comprising 13 % sucrose in each of such reagent bottles was allowed to settle for 5 days at a fill volume of 10 mL each, followed by mechanical resuspension with a 50/50 step protocol with assymetric pauses. As shown in Figure 35 (where V = rectangular bottle and R = Abbott IMx bottle), resuspension of the microparticle reagent was complete within 2.5 minutes in the rectangular bottle and within 1 minute in the Abbott IMx bottle.

(b) Effect Of Assay Reagent Container Fill Volume On Microparticle Settling

The effect of fill volume of an assay reagent container was demonstrated employing small Abbott IMx reagent bottles. A TSH microparticle reagent comprising 13% sucrose in two small bottles (fill volumes of 5 mL and 10 mL) was allowed to settle for a 5 week period at 4°C, followed by mechanical resuspension (50/50 step protocol) with assymetric pauses. As shown in Figure 36, resuspension was completed at 5 minutes in the bottle containing 5 mL and at 20 minutes in the bottle containing 10 mL.

(c) Resuspension Of Microparticles With Automated Agitation

TSH microparticle reagents comprising sucrose concentrations of 0%, 6.5% and 13% were allowed to settle for 6 days at 4°C in large Abbott IMx reagent bottles at fill volumes of 10 mL, followed by mechanical resuspension (50/50 step protocol) with asymmetric pauses. As shown in Figure 37, resuspension of each concentration was completed

completed at 2.5 minutes with the large reagent bottle and at 20 minutes with the small reagent bottle.

(ii) The effect of reagent container configuration on the resuspension of a microparticle reagent was demonstrated employing a rectangular 20 mL reagent bottle and a large Abbott IMx reagent bottle. A TSH microparticle reagent comprising 13 % sucrose in each of such reagent bottles was allowed to settle for 5 days at a fill volume of 10 mL each, followed by mechanical resuspension with a 50/50 step protocol with assymetric pauses. As shown in Figure 35 (where V = rectangular bottle and R = Abbott IMx bottle), resuspension of the microparticle reagent was complete within 2.5 minutes in the rectangular bottle and within 1 minute in the Abbott IMx bottle.

(b) Effect Of Assay Reagent Container Fill Volume On Microparticle Settling

The effect of fill volume of an assay reagent container was demonstrated employing small Abbott IMx reagent bottles. A TSH microparticle reagent comprising 13% sucrose in two small bottles (fill volumes of 5 mL and 10 mL) was allowed to settle for a 5 week period at 4°C, followed by mechanical resuspension (50/50 step protocol) with assymetric pauses. As shown in Figure 36, resuspension was completed at 5 minutes in the bottle containing 5 mL and at 20 minutes in the bottle containing 10 mL.

(c) Resuspension Of Microparticles With Automated Agitation

TSH microparticle reagents comprising sucrose concentrations of 0%, 6.5% and 13% were allowed to settle for 6 days at 4°C in large Abbott IMx reagent bottles at fill volumes of 10 mL, followed by mechanical resuspension (50/50 step protocol) with asymmetric pauses. As shown in Figure 37, resuspension of each concentration was completed

within 2.5 minutes. Although the optimal sucrose concentration for such microparticle reagent is approximately 16%, rapid resuspension of such other TSH microparticles reagents was demonstrated.

5

(d) Effect Of Range Of Automated Agitation Of Microparticles

10 The effect of the range of carousel movement to provide rapid mixing of a TSH microparticle reagent comprising 13% sucrose was demonstrated by allowing such microparticle reagent to settle for 5 days at 4°C in large Abbott IMx reagent bottles, followed by mechanical resuspension thereof utilizing alternating carousel movements of 25, 50, 100 and 150 steps with assymetric pauses. As shown in Figure 38, resuspension of the TSH microparticle reagent was complete within 2 minutes.

20 Example 5
Description Of Kitting And Process Area Activities For Performing An FPIA FPIA On The Continuous And Random Access Analytical System

25 SYSTEM DESCRIPTION OF KITTING AREA FOR PHENOBARBITAL ASSAY

A. ASSUMPTIONS

- 30 1. Analyzer is in Standby/Ready mode when sample is loaded. System has been previously initialized (All motors are homed, syringe and pumps are purged, all electronics and sensors are checked.)
2. Waste has been emptied, Diluent, MEIA buffer, MUP, and Quat bulk liquid consumables have been checked for sufficient volume.
- 35 3. All Consumable inventory files have been updated.

B. PREPARATION STEPS

1. User loads empty Reaction Vessel (RV) into RV carousel.
- 5 2. To load a reagent pack(s), the user must first pause the front end carousels. The system will complete kitting of the current test and transfer the test to the process area.
- 10 3. User opens the reagent carousel cover, loads reagent pack(s) into reagent carousel, closes the reagent carousel cover, then resumes the front-end.
4. Instrument automatically scans all reagent packs onboard to verify reagent status.
 - 15 (a) Each reagent pack is positioned in front of the reagent pack barcode reader by rotation of the reagent carousel.
 - (b) Reagent pack barcode reader reads barcode to identify assay type and carousel location.
 - 20 (c) If the barcode is unreadable, the system will request a barcode override.
 - (d) If the barcode is good or override complete, the system will check the system inventory. The user will be notified if the pack is found to be empty, invalid or outdated. Once the
25 reagent pack is found to be good, it is ready to use.
 - (e) If any resident reagent packs have been added since previous read, carousel will
30 agitate.

C. REQUESTING A TEST

1. User has two options for requesting a test or group of tests on one or more patient samples.
 - 35 (a) User may download the test request loadlist from a host computer to create an order list.
 - (b) User enters test request or creates an order

list on the System directly.

2. If sample cups (no barcode) are used, the following scenario occurs:

- 5 (a) User refers to order list for segment ID and position number to place sample.
- (b) User loads a sample cup into referenced position in segment.
- (c) User transfers patient sample from blood collection tube into sample cup.
- 10 (d) Segment is placed into sample carousel.
- (e) Indication is made to instrument that samples have been loaded.
- (f) Instrument checks consumable inventories, waste status, cal status, etc.
- 15 (g) Sample carousel rotates segment to segment identification reader.
- (h) Instrument reads segment identification.

3. If primary tubes (with barcode) are used, the following scenario occurs (two types of carriers are used for primary tubes: one for tubes with heights of 75 mm and a second for tubes with heights of 100 mm.):

- 20 (a) User loads primary tube into next available segment location on sample carousel.
- (b) Indication is made to instrument that samples are available to be run.
- 25 (c) Instrument checks consumable inventories, waste status, cal status, etc.

30 D. SCHEDULING A TEST

1. When the sample is presented to the pipettor, the System attempts to schedule the tests ordered on that sample for processing. Each test ordered for the sample will be scheduled separately.

- 35 (b) The System checks for adequate inventory (reagent packs, cartridges, buffer, MUP),

system resources, sample time to complete the test.

- 5 (c) The System checks for valid calibration or orders for them on the order list.
- (d) If all test requirements are met, the test is scheduled for processing.
- 10 (e) If all test requirements are not met, the test request is moved to the exception list. Once the test requirements have been met, the test request is moved back to the order list by the user.

2. When a test has been scheduled, the System moves it to the processing list and attempts to schedule other tests ordered for that sample.

- 15 3. When all tests for the current sample have been kitted, the System advances to the next sample on the sample carousel.

20 E KITTING A TEST

1. Once a test is scheduled, it is immediately kitted. (No tests are kitted until the scheduler ensures that the test can be transferred onto the process carousel immediately and processed within the timing requirements of the assay.)
- 25 2. RV carousel is rotated clockwise until an RV is detected in pipette axis position.
3. Reagent pack carousel is rotated until reagent pack for test ordered is at the actuator position. The actuator opens the reagent cartridge caps and the reagent pack carousel is then rotated until a reagent pack for test ordered is in the pipette axis position. After all pipetting steps have been completed, the reagent pack carousel is rotated back to the actuator position where the reagent cartridge caps are closed.
- 30 4. Sample carousel is rotated until sample cup (or primary tube) is in pipette axis position.
- 35

5. Pipette is always at "HOME" position (Pipette R-axis is parked over wash station and Pipette Z-axis is at the Z-clear position) when not in use.

6. Sample kitting.

5

(a) Sample aspirate.

10

15

20

25

30

35

(i) Syringe aspirates "X" uL of air at a rate of "X" ul/sec.

(ii) Pipette R-axis is moved over sample cup.

(iii) Pipette Z-axis is moved down to the Z-above position.

(iv) LLS is enabled to ensure that no liquid is currently detected.

(v) Pipette Z-axis is moved down at constant speed until fluid is detected or until Z-Asp limit has been reached (It will be assumed that fluid is detected)

(vi) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present in the well, the aspiration sequence is initiated (If insufficient volume is present, the test is aborted and the test request moved to the exception list. The exception list provides notice to an operator of tests which cannot be completed).

(vii) The following occur simultaneously until the total volume of sample required is aspirated:

(1) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.

- 5
- (2) Syringe motor aspirates "X" uL at a rate of "X" ul/sec.
- (3) LLS is checked to ensure probe still in liquid. Liquid Level Sense (LLS) is disabled. Pipette Z-axis is moved up to Z-clear position.
- (4) Pipette R-axis is moved over the RV sample well.
- 10
- (5) Pipette Z-axis is moved down to the dispense position within the RV sample well.
- (6) Syringe dispenses "X" uL of sample at a rate of "X" ul/sec.
- 15
- (7) Pipette Z-axis is moved up to Z-clear position.

(b) Probe Post-Wash

20

The probe is washed to ensure that it is free from contamination. It is to be understood that all pipette activities (in both kitting and process areas) are followed with a probe post-wash to minimize carryover from one fluid aspirate to another. In some cases, pipette activities may be preceded with a probe prewash if necessary to

25

guarantee the validity of the next fluid aspirate. For this assay description, it will be assumed that only a post-wash is used.

- (i) The inside of the probe is cleaned first.
- 30
- (1) Pipette R-axis is moved over waste area.
- (2) Pipette Z-axis is moved down to appropriate position within the waste area.
- 35
- (3) The wash valve is opened for the amount of time specified in the assay protocol.
- (4) Wash valve is closed.

- (5) Pipette Z-axis is moved up to the Z-clear position.
- (ii) The outside of the probe is cleaned next.
- 5 (1) Pipette R-axis is moved over wash cup.
- (2) Pipette Z-axis is moved down to wash position within the wash cup.
- 10 (3) The wash valve is opened for the amount of time specified in the assay protocol.
- (4) Wash valve is closed.
- (iii) Pipette is returned to "HOME" position.
- 15 7. Popper kitting ("Popper" is defined as a substance which eliminates in general interfering substances in assays such as, for example, those discussed and claimed in U.S. Patent 4,492,762 issued January 8, 1985 and hereby incorporated by reference).
- 20 (a) Popper aspirate.
- (i) Syringe aspirates "X" uL of air at a rate of "X" ul/sec.
- (ii) Pipette R-Axis is moved over the popper reagent bottle in the Reagent Pack.
- 25 (iii) Pipette Z-axis is moved down to the Z-above position.
- (iv) LLS is enabled to ensure no liquid currently detected.
- 30 (v) Pipette Z-axis is moved down at constant speed until fluid is detected or until the Z-aspiration-lower (Z-Asp) limit is reached (it will be assumed that fluid is detected).
- 35 (vi) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System

- 5 calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present in the well, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 10 (vii) The following occur simultaneously until the total volume of popper required is aspirated:
- 15 (1) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
- (2) Syringe aspirates "X" uL at a rate of "X" ul/sec.
- (3) LLS is checked to ensure probe still in liquid.
- (4) LLS is disabled.
- 20 (5) Pipette Z-axis is moved up to Z-clear position.
- (6) Pipette R-axis is moved over the RV reagent 1 well.
- (7) Pipette Z-axis is moved down to the dispense position within the RV reagent 1 well.
- 25 (8) Syringe dispenses "X" uL of popper at a rate of "X" ul/sec.
- (9) Pipette Z-axis is moved up to Z-clear position.
- 30 (b) Probe post-wash.
The probe is again washed to ensure that it is free from contamination as described in section 6 (Sample Kitting).
- 35 8. Antiserum kitting
(a) Antiserum aspirate

- 5
- (i) Syringe aspirates "X" uL of air at a rate of "X" ul/sec.
- (ii) Pipette R-Axis is moved over the antiserum reagent bottle in the Reagent Pack.
- 10
- (iii) Pipette Z-axis is moved down to the Z-above position.
- (iv) LLS is enabled to ensure no liquid currently detected.
- 15
- (v) Pipette Z-axis is moved down at constant speed until fluid is detected or until the Z-Asp limit is reached (it will be assumed that fluid is detected).
- 20
- (vi) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present in the well, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 25
- (vii) The following occur simultaneously until the total volume of antiserum required is aspirated:
- 30
- (1) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
- (2) Syringe aspirates "X" micro liter (uL) at a rate of "X" ul/sec. LLS is checked to ensure probe still in liquid.
- 35
- (3) LLS is disabled.
- (4) Pipette Z-axis is moved up to Z-clear position.

- 5
- (5) Pipette R-axis is moved over the RV reagent 2 well.
 - (6) Pipette Z-axis is moved down to the dispense position within the RV reagent 2 well.
 - (7) Syringe dispenses "X" uL of antiserum at a rate of "X" ul/sec.
 - (8) Pipette Z-axis is moved up to Z-clear position.
- 10
- (b) Probe post-wash.
The probe is again washed to ensure that it is free from contamination as described in section 6 (Sample Kitting).
- 15
9. Tracer kitting.
- (a) Tracer aspirate.
 - (i) Syringe aspirates "X" uL of air at a rate of "X" ul/sec.
 - (ii) Pipette R-Axis is moved over the tracer reagent bottle in the Reagent Pack.
 - (iii) Pipette Z-axis is moved down to the Z-above position.
 - (iv) LLS is enabled to ensure no liquid currently detected.
 - (v) Pipette Z-axis is moved down at constant speed until fluid is detected or until the Z-Asp limit is reached (it will be assumed that fluid is detected).
 - (vi) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present in the well, the aspiration sequence is initiated.(if sufficient volume not is
- 20
- 25
- 30
- 35

present, the test is aborted and the test request moved to the exception list).

- 5 (vii) The following occur simultaneously until the total volume of tracer required is aspirated:
- 10 (1) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
 - (2) Syringe aspirates "X" uL at a rate of "X" ul/sec.
 - (3) LLS is checked to ensure probe still in liquid.
 - (4) LLS is disabled.
 - 15 (5) Pipette Z-axis is moved up to Z-clear position.
 - (6) Pipette R-axis is moved over the RV reagent 3 well.
 - (7) Pipette Z-axis is moved down to the dispense position within the RV reagent 2 well.
 - 20 (8) Syringe dispenses "X" uL of tracer at a rate of "X" ul/sec.
 - (9) Pipette Z-axis is moved up to Z-clear position.
 - 25 (b) Probe post-wash.
- The probe is again washed to ensure that it is free from contamination as described in section 6 (Sample Kitting).

30 F. TRANSFER OF REACTION VESSEL (RV) INTO PROCESS AREA

1. RV carousel is rotated to transfer station.
2. Process carousel is rotated so that the empty position is aligned with the transfer station.
- 35 3. Transfer mechanism 0-axis is rotated to sample entry area.

4. Transfer mechanism R-axis grabs the RV and pulls it into the transfer mechanism.
5. Transfer mechanism O-axis is rotated so that RV is aligned with the empty position on the process carousel.
- 5 6. RV is loaded onto process carousel.

SYSTEM DESCRIPTION OF FPIA PROCESS AREA FOR
PHENOBARBITAL

- 10 A. Wait for temperature equilibration time and evaporation window to expire.
- B. FIRST PIPETTE ACTIVITY (preparation of sample blank comprising diluted sample and popper).
 - 15 1. Incubation timer is set according to assay file specifications.
 2. Precision diluent aspirate. The following activities are performed simultaneously:
 - 20 (a) Syringe aspirates "X" uL at a rate of "X" ul/sec.
 - (b) Wash valve is opened.
 - (c) Wait "n" seconds.
 - (d) Wash valve is closed.
 3. Sample aspirate.
 - 25 (a) Pipette R-axis is moved over the RV sample well.
 - (b) LLS is enabled to ensure no liquid currently detected.
 - (c) Pipette Z-axis is moved down at constant speed until fluid is detected OR until the Z-Asp limit is reached (it will be assumed that fluid is detected).
 - 30 (d) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description.
 - 35

- 5 If sufficient volume is present, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 10 (e) The following occur simultaneously until the total volume of sample required is aspirated:
- (i) Pipettor Z-axis motor is moved down at a rate of "X" steps/sec.
 - (ii) Syringe aspirates "x" uL of sample at a rate of "X" ul/sec.
 - (iii) LLS is checked to ensure probe still in liquid.
 - (iv) LLS is disabled.
 - 15 (v) Pipette Z-axis is moved up to Z-above position.
4. Diluent/sample dispensed to the RV predilute well.
- (a) Pipette R-axis is moved over the RV predilute well.
 - 20 (b) Pipette Z-axis is moved down to the dispense position within the RV predilute well.
 - (c) Syringe dispenses "X" uL of diluent/sample at a rate of "X" ul/sec.
 - (d) Pipette Z-axis is moved up to Z-clear position.
 - 25
5. Probe post-wash.
- The probe is again washed to ensure that it is free from contamination as described in section 6 (Sample kitting).
- 30 6. Precision diluent aspirate. The following activities are performed simultaneously:
- (a) Syringe aspirates "X" uL at a rate of "X" ul/sec.
 - (b) Wash valve is opened.
 - 35 (c) Wait "n" seconds.
 - (d) Wash valve is closed.
7. Popper aspirate.

- 5 (a) Pipette R-axis is moved over the RV Reagent (popper) well.
- (b) LLS is enabled to ensure no liquid currently detected.
- 10 (c) Pipette Z-axis is moved down at constant speed until fluid is detected or until the Z-Asp limit is reached (it will be assumed that fluid is detected).
- (d) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 15 (e) The following occur simultaneously until the total volume of popper required is aspirated:
- 20 (i) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
- (ii) Syringe aspirates "X" uL at a rate of "x" ul/sec.
- 25 (iii) LLS is checked to ensure probe still in liquid.
- (iv) LLS is disabled.
- (v) Pipette Z-axis is moved up to the Z-above position.
- 30 8. Diluted sample aspirate.
- (a) Pipette R-axis is moved over the RV predilute well.
- (b) LLS is enabled to ensure no liquid currently detected.
- 35 (c) Pipette Z-axis is moved down at constant speed until fluid is detected or until the Z-Asp limit is reached (it will be assumed that fluid is detected).

- 5 (d) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 10 (e) The following occur simultaneously until the total volume of diluted sample required is aspirated:
- 15 (i) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
- (ii) Syringe aspirates "X" uL at a rate of "x" uL/sec.
- (iii) LLS is checked to ensure probe still in liquid.
- 20 (iv) LLS is disabled.
- (v) Pipette Z-axis is moved up to the Z-above position.
11. Diluted sample/popper diluent dispensed to RV cuvette.
- 25 (a) Pipette R-axis is moved over to the RV cuvette position.
- (b) Pipette Z-axis is moved down to the dispense position in the RV cuvette.
- 30 (c) Syringe dispenses "X" uL of diluted sample/popper/diluent at a rate of "X" uL/sec.
- (d) Pipette Z-axis is moved up to the Z-above position.
12. Probe post-wash.
- 35 The probe is again washed to ensure that it is free from contamination as described in section 6 (sample kitting) to complete first pipette activity

C. BLANK READ PREPARATION

When incubation timer expired, the following activities are started:

- 5 1. The FPIA reader is prepared to take a read; lamp intensity is brought from simmer state to burn state.
2. Photomultiplier tube (PMT) gain is set.

D. BLANK READ (BACKGROUND)

- 10 1. Incubation timer is set according to assay file specifications.
2. Process carousel is rotated so that the RV is at the read station.
- 15 3. Horizontal intensity is read for "X.XX" seconds.
4. The crystal is flipped for the vertical read.
5. Wait "n" seconds until the crystal settles.
6. Vertical intensity is read for "X.XX" seconds.
7. The raw reads are converted to normalized reads
- 20 (light intensity hitting detector/lamp intensity) by the optics microprocessor.
8. Background reads are stored.
9. System calculates BLANK I to complete blank read.
10. Next activity started when incubation timer
- 25 expires.

E SECOND PIPETTE ACTIVITY (for reaction between diluted sample, popper, tracer and antiserum).

- 30 1. Incubation timer is set according to assay file specifications.
2. Precision diluent aspirate.
 - (a) The following activities are performed simultaneously:
 - 35 (i) Syringe aspirates "X" uL at a rate of "X" ul/sec.
 - (ii) Wash valve is opened.

- (iii) Wait "n" seconds.
- (iv) Wash valve is closed.
- 3. Antiserum aspirate.
 - (i) Pipette R-axis is moved over the RV Reagent 2 (antiserum) well.
 - (ii) LS is enabled to ensure no liquid currently detected.
 - (iii) Pipette Z-axis is moved down at constant speed until fluid is detected OR until the Z-Asp limit is reached (it will be assumed that fluid is detected).
 - (iv) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present, the aspiration sequence is initiated. (If sufficient volume is not present, the test is aborted and the test request moved to the exception list.)
 - (v) The following occur simultaneously until the total volume of antiserum required is aspirated:
 - (1) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
 - (2) Syringe aspirates "X" uL at a rate of "X" ul/sec.
 - (3) LLS is checked to ensure probe still in liquid.
 - (4) LLS is disabled.
 - (5) Pipette Z-axis is moved up to the Z-above position.
- 4. Tracer aspirate.
 - (a) Syringe aspirates "X" uL of air at a rate of "X" ul/sec.
 - (b) Pipette R-axis is moved over the RV Reagent 3 (tracer) well.

- 5 (c) LLS is enabled to ensure no liquid currently detected.
- (d) Pipette Z-axis is moved down at constant speed until fluid is detected OR until the Z-
Asp limit is reached (it will be assumed that fluid is detected).
- 10 (e) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 15 (f) The following occur simultaneously until the total volume of tracer required is aspirated:
- (i) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
- 20 (ii) Syringe aspirates "X" uL at a rate of "X" ul/sec.
- (iii) LLS is checked to ensure probe still in liquid.
- (v) LLS is disabled.
- 25 (vi) Pipette Z-axis is moved up to the Z-above position.
5. Diluted sample aspirate.
- (a) Pipette R-axis is moved over the RV predilute well.
- 30 (b) LLS is enabled to ensure no liquid currently detected.
- (c) Pipette Z-axis is moved down at constant speed until fluid is detected OR until the Z-
Asp limit is reached (it will be assumed that fluid is detected).
- 35 (d) Based on the Z-height position at which fluid is detected and the Z-height/volume table,

the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list.)

(e) The following occur simultaneously until the total volume of diluted sample required is aspirated:

(1) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.

(2) Syringe aspirates "X" uL at a rate of "X" uL/sec.

(3) LLS is checked to ensure probe still in liquid.

(4) LLS is disabled.

(5) Pipette Z-axis is moved up to the Z-above position.

6. Diluted sample/tracer/aspirate/antiserum/diluent dispensed to RV cuvette.

(a) Pipette R-axis is moved over to the RV cuvette position.

(b) Pipette Z-axis is moved down to the dispense position in the RV cuvette.

(c) Syringe dispenses "X" uL of diluted sample/tracer/air/antiserum/diluent at a rate of "X" ul/sec.

(d) Pipette Z-axis is moved up to the Z-above position.

7. Probe post-wash.

The probe is again washed to ensure that it is free from contamination as described in section 6 (Sample kitting) to complete the second pipette activity.

8. Next activity started when incubation timer expires.

E FINAL READ PREPARATION

1. The FPIA reader is prepared to take a read; lamp intensity is brought from simmer state to burn state.
- 5 2. PMT gain is set.

F FINAL READ

- 10 1. Process carousel is rotated so that the RV is at the read station.
2. Horizontal intensity is read for "X.XX" seconds.
3. The crystal is flipped for the vertical read.
4. The System delays "n" seconds until the crystal settles.
- 15 5. Vertical intensity is read for "X.XX" seconds.
6. The raw reads are converted to normalized reads (light intensity hitting detector/lamp intensity) by the optics microprocessor.
7. Reads are stored.
- 20 8. System calculates NET intensity (I) and millipolarization (mP).
9. mP value is fitted to calibration curve to yield a concentration result.

25 G RV UNLOAD (this activity occurs when resources are not in use. The following are performed simultaneously:

1. Process carousel is rotated so that the empty position is at the transfer station. Transfer mechanism 0-axis is moved to process carousel.
- 30 2. RV is grabbed with the transfer mechanism R-axis and pulled into the transfer mechanism.
3. Transfer mechanism 0-axis is rotated so that RV is aligned with the waste container.
4. RV is pushed into the waste container.
- 35

Example 6
Description Of Kitting And Process Area Activities For
Performing An MEIA On The Continuous And Random Access
Analytical System

5

SYSTEM DESCRIPTION OF KITTING AREA FOR CEA ASSAY

A. ASSUMPTIONS

- 10 1. Analyzer is in Standby/Ready mode when sample is loaded. System has been previously initialized (All motors are homed, syringe and pumps are purged, all electronics and sensors are checked).
- 15 2. Waste has been emptied, dilution, MEIA buffer, MUP, and Quat bulk liquid consumables have been checked for sufficient volume.
- 20 3. Cartridges have been placed into hopper and are available for loading onto auxiliary carousel when needed (for MEIA assays only).
- 20 4. All Consumable inventory files have been updated.

B. PREPARATION STEPS

- 25 1. User loads empty RVs into RV carousel.
- 25 2. To load a reagent pack(s), the user must first pause the front end carousels. The system will complete kitting of the current test and transfer the test to the process area.
- 30 3. User opens the reagent carousel, loads reagent pack(s) into reagent carousel, closes the reagent carousel cover, then resumes the front-end.
- 30 4. Instrument automatically scans all reagent packs onboard to verify reagent status.
- 35 5. Each reagent pack is positioned in front of the reagent pack barcode reader by rotation of the reagent carousel.

6. Reagent pack barcode reader reads barcode to identify assay type and carousel location. If the barcode is unreadable, the system will request a barcode override.

5 7. If the barcode is good or override complete, the system will check the system inventory. The user will be notified if the pack is found to be empty, invalid or outdated. Once the reagent pack is found to be good, it is ready to use.

8. If any resident reagent packs have been added since previous read, carousel will agitate.

10

C. REQUESTING A TEST

1. User has two options for requesting a test or group of tests on one or more patient samples.

- 15 (a) User may download the test request loadlist from a host computer to create an order list.
- (b) User enters test request or creates an order list on the System directly.

20 2. If sample cups (no barcode) are used, the following scenario occurs:

- (a) User refers to order list for segment ID and position number to place sample.
- (b) User loads a sample cup into referenced position in segment.
- 25 (c) User transfers patient sample from blood collection tube into sample cup.
- (d) Segment is placed into sample carousel.
- (e) Indication is made to instrument that samples have been loaded.
- 30 (f) Instrument checks consumable inventories, waste status, assay calibration, etc.
- (g) Sample carousel rotates segment to segment identification reader.
- (h) Instrument reads segment identification.

35 3. If primary tubes (with barcode) are used, the following scenario occurs:

- 5 (a) User loads primary tube into next available segment location on sample carousel (two types of carriers are used for primary tubes: one for tubes with heights of 75 mm and a second for tubes with heights of 100 mm.).
- (b) Indication is made to instrument that samples are available to be run.
- 10 (c) Sample carousel rotates segment to segment identification reader.

D. SCHEDULING A TEST

- 15 1. When the sample is presented to the pipettor, the System attempts to schedule the tests ordered on that sample for processing. Each test ordered for the sample will be scheduled separately.

- 20 (a) The System checks for adequate inventory (reagent packs, cartridges, buffer, MUP), system resources, sample time to complete the test.
- (b) The System checks for valid calibration or orders for them on the order list.
- (c) If all test requirements are met, the test is scheduled for processing.
- 25 (d) If all test requirements are not met, the test request is moved to the exception list. Once the test requirements have been met, the test request is moved back to the order list by the user.

- 30 2. When a test has been scheduled, the system moves it to the processing list and attempts to schedule other tests ordered for that sample.

- 35 3. When all tests for the current sample have been kitted, the System advances to the next sample on the sample carousel.

E. KITTING A TEST

1. Once a test is scheduled, it is immediately kitted.
(no tests are kitted until the scheduler ensures that the test can be transferred onto the process carousel immediately and processed within the timing requirements of the assay).
- 5 2. RV carousel is rotated clockwise until an RV is detected in pipette axis position.
3. Reagent pack carousel is rotated until reagent pack for test ordered is at the actuator position. The actuator opens the reagent cartridge caps and the reagent pack carousel
10 is then rotated until reagent pack for test ordered is in the pipette axis position. After all pipetting steps have been completed, the reagent pack carousel is rotated back to the actuator position where the reagent cartridge caps are closed.
4. Sample carousel is rotated until sample cup (or
15 primary tube) is in pipette axis position.
5. Pipette is always at HOME position (Pipette R-axis is parked over wash station and Pipette Z-axis is at the Z-clear position) when not in use.
6. Sample kitting.
20 (a) Sample aspirate.
 - (i) Syringe aspirates "X" uL of air at a rate of "X" ul/sec.
 - (ii) Pipette R-axis is moved over sample cup.
 - 25 (iii) Pipette Z-axis is moved down to the Z-above position.
 - (iv) Pipette Z-axis is moved down to the Z-LLS position.
 - (v) LLS is enabled to ensure that no liquid
30 is currently detected.
 - (vi) Pipette Z-axis is moved down at constant speed until fluid is detected or until Z-Asp limit has been reached (it will be assumed that fluid is
35 detected).
 - (vii) Based on the Z-height position at which fluid is detected and the Z-

- 5 height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present in the well, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 10 (viii) The following occur simultaneously until the total volume of sample required is aspirated:
- 15 (1) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
(2) Syringe aspirates "X" uL at a rate of "X" ul/sec.
(3) LLS is checked to ensure probe still in liquid.
- 20 (4) LLS is disabled.
(5) Pipette Z-axis is moved up to Z-clear position.
(6) Pipette R-axis is moved over the RV sample well.
- 25 (7) Pipette Z-axis is moved down to the dispense position within the RV sample well.
- 30 (8) Syringe dispenses "X" uL of sample at a rate of "X" ul/sec.
(9) Pipette Z-axis is moved up to Z-clear position.
- 35 (b) Probe post-wash.
The probe is washed to ensure that it is free from contamination. It is to be understood that pipette activities in both kitting and

process areas are generally followed with a probe post-wash to minimize carryover from one fluid aspirate to another. In some cases, pipette activities may be preceded with a probe prewash if necessary to guarantee the validity of the next fluid aspirate. For this assay description, it will be assumed that only a post-wash is used .

5

(i) The inside of the probe is cleaned first.

10

(1) Pipette R-axis is moved over waste area.

(2) Pipette Z-axis is moved down to appropriate position within the waste area.

15

(3) The wash valve is opened for the amount of time specified in the assay protocol.

(4) Wash valve is closed.

20

(ii) Pipette Z-axis is moved up to the Z-clear position.

(iii) The outside of the probe is cleaned next.

(1) Pipette R-axis is moved over wash cup.

25

(2) Pipette Z-axis is moved down to wash position within the wash cup.

(3) The wash valve is opened for the amount of time specified in the assay protocol.

30

(4) Wash valve is closed.

(5) Pipette is returned to "HOME" position.

7. Microparticle kitting.

35

(a) Microparticle aspirate (microparticles are pipetted directly into the RV incubation well

to save on volume, as this is the most costly MEIA reagent).

- (i) Syringe aspirates "X" uL of air at a rate of "X" ul/sec.
- 5 (ii) Pipette R-Axis is moved over the microparticle reagent bottle in the Reagent Pack.
- (iii) Pipette Z-axis is moved down to the Z-above position.
- 10 (iv) Pipette Z-axis is moved down to the Z-LLS position.
- (v) LLS is enabled to ensure no-liquid currently detected.
- 15 (vi) Pipette Z-axis is moved down at constant speed until fluid is detected or until the Z-Asp limit is reached (it will be assumed that fluid is detected)
- (vii) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present in the well, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 20
- 25
- 30 (viii) The following occur simultaneously until the total volume of microparticles required is aspirated:
 - (1) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
 - 35 (2) Syringe aspirates "X" uL at a rate of "X" ul/sec.

- (3) LLS is checked to ensure probe still in liquid.
- (ix) LLS is disabled.
- (x) Pipette Z-axis is moved up to Z-clear position.
- (xi) Pipette R-axis is moved over the RV incubation well.
- (xii) Pipette Z-axis is moved down to the dispense position within the RV incubation well.
- (xiii) Syringe dispenses "X" uL of microparticles at a rate of "X" ul/sec. Pipette Z-axis is moved up to Z-clear position.
- (b) Probe post-wash.
- The probe is again washed to ensure that it is free from contamination as described in section 6 (Sample kitting).
8. Conjugate kitting.
- (a) Conjugate aspirate (conjugate, special wash fluid, and/or specimen diluent are pipetted into either RV reagent wells or RV predilution well, depending on volume requirements).
- (i) Syringe aspirates "X": uL of air at a rate of "X" ul/sec.
- (ii) Pipette R-Axis is moved over the conjugate reagent bottle in the Reagent Pack.
- (iii) Pipette Z-axis is moved down to the Z-above position.
- (iv) Pipette Z-axis is moved down to the Z-LLS position.
- (v) LLS is enabled to ensure no liquid currently detected.
- (vi) Pipette Z-axis is moved down at constant speed until fluid is detected

- or until the Z-Asp limit is reached (it will be assumed that fluid is detected).
- 5 (vii) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present in the well, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 10
- 15 (viii) The following occur simultaneously until the total volume of conjugate required is aspirated:
- 20 (1) Pipette Z-axis motor is moved down at a rate of "x" steps/sec.
- (2) Syringe aspirates "X" uL at a rate of "X" ul/sec.
- (3) LLS is checked to ensure probe still in liquid.
- 25 (ix) LLS is disabled.
- (x) Pipette Z-axis is moved up to Z-clear position.
- (xi) Pipette R-axis is moved over the RV reagent well.
- 30 (xii) Pipette Z-axis is moved down to the dispense position within the RV r reagent well.
- (xiii) Syringe dispenses "X" uL of conjugate at a rate of "X" ul/sec.
- 35 (xiv) Pipette Z-axis is moved up to Z-clear position.
- (b) Probe post-wash.

The probe is again washed to ensure that it is free from contamination as described in section 6 (Sample kitting).

9. MEIA Buffer Kitting.

- 5 (a) RV Carousel is rotated until RV buffer well is under the MEIA buffer dispenser at buffer kitting station.
- (b) "X" uL of MEIA buffer is dispensed into the buffer well at a rate of "X" ul/sec

10

F. TRANSFERRING RV INTO PROCESS AREA

1. RV carousel is rotated to transfer station.
2. Process carousel is rotated so that the empty
- 15 position is aligned with the transfer station.
3. Transfer mechanism O-axis is rotated to sample entry area.
4. Transfer mechanism R-axis grabs the RV and pulls it into the transfer mechanism.
- 20 5. Transfer mechanism O-axis is rotated so that RV is aligned with the empty position on the process carousel.
6. RV is loaded onto process carousel.

SYSTEM DESCRIPTION OF MEIA PROCESS AREA FOR CEA

25

A. System waits for temperature equilibration time and evaporation window to expire.

B. FIRST PIPETTE ACTIVITY (microparticle/sample

30 reaction)

1. Incubation timer is set according to assay file specifications.
2. MEIA buffer aspirate.
- 35 (a) The process carousel is moved so that the RV is at the pipetting station.

- 5 (b) Syringe aspirates "X" uL of air at a rate of "X" ul/sec.
- (c) Pipette R-axis is moved over the RV buffer well.
- (d) Pipette Z-axis is moved down to the Z-above position over the RV buffer well.
- (e) Pipette Z-axis is moved down to the Z-LLS position.
- 10 (f) LLS is enabled to ensure no liquid currently detected.
- (g) Pipette Z-axis is moved down at constant speed until fluid is detected or until the Z-Asp limit is reached (it will be assumed that fluid is detected).
- 15 (h) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If
- 20 sufficient volume is present, the aspiration sequence is initiated (If sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 25 (i) The following occur simultaneously until the total volume of MEIA buffer required is aspirated:
- (1) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
- (2) Syringe aspirates "X" uL at a rate of "X" ul/sec.
- 30 (j) LLS is checked to ensure probe still in liquid.
- (k) LLS is disabled.
- (l) Pipette Z-axis is moved up to Z-above position.
- 35 3. Sample aspirate
- (a) Pipette R-axis is moved over the RV sample well.

- 5 (b) Pipette Z-axis is moved down to the Z-LLS position.
- (c) LLS is enabled to ensure no liquid currently detected.
- 10 (d) Pipette Z-axis is moved down at constant speed until fluid is detected or until the Z-Asp limit is reached (it will be assumed that fluid is detected).
- 15 (e) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the system calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 20 (f) The following occur simultaneously until the total volume of sample required is aspirated:
- (1) Pipettor Z-axis motor is moved down at a rate of "X" steps/sec.
- (2) Syringe aspirates "X" uL at a rate of "X" ul/sec.
- 25 (g) LLS is checked to ensure probe still in liquid.
- (h) LLS is disabled.
- (i) Pipette Z-axis is moved up to Z-above position.
- 30 4. MEIA buffer and sample are added to microparticles in incubation well.
- (a) Pipette Z-axis is moved down to the dispense position within the RV incubation well.
- (b) Syringe dispenses "X" uL of MEIA buffer and sample at a rate of "X" ul/sec.
- 35 (c) Pipette Z-axis is moved up to Z-clear position.
5. Probe post-wash.

The probe is again washed to ensure that it is free from contamination as described in section 6 (Sample kitting).

5 C. CARTRIDGE LOAD (This activity occurs when resources are not in use)

1. Move the auxiliary carousel so that reserved position is under feeder.
2. Cycle trap-door mechanism to load flashlight into
10 carousel.
3. Cycle shuttle mechanism to place another MEIA cartridge on trap door (for next tab load).
4. Check incubation timer. When expires start next
15 pipetting.

D. SECOND PIPETTE ACTIVITY (transfer of reaction mixture to matrix)

- 20 1. Incubation timer is set according to assay file specifications.
2. Buffer aspirate.
 - (a) The process carousel is moved so that the RV is at the pipetting station.
 - 25 (b) Syringe aspirates "X" uL of air at a rate of "X" ul/sec.
 - (c) Pipette R-axis is moved over the RV buffer well.
 - (d) Pipette Z-axis is moved down to the Z-above
30 position.
 - (e) Pipette Z-axis is moved down to the Z-LLS position.
 - (f) LLS is enabled to ensure no liquid currently detected.
 - 35 (g) Pipette Z-axis is moved down at constant speed until fluid is detected or until the Z-

- Asp limit is reached (it will be assumed that fluid is detected).
- 5 (h) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the system calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 10 (i) The following occur simultaneously until the total volume of buffer required is aspirated:
- 15 (1) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
(2) Syringe aspirates "X" uL at a rate of "X" ul/sec.
- (j) LLS is checked to ensure probe still in liquid.
(k) LLS is disabled.
- 20 (l) Pipette Z-axis is moved up to the Z-above position.
3. Reaction mixture aspirate.
- (a) Pipette R-axis is moved over the RV incubation well.
- 25 (b) Pipette Z-axis is moved down to the Z-LLS position.
- (c) LLS is enabled to ensure no liquid currently detected.
- (d) Pipette Z-axis is moved down at constant speed until fluid is detected or until the Z-Asp limit is reached (it will be assumed that fluid is detected).
- 30 (e) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the system calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If
- 35

- sufficient volume is present, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).
- 5 (f) The following occur simultaneously until the total volume of reaction mixture required is aspirated:
- (1) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
- 10 (2) Syringe aspirates "X" uL at a rate of "X" ul/sec.
- (g) LLS is checked to ensure probe still in liquid.
- (h) LLS is disabled.
- (i) Pipette Z-axis is moved up to the Z-clear position.
- 15 4. Reaction mixture dispense onto matrix.
- (a) The following are performed simultaneously and concurrently with the reaction mixture aspirate (above):
- (i) The auxiliary carousel is moved so that the cartridge is at the pipetting station.
- 20 (ii) Pipette R-axis is moved over the MEIA cartridge (matrix) surface.
- (iii) Pipette Z-axis is moved down to the matrix dispense position.
- 25 (iv) Syringe dispenses "X" uL of reaction mixture at a rate of "X" ul/sec.
- (v) System delays "X" seconds until reaction mixture has been absorbed by matrix.
- 30 5. Buffer wash of matrix.
- (a) Syringe dispenses "X" uL of buffer at a rate of "X" ul/sec.
- (b) Pipette Z-axis is moved up to the Z-clear position.
- 35 6. Probe post-wash.

The probe is again washed to ensure that it is free from contamination as described in section 6 (Sample kitting).

7. When incubation timer expires, next pipette activity begins.

5

E THIRD PIPETTE ACTIVITY (conjugate addition)

1. Incubation timer is set according to assay file specifications.

10

2. Conjugate aspirate.

(a) The process carousel is moved so that the RV is at the pipetting station.

(b) Syringe aspirates "X" uL of air at a rate of "X" ul/sec.

15

(c) Pipette R-axis is moved over the RV reagent 1 (conjugate) well.

(d) Pipette Z-axis is moved down to the Z-above position.

20

(e) LLS is enabled to ensure no liquid currently detected.

(f) Pipette Z-axis is moved down at constant speed until fluid is detected or until the Z-Asp limit is reached (it will be assumed that fluid is detected).

25

(g) Based on the Z-height position at which fluid is detected and the Z-height/volume table, the System calculates the volume of fluid in the well and compares it to the volume specified in the pipetting description. If sufficient volume is present, the aspiration sequence is initiated (if sufficient volume is not present, the test is aborted and the test request moved to the exception list).

30

(h) The following occur simultaneously until the total volume of conjugate required is aspirated:

35

- (i) Pipette Z-axis motor is moved down at a rate of "X" steps/sec.
 - (ii) Syringe aspirates "X" uL at a rate of "X" ul/sec.
- 5 (i) LLS is checked to ensure probe still in liquid.
- (j) LLS is disabled.
- (k) Pipette Z-axis is moved up to the Z-clear position.
- 3. Conjugate dispense (performed simultaneously).
- 10 (a) The auxiliary carousel is moved so that the cartridge is at the pipetting station.
- (b) Pipette R-axis is moved over the cartridge (matrix) surface.
- (c) Pipette Z-axis is moved down to the matrix dispense position.
- 15 (d) Syringe dispenses "X" uL of conjugate at a rate of "X" ul/sec.
- (e) Pipette Z-axis is moved up to the Z-clear position.
- 20 (f) Wait "X" seconds until reaction mixture has been absorbed by matrix.
- 4. Probe post-wash.
 The probe is again washed to ensure that it is free from contamination as described in section 6 (Sample kitting).
- 25 F. RV UNLOAD (This activity occurs when resources are not in use)
- 1. The following are performed simultaneously:
 - 30 (a) Process carousel is rotated so that the empty position is at the transfer station.
 - (b) Transfer mechanism O-axis is moved to process carousel.
- 2. RV is grabbed with the transfer mechanism R-axis and pulled into the transfer mechanism.
- 35 3. Transfer mechanism O-axis is rotated so that RV is aligned with the waste container.

101

4. RV is pushed into the waste container.
5. Check incubation timer. When expires start next activity.

5 G MEIA READ PREPARATION

1. Lamp intensity is brought from simmer state to burn state.
2. PMT gain is set.

10

H MATRIX WASH

1. Auxiliary carousel is rotated so that the cartridge
15 is at the matrix wash station.
2. The following steps are repeated until all the
buffer specified in the assay file for cartridge wash has been
dispensed.
 - (a) "X" uL of heated MEIA buffer are dispensed in
20 50uL cycles at a rate of "X" uL/sec onto the
matrix.
 - (b) Wait "n" seconds.

I. MUP DISPENSE

25

1. Auxiliary carousel is rotated so that the cartridge
is at the MUP station.
2. 50uL of heated MUP are dispensed at a rate of "X"
uL/sec onto the matrix.
- 30 3. Wait "n" seconds.

J. MEIA READ

1. Auxiliary carousel is rotated so that the cartridge
35 is at the read station.

2. The following steps are repeated until the number of micro-reads specified in the assay file have been taken (usually 8).
 - (a) Read for "X.XX" seconds.
 - 5 (b) Wait "X.XX" seconds.
3. The reader is returned to its idle state.
 - (a) Lamp intensity is turned to simmer state.
 - (b) PMT gain is set.
4. The raw reads are converted to normalized reads
10 (light intensity hitting detector/lamp intensity) by the optics microprocessor.
5. A rate is calculated by the System from the normalized reads vs time.
6. For quantitative assays, the rate is fitted to a
15 calibration curve to yield a concentration result.
7. For qualitative assays, the sample rate is compared to an index or cutoff rate to determine if the sample is positive or negative (or reactive or nonreactive).
- 20 K. CARTRIDGE UNLOAD (This activity occurs when resources are not in use)
 1. Auxiliary carousel is rotated so that cartridge is at the ejector station.
 - 25 2. Ejector is cycled to place cartridge into waste container.

Schematic reaction sequences are presented in FIGURES-26, 27 and 28 which are typical of assays that can be handled
30 by the automated immunoassay analytical system of the invention. In FIGURE 26, a T4 assay, FPIA sequence 420, is presented wherein Step 1, T4 bound by thyroxine binding protein (TBP) 424, is reacted with T4 displacing agent 426 to yield TBP 428 plus unbound T4 (430). In step 2, the T4 (430) is
35 added to T4 antibody 432 which yields a reaction product 434 (T4 antibody-T4 complex). In Step 3, the T4 antibody-T4 complex 434 is treated with T4 tracer (fluorescent) 436

which yields a fluorescent polarization measurable reaction product 438.

In FIGURE 27, a schematic reaction sequence 440 for a 1-step sandwich MEIA determination (ferritin) is presented. In Steps 1 and 2 an anti-ferritin alkaline phosphatase conjugate is mixed with ferritin sample 444 and anti-ferritin microparticles 446 of a homogeneous microparticle assay reagent of the present invention to yield a ferritin antibody-antigen-antibody complex 448. In step 3, the antibody-antigen-antibody complex 448 is reacted with 4-methylumbelliferyl phosphate (MUP) 450 which yields methylumbelliferone (MU) which is fluorescent. The rate of MU production is measured.

In FIGURE 28, the schematic reaction sequence 456 for a 2-step sandwich MEIA is provided for HTSH assay. Anti-hTSH specific microparticles 458 of a homogeneous microparticle assay reagent of the present invention are added to the HTSH sample 460 which provides a reaction product HTSH antibody-antigen complex 462. In Steps 2 through 4, the complex 462 is combined with an anti-hTSH alkaline phosphatase 464 yielding hTSH antibody-antigen-antibody complex 466. In step 5, the complex 466 is reacted with MUP 450 to yield MU which is fluorescent. The rate of MU production is measured. In accordance with the embodiments, the automated immunoassay analytical system provides apparatus, software, hardware and process technology for performing a multitude of assays continuously and with random access being available to the operator. The utilization of carousel pipettor technology for kitting and pipetting operations at either the main carousel or the process carousel, depending on the scheduled test, provides scheduling flexibilities heretofore unachievable. The inventive system allows for a commonality of kitting and pipetting for either immuno precipitation or competitive immunoassay technologies utilizing a common main carousel, transfer station, first kitting and pipetting probe and process carousel as well as a second pipetting probe before separating into respective apparatus and process

requirements. Also shared is the commonality of cabinetry disposal and supply materials as well as a common computer network for scheduling, testing, kitting and pipetting.

It will be seen that multiple assays can be performed
5 with a minimum of operator input or handling on the system and the system can be utilized for other processes and assays which have not been directly discussed but will be readily apparent to one practiced in the art in view of the above invention disclosure and the claims. It will also be
10 appreciated that although particular embodiments of the present invention have been disclosed, various changes and adaptations to the apparatus and methods can be made without departing from the teachings of the specification and scope of the invention as set out in the following claims.

What is claimed is:

1. A method of providing a substantially homogeneous liquid assay reagent comprising a liquid component and a particulate assay component, said method comprising the steps of:
 - (a) providing an assay reagent comprising a liquid component and a particulate assay component, wherein the density of said liquid component and the density of said particulate assay component are substantially different; and
 - (b) adding an inert reagent to said liquid assay reagent in an amount whereby the density of said liquid component and the density of said particulate assay component are substantially the same to thereby provide a homogeneous suspension of said particulate assay component.
2. The method of claim 1 wherein said inert material is selected from the group consisting of sucrose, metrizamide, and metrioic acid.
3. The method of claim 1 wherein said particulate assay component comprises a binding protein immobilized thereto.
4. The method of claim 3 wherein said particulate material is selected from the group consisting of beads, particles, and microparticles.
5. The method of claim 3 wherein said binding protein is an antibody or fragment thereof.
6. The method of claim 3 wherein said binding protein is an analyte or analog thereof.
7. The method of claim 1 wherein said homogeneous liquid assay reagent is manually agitated upon expiration of a predetermined period of time or sooner.
8. The method of claim 1 wherein said homogeneous liquid assay reagent is agitated with automated means upon expiration of a predetermined period of time or sooner.

9. A method of providing a substantially homogeneous liquid immunoassay reagent for use in a heterogeneous immunoassay, said liquid immunoassay reagent comprising a liquid component and an immobilized binding protein component, said method comprising the steps of:

(a) providing an immunoassay reagent comprising a liquid component and an immobilized protein component, said immobilized protein component comprising a particulate material, wherein the density of said liquid component and the density of said immobilized protein component are substantially different; and

(b) adding an inert reagent to said liquid assay reagent in an amount whereby the density of said liquid component and the density of said immobilized protein component are substantially the same to thereby provide a substantially homogeneous suspension of said immobilized binding protein component.

10 The method of claim 9 wherein said inert material is selected from the group consisting of sucrose, metrizamide, and metrcic acid.

11. The method of claim 9 wherein said binding protein is an antibody or fragment thereof.

12. The method of claim 9 wherein said binding protein is an analyte or analog thereof.

13. The method of claim 9 wherein said particulate material is selected from the group consisting of beads, particles, and microparticles.

14. The method of claim 9 wherein said heterogeneous immunoassay is a competitive heterogeneous immunoassay.

15. The method of claim 9 wherein said heterogeneous immunoassay is a sandwich immunoassay.

16. The method of claim 9 wherein said homogeneous liquid assay reagent is manually agitated upon expiration of a predetermined period of time or sooner.

17. The method of claim 9 wherein said homogeneous liquid assay reagent is agitated with automated means upon expiration of a predetermined period of time or sooner.

18. A method for performing a heterogeneous immunoassay for determining the presence or amount of an analyte or an antibody present in a test sample, said method comprising the steps:

5 (a) contacting said test sample with a labeled reagent comprising the one or the other of (i) said analyte or analog thereof or (ii) said antibody, labeled with a detectable moiety, to form a free species and a bound species thereof; and

(b) separating said free species from said bound species
10 with an immunoassay reagent comprising a liquid component and an immobilized component, said immobilized component comprising (i) the one or the other of said analyte or analog thereto or said antibody immobilized to a particulate material and (ii) an inert reagent,

15 wherein said inert reagent is present in said immunoassay reagent in an amount whereby the density of said liquid component and the density of said immobilized component are substantially the same to thereby provide a homogeneous suspension of said immobilized component.

20 19. The method of claim 18 wherein said inert material is selected from the group consisting of sucrose, metrizamide, and metroic acid.

25 20. The method of claim 18 wherein said particulate material is selected from the group consisting of beads, particles, and microparticles.

21. The method of claim 18 wherein said homogeneous liquid assay reagent is manually agitated upon expiration of a predetermined period of time or sooner.

30 22. The method of claim 19 wherein said homogeneous liquid assay reagent is agitated with automated means upon expiration of a predetermined period of time or sooner.

23. A method of operating an automated, continuous and random access analytical system capable of simultaneously effecting multiple assays of a plurality of liquid samples, comprising:
- 5 a. scheduling various assays of a plurality of liquid samples;
 - b. creating one or more unit dose disposables by separately transferring a first of said liquid samples and one or more reagents to a reaction vessel without
10 initiation of an assay reaction sequence, wherein said one or more reagents is a homogeneous liquid assay reagent comprising a liquid component and a particulate assay component and (ii) an inert reagent, wherein said inert reagent is present in said particulate assay reagent in an amount
15 whereby the density of said liquid component and the density of said particulate assay component are substantially the same to thereby provide a homogeneous suspension of said particulate component;
 - c. transferring one or more of said unit dose
20 disposables to a processing workstation;
 - d. mixing an aliquot of said first liquid sample with one or more of said reagents at different times in said reaction vessel to form a first reaction mixture;
 - e. mixing aliquots of the same or different one or
25 more samples with one or more of said reagents at different times in different reaction vessels to form multiple, independently scheduled reaction mixtures;
 - f. incubating said multiple reaction mixtures simultaneously and independently;
 - 30 g. Performing more than one scheduled assay on said reaction mixtures in any order in which more than scheduled assays are presented; and
 - h. analyzing said incubated reaction mixtures independently and individually by at least two assay
35 procedures.

24. The method of claim 23 wherein said inert material is selected from the group consisting of sucrose, metrizamide, and metrcic acid.

5 24. The method of claim 23 wherein said particulate assay reagent comprises a binding protein immobilized to a particulate material.

25. The method of claim 24 wherein said particulate material is selected from the group consisting of beads, particles, and microparticles.

10 26. The method of claim 24 wherein said binding protein is an antibody or fragment thereof.

27. The method of claim 24 wherein said binding protein is an analyte or analog thereof.

15 28. The method of claim 23 wherein said automated, continuous and random access analytical system comprises a front end carousel rotatably disposed and capable of bidirectional rotational motion to provide a bidirectional shaking motion for agitating said homogeneous liquid assay reagent upon expiration of a predetermined period of time
20 period or sooner.

29. The method of claim 23 wherein at least two different assays are scheduled to be performed on the system for a plurality of liquid samples, said method provides scheduling of said assays in advance of performance thereof,
25 each assay test definition containing several timing parameters with each activity of the assay test containing time values which scheduling uses to determine which resources of the system and activity required by each of said assays and the time period that said resources need.

30 30. The method of claim 23 wherein the scheduling process includes scheduling of each activity to be performed in an assay before the assay is kitted, and scheduling performance of each assay activity prior to its originally scheduled execution time, thus minimizing resource idle time.

35 31. The method of claim 23 wherein assay throughput is increased in the system.

32. The method of claim 23 wherein

operating the automated continuous and random access analytical system includes kitting a unit dose disposable by separately transferring assay samples and reagents to a reaction vessel without initiation of an assay reaction sequence.

33. The method of claim 32 wherein the system is capable of allowing special priority handling through stat procedure scheduling of specific samples, and wherein said stat procedure scheduling interrupts prior scheduling, allowing the system to finish preparing assays on a current sample and then prepare to perform an assay on the sample through a modification of the scheduling.

34. The method of claim 32 wherein scheduling for performing assays maximizes the number of assays the system is capable of processing per unit time by allowing sufficient time gaps between assay protocol steps to enable other assay protocol steps to be performed within such time gaps.

35. The method of claim 34 wherein calibration procedures scheduling is scheduled as a stat procedure.

36. The method of claim 23 wherein the assay performed on said reaction mixture in said reaction vessel is a homogeneous assay.

37. The method of claim 23 wherein the assay performed on said reaction mixture in said reaction vessel is a heterogeneous assay.

38. The method of claim 23 wherein at least two assays are immunoassays.

39. The method of claim 38 wherein said immunoassays are comprised of MEIA and FPIA assays.

40. The method of claim 23 wherein said analyzing step includes optically monitoring said reaction mixtures.


41. The method of claim 23 wherein said reaction mixtures are monitored by turbidimetric, colorimetric, fluorometric or luminescent means.


42. The method of claim 23 wherein partial initiation of an assay reaction sequence is achieved simultaneously with the creating of the unit dose disposable.

43. The method of claim 23 wherein the system
5 achieves simultaneous creation of unit dose disposables, transferring of a unit dose disposable reaction vessel and mixing of a reaction mixture while incubating multiple reaction mixtures and performing at least one scheduled assay and analysis simultaneously.

10

1/38

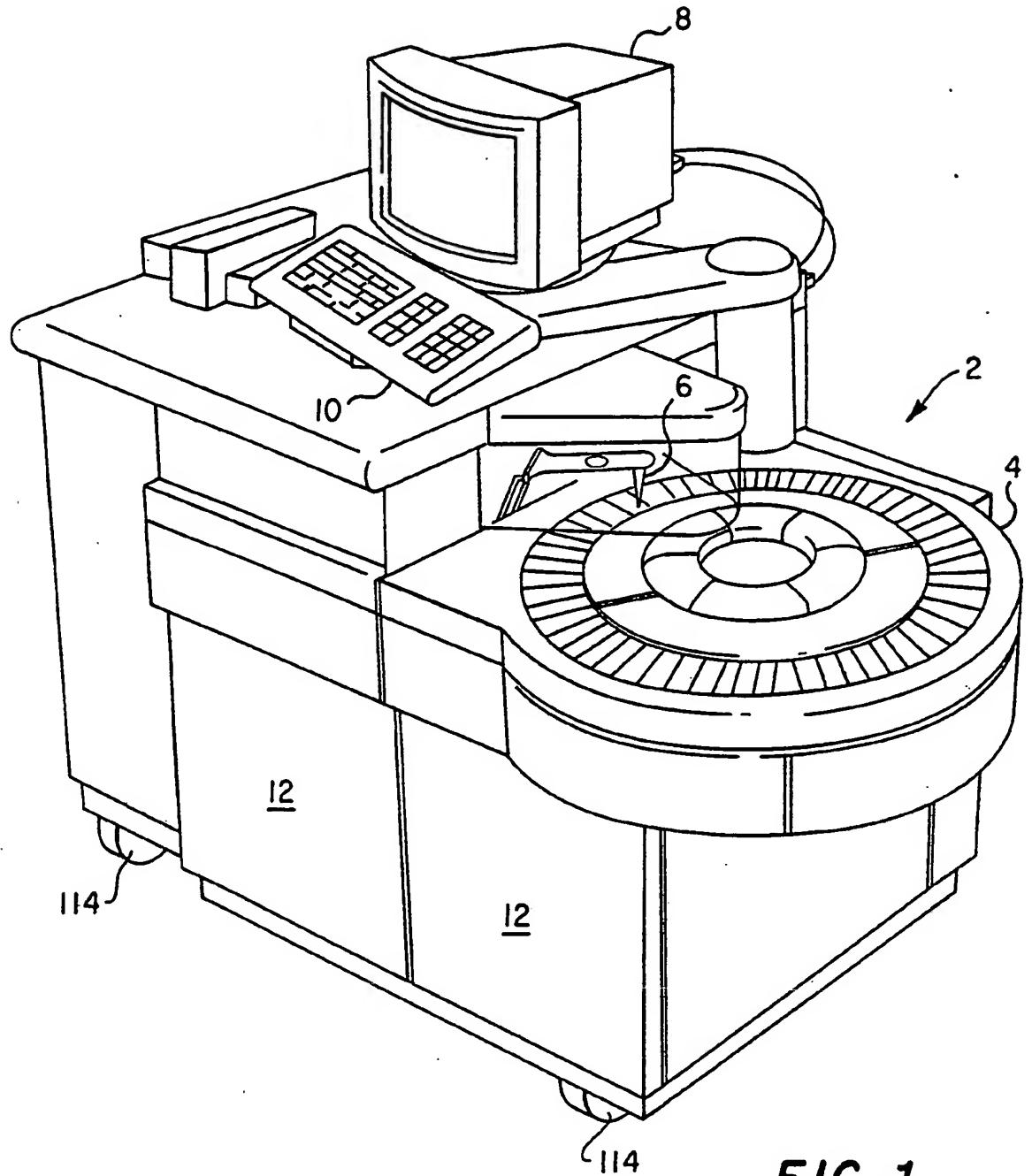**FIG. 1**

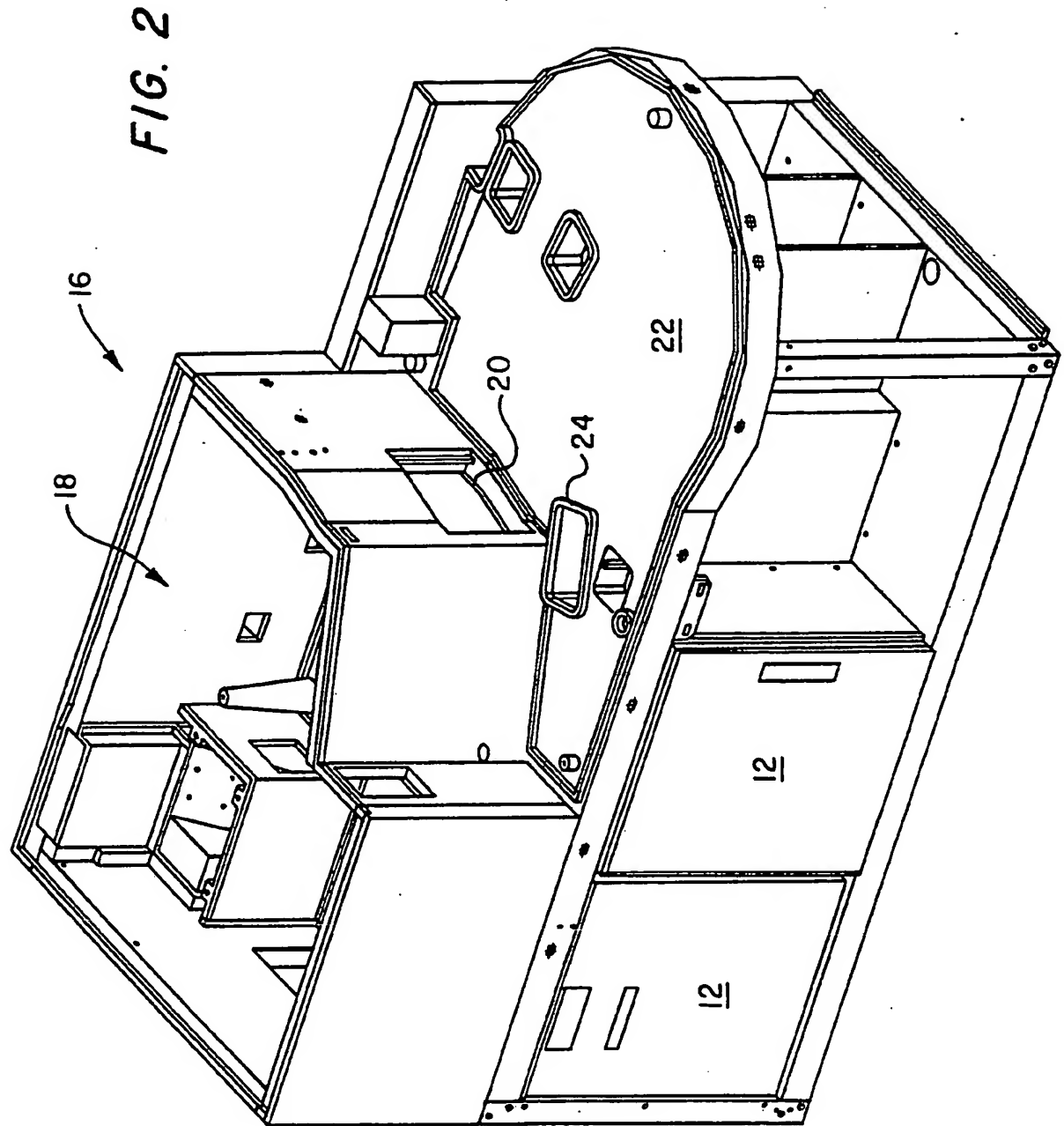

3/38

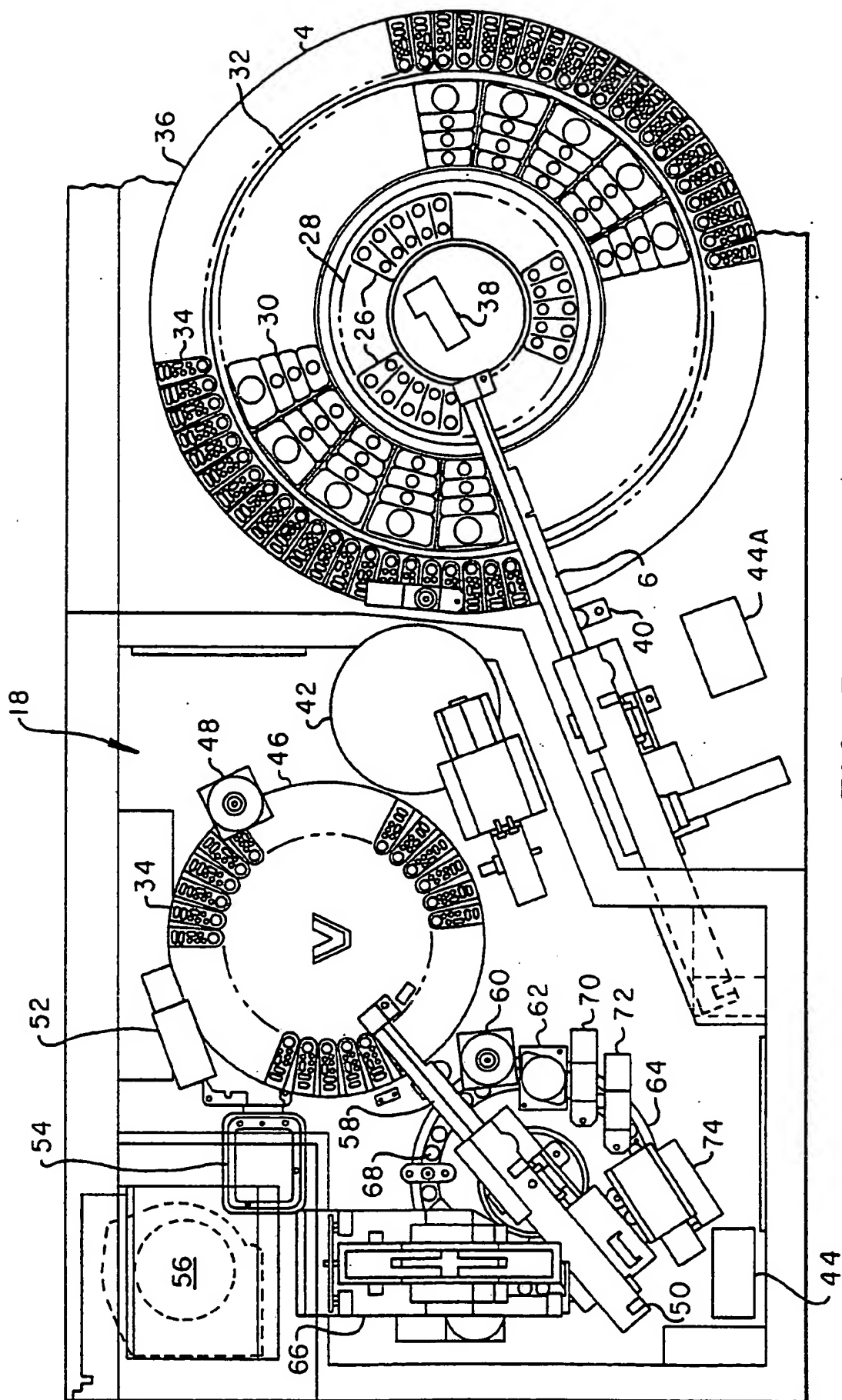


FIG. 3

4 / 38

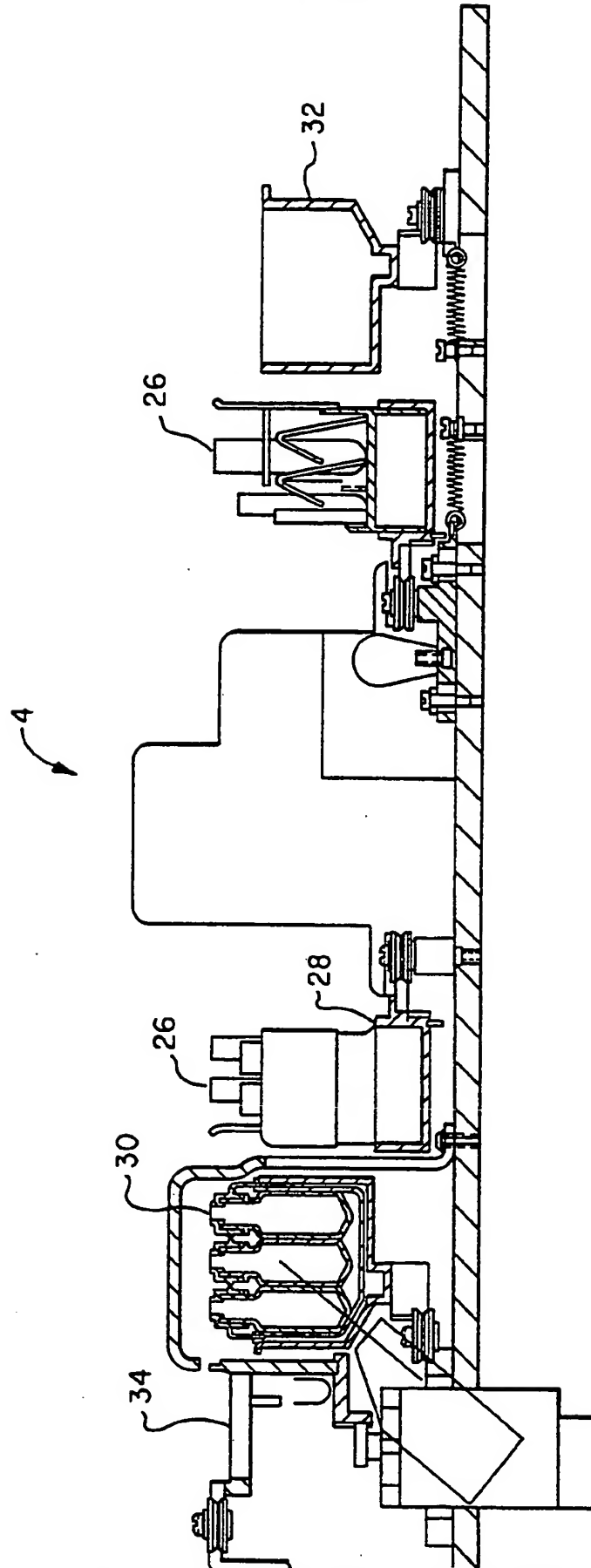


FIG. 4

5/38

FIG. 4A

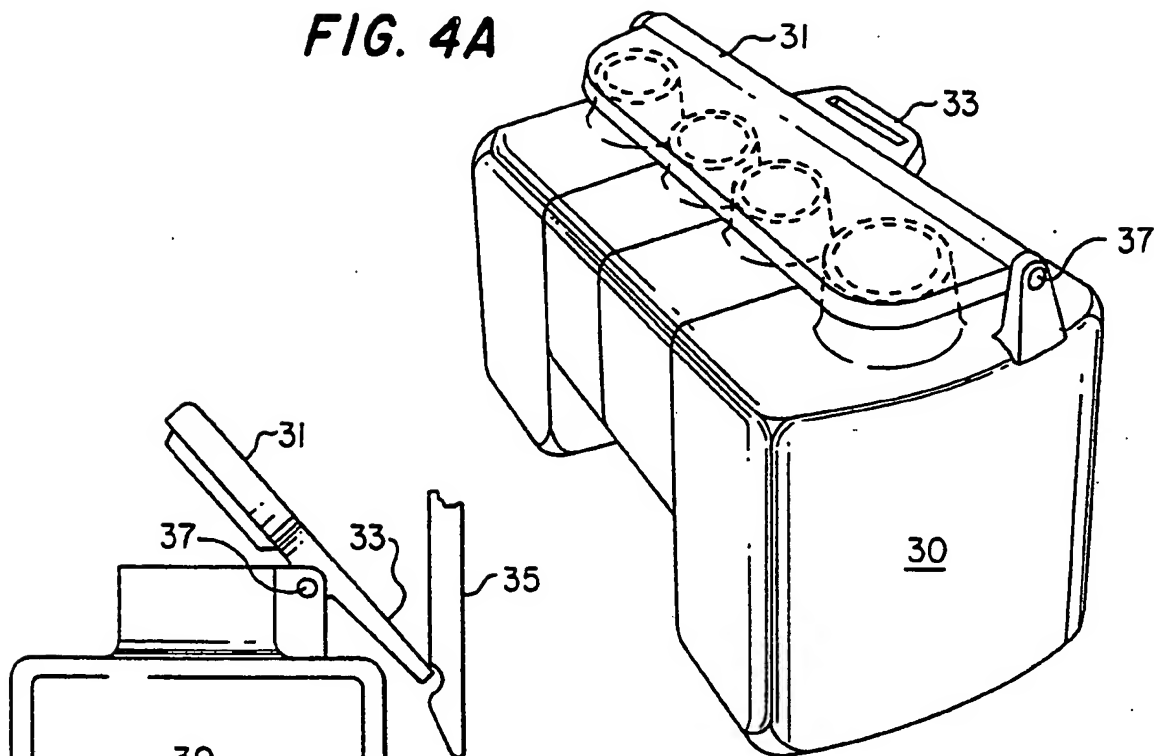


FIG. 4B



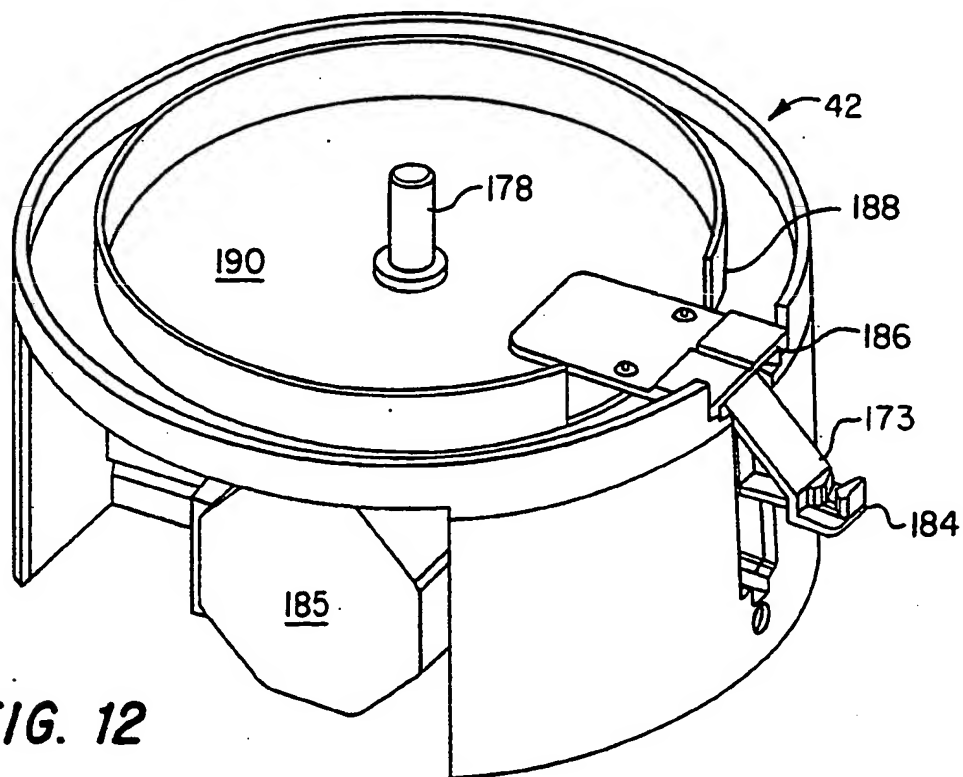



FIG. 12

6/38


FIG. 5

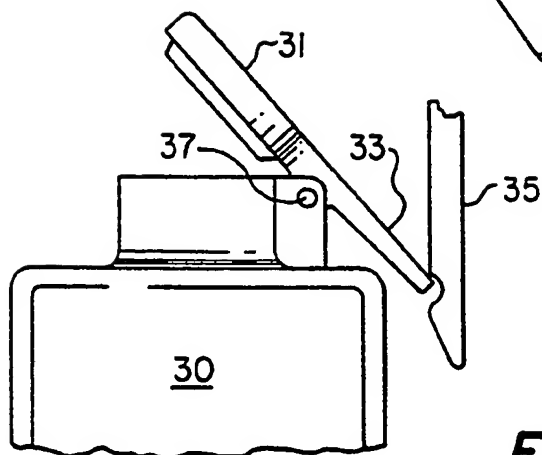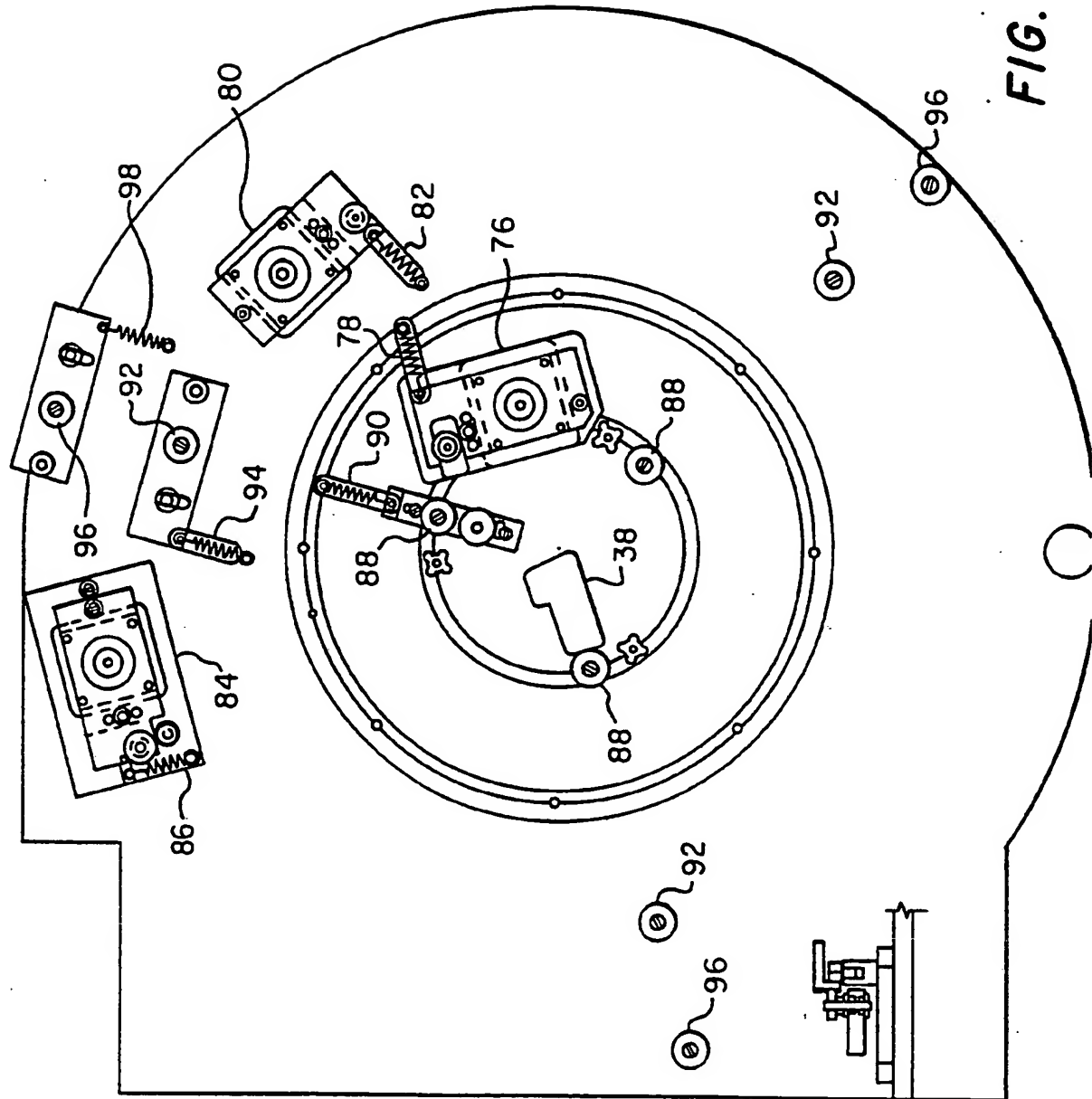

7/38

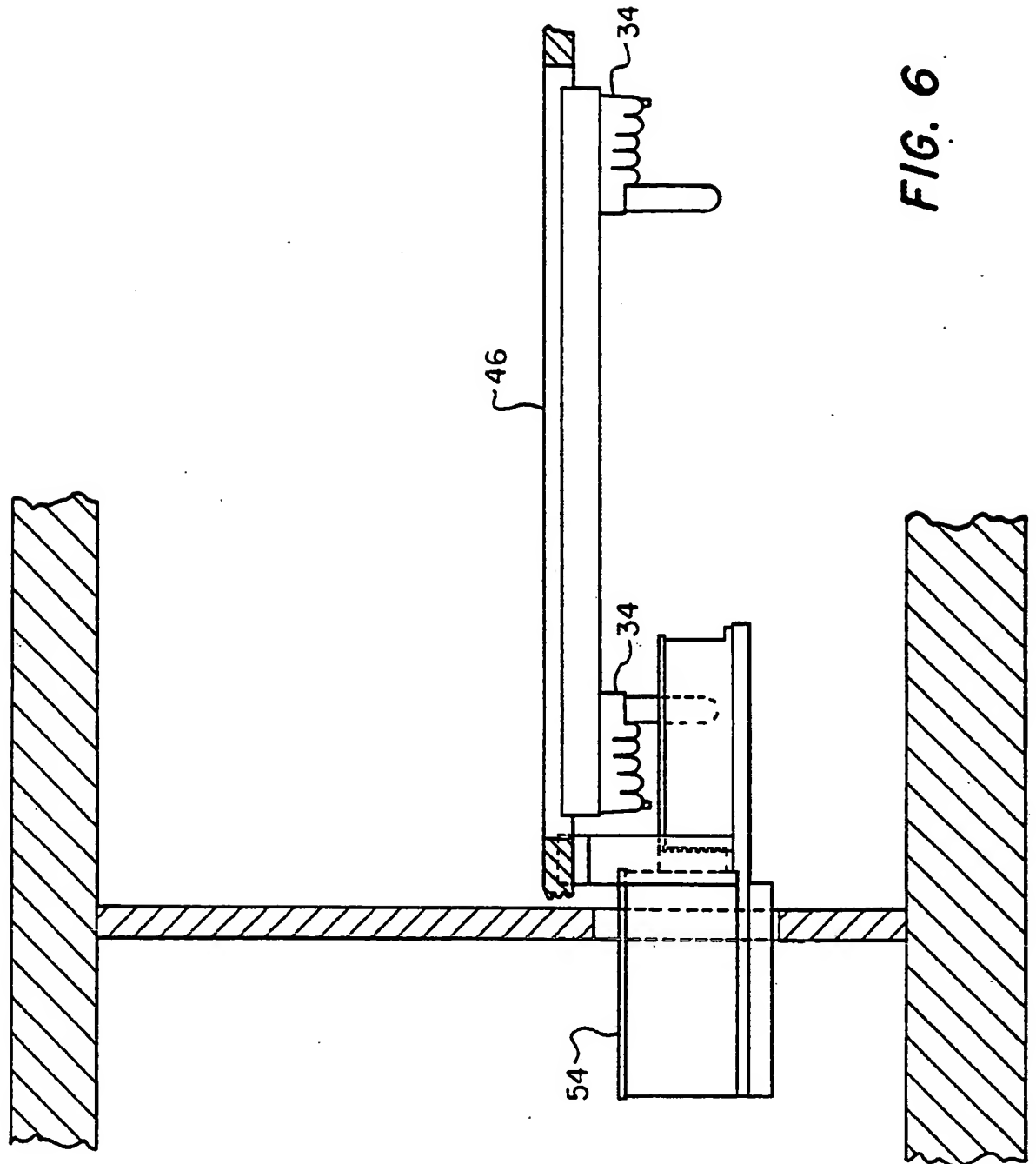

8 / 38

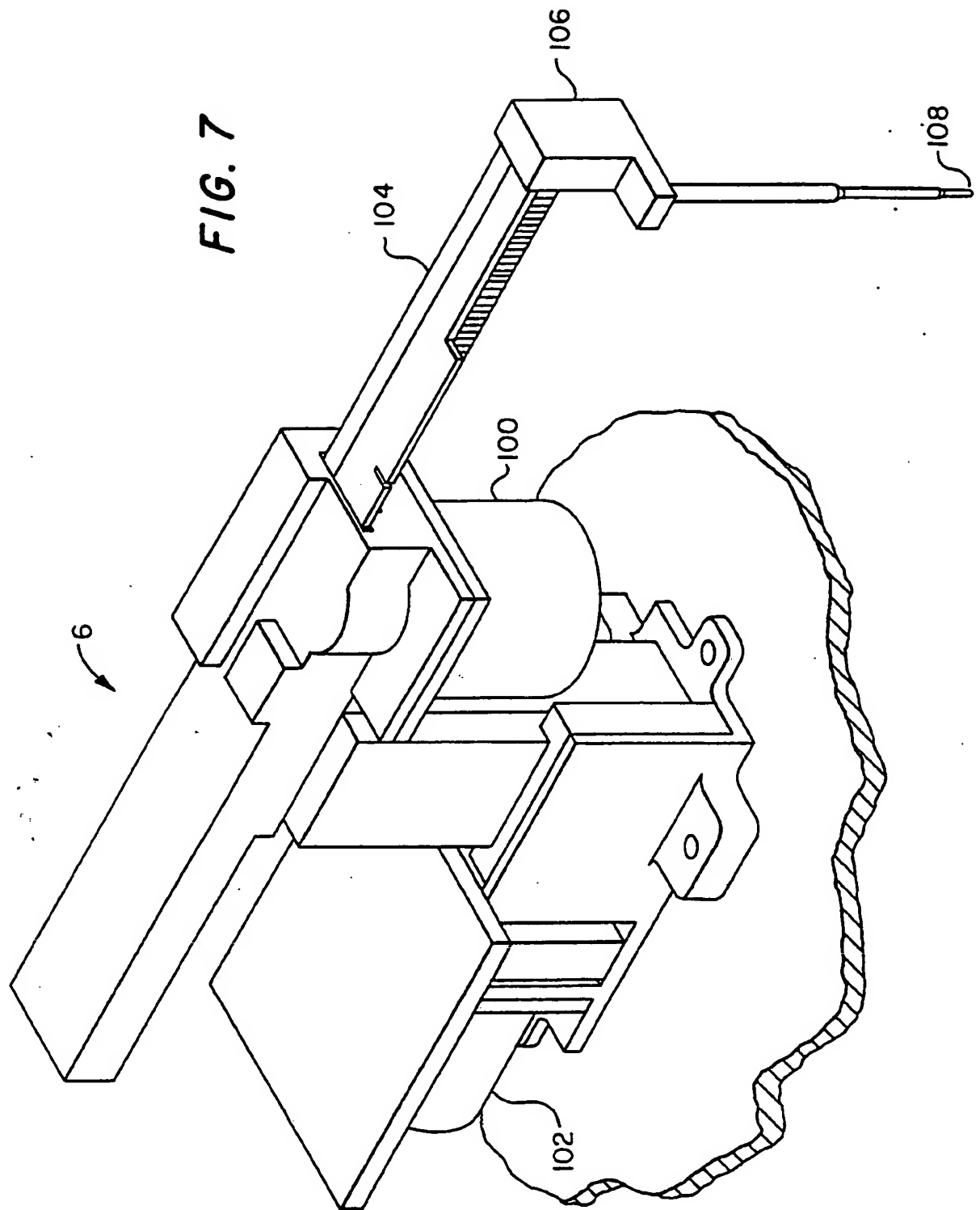

9/38

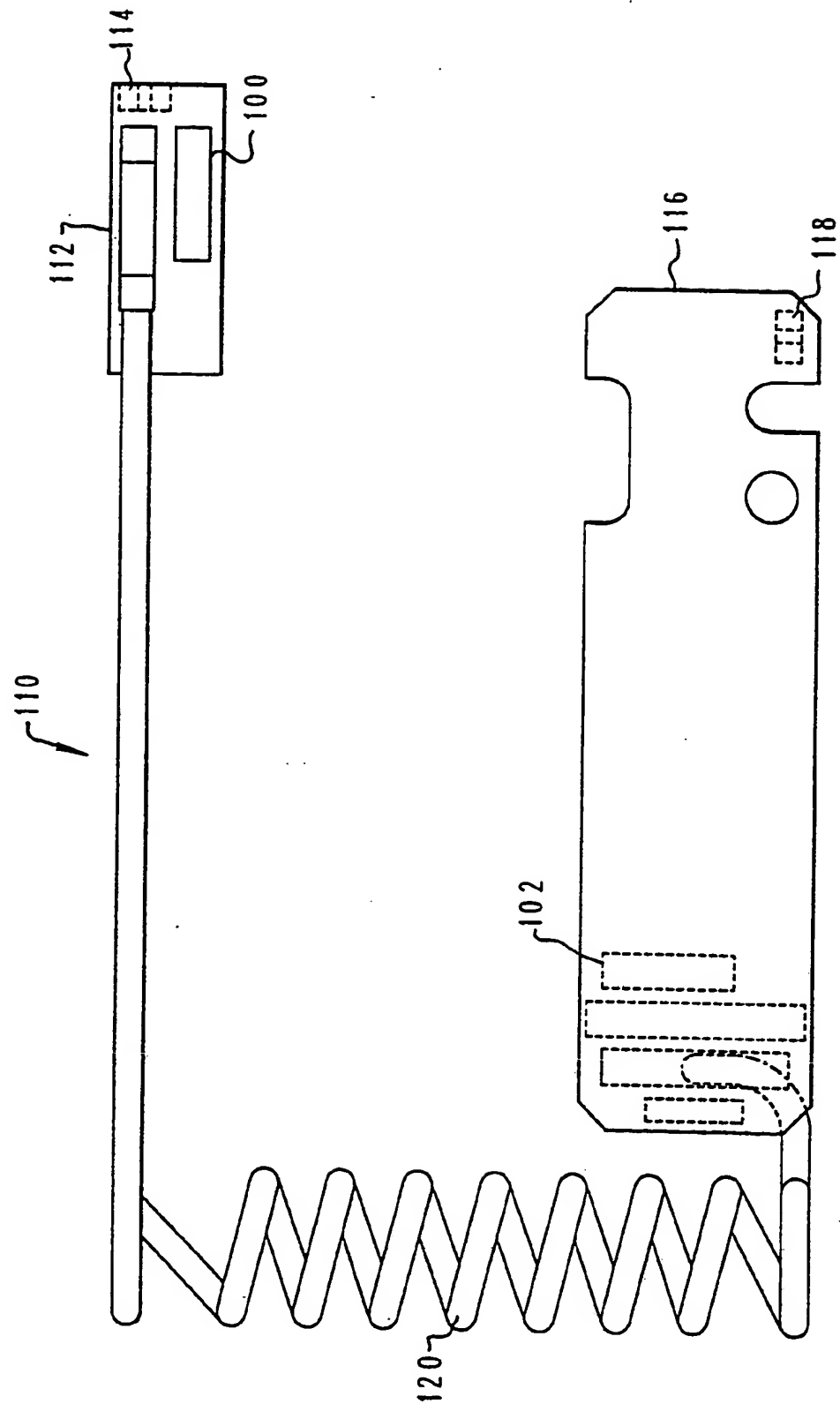


Fig. 8

10/38

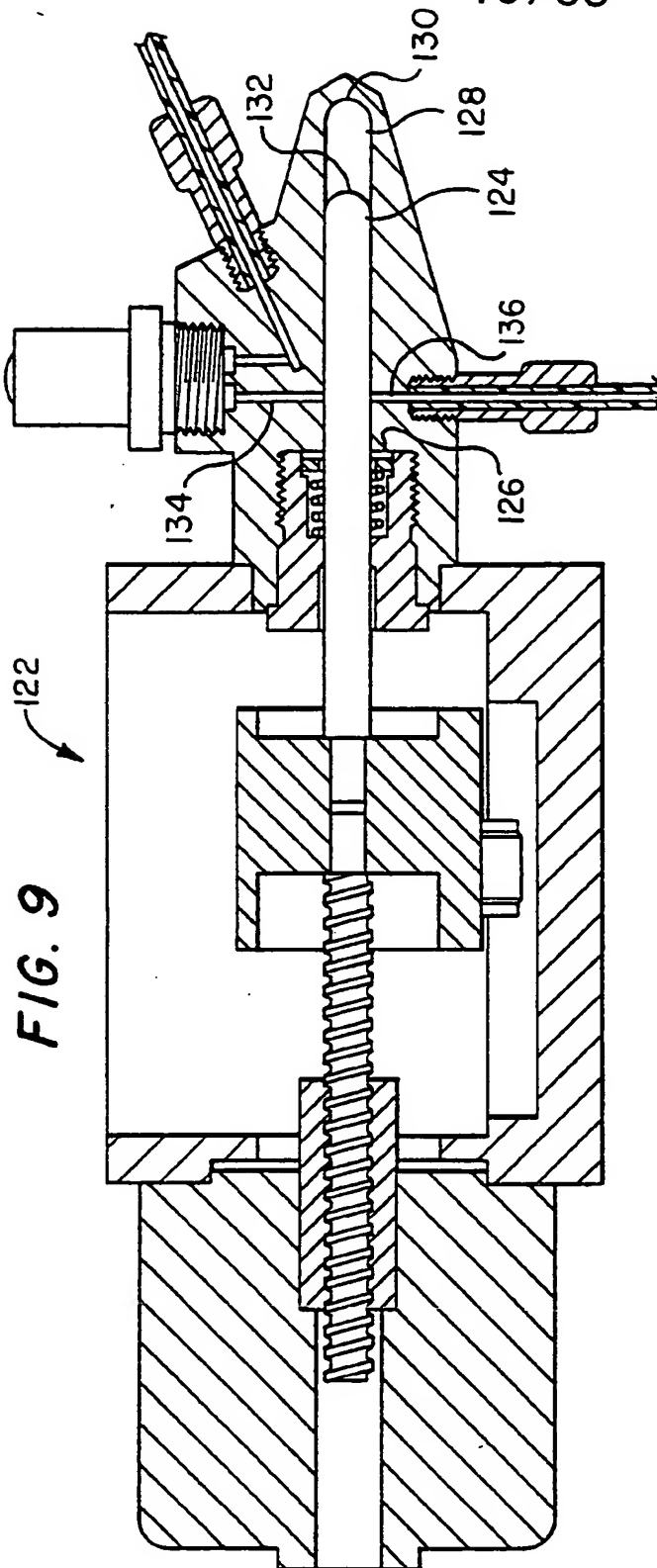


FIG. 9

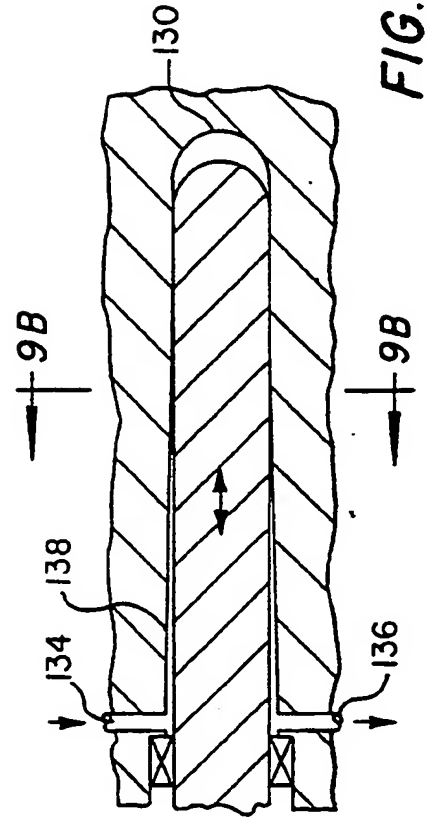


FIG. 9A

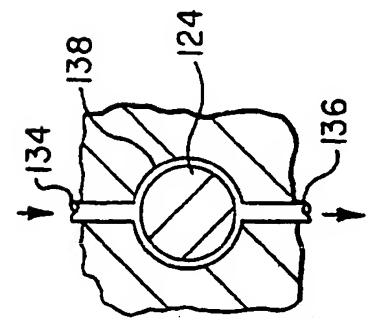


FIG. 9B

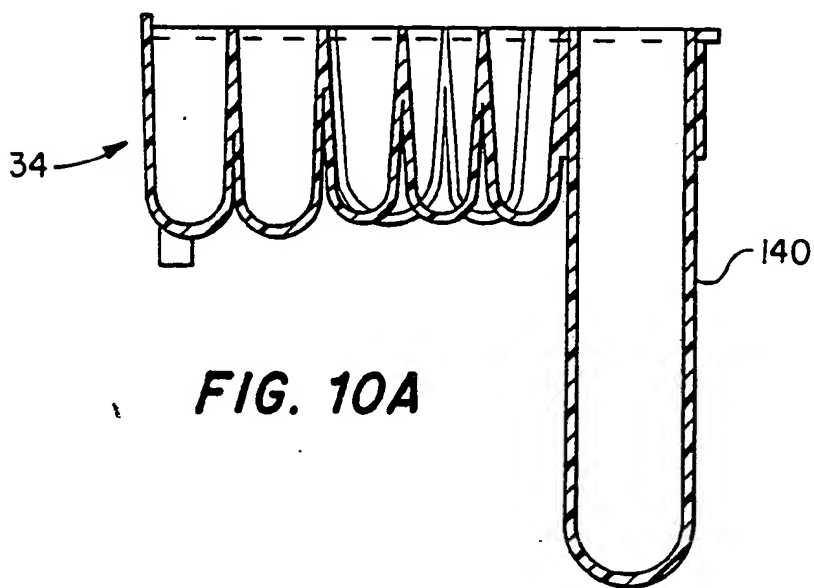
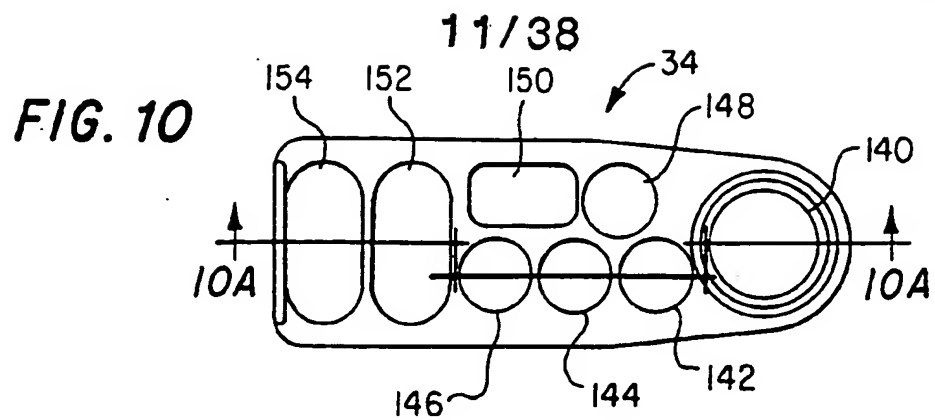


FIG. 10A

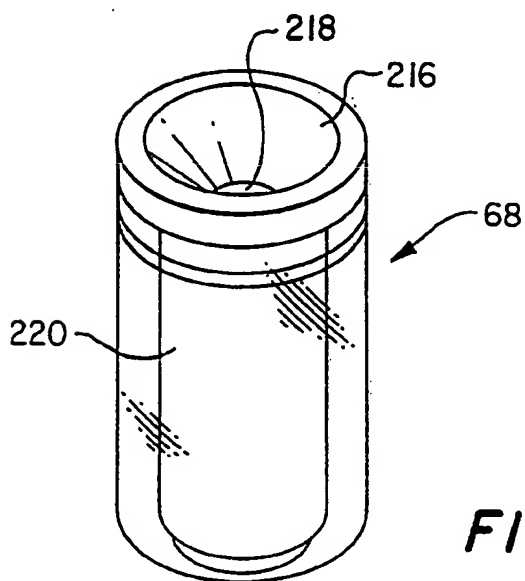


FIG. 16

12 / 38

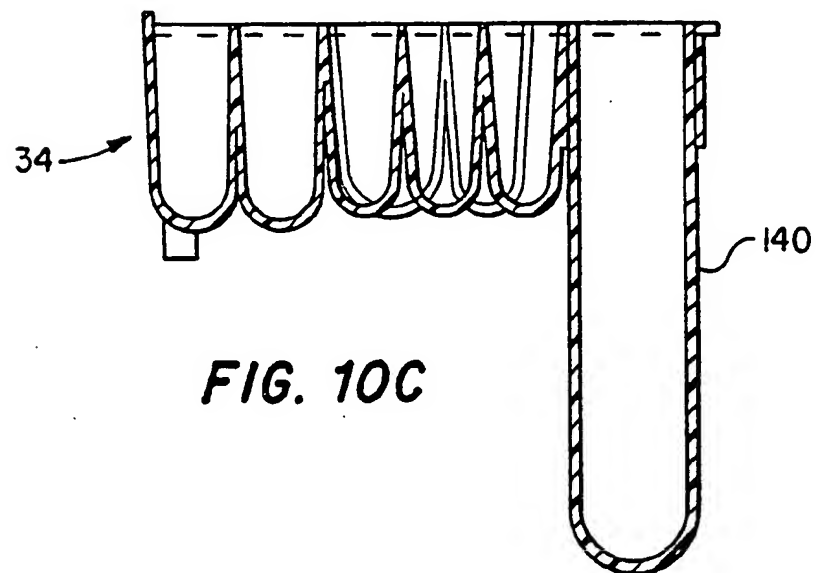
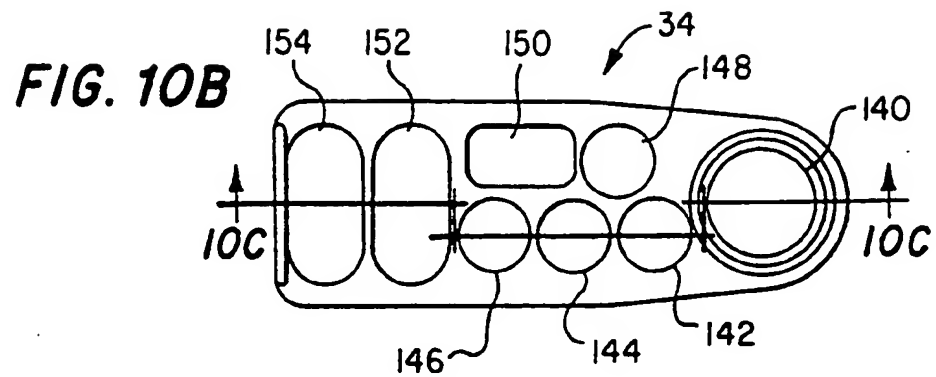


FIG. 10C

13/38

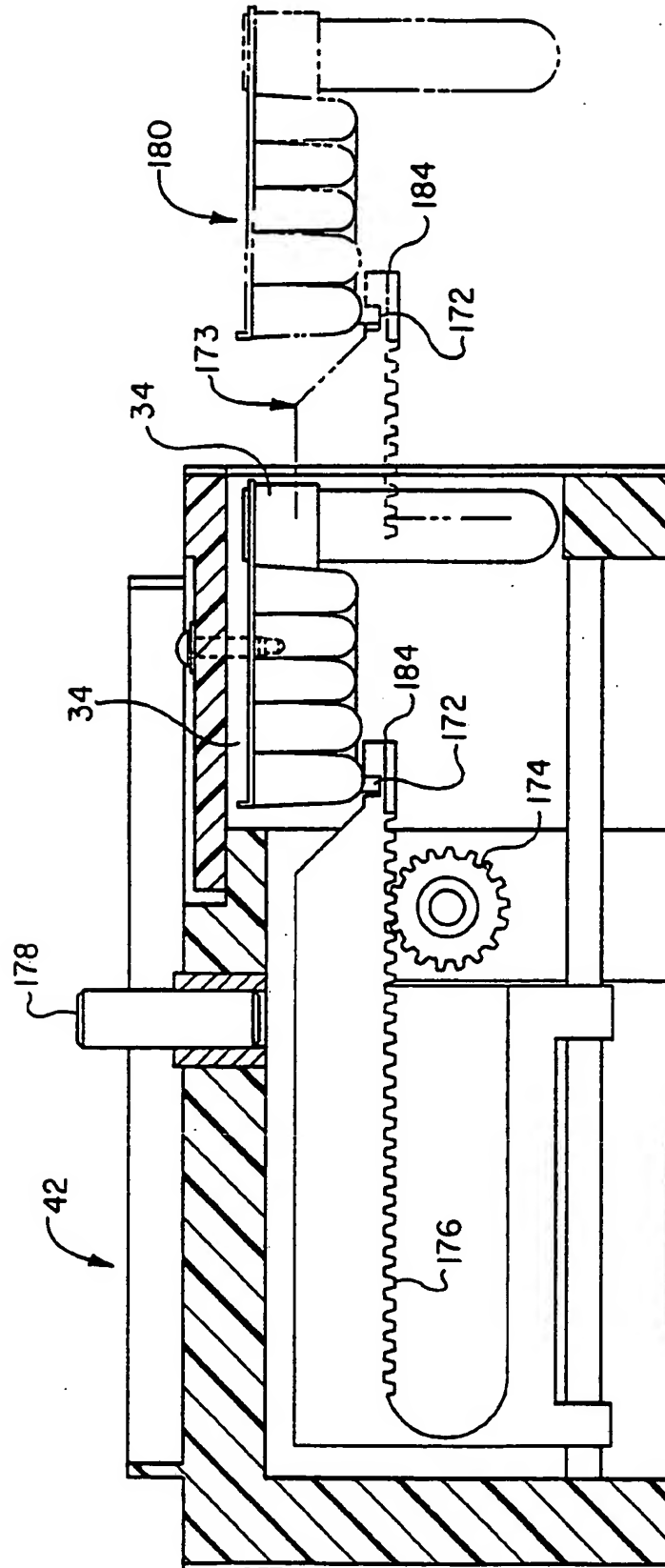


FIG. 11

14/38

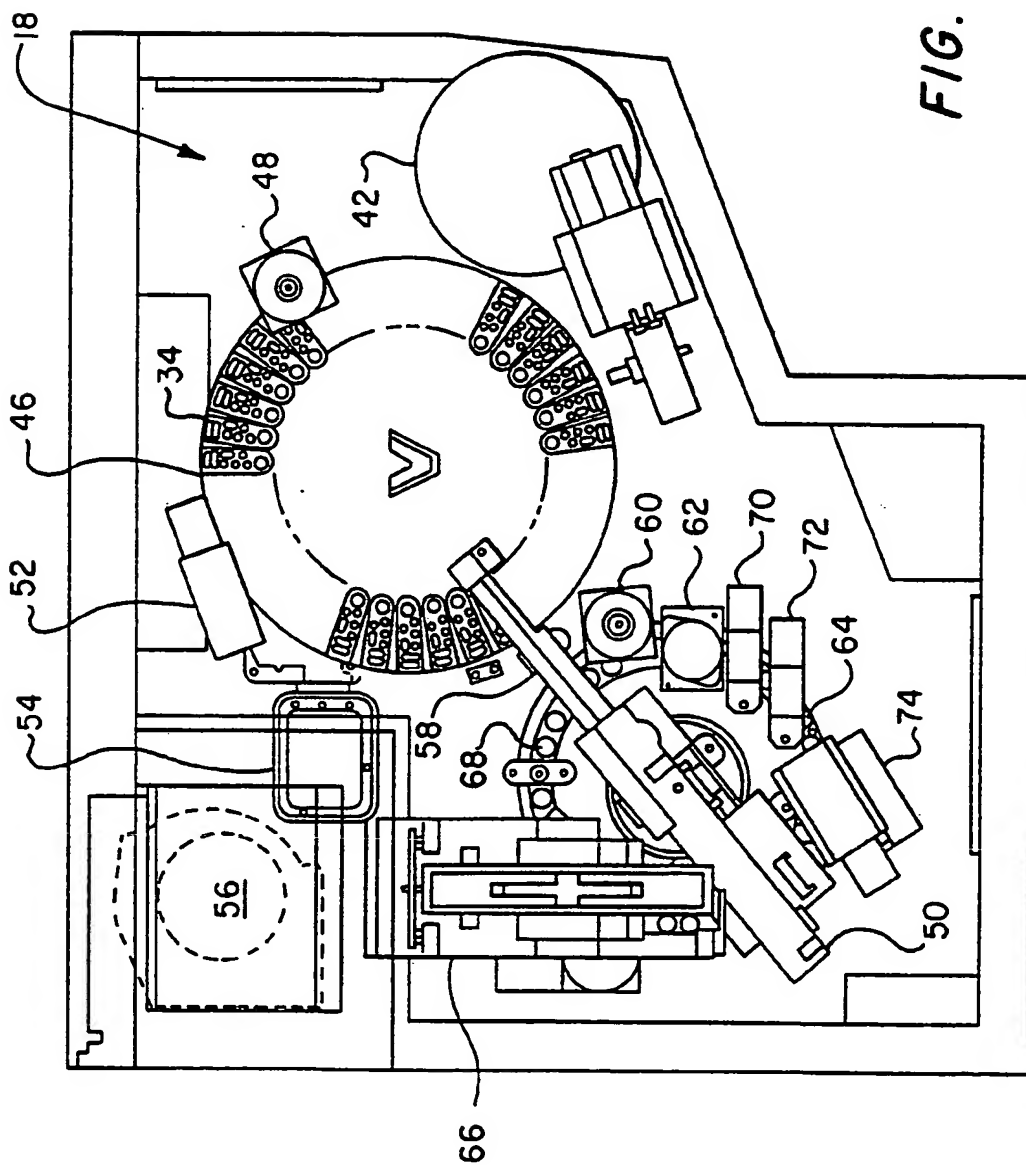
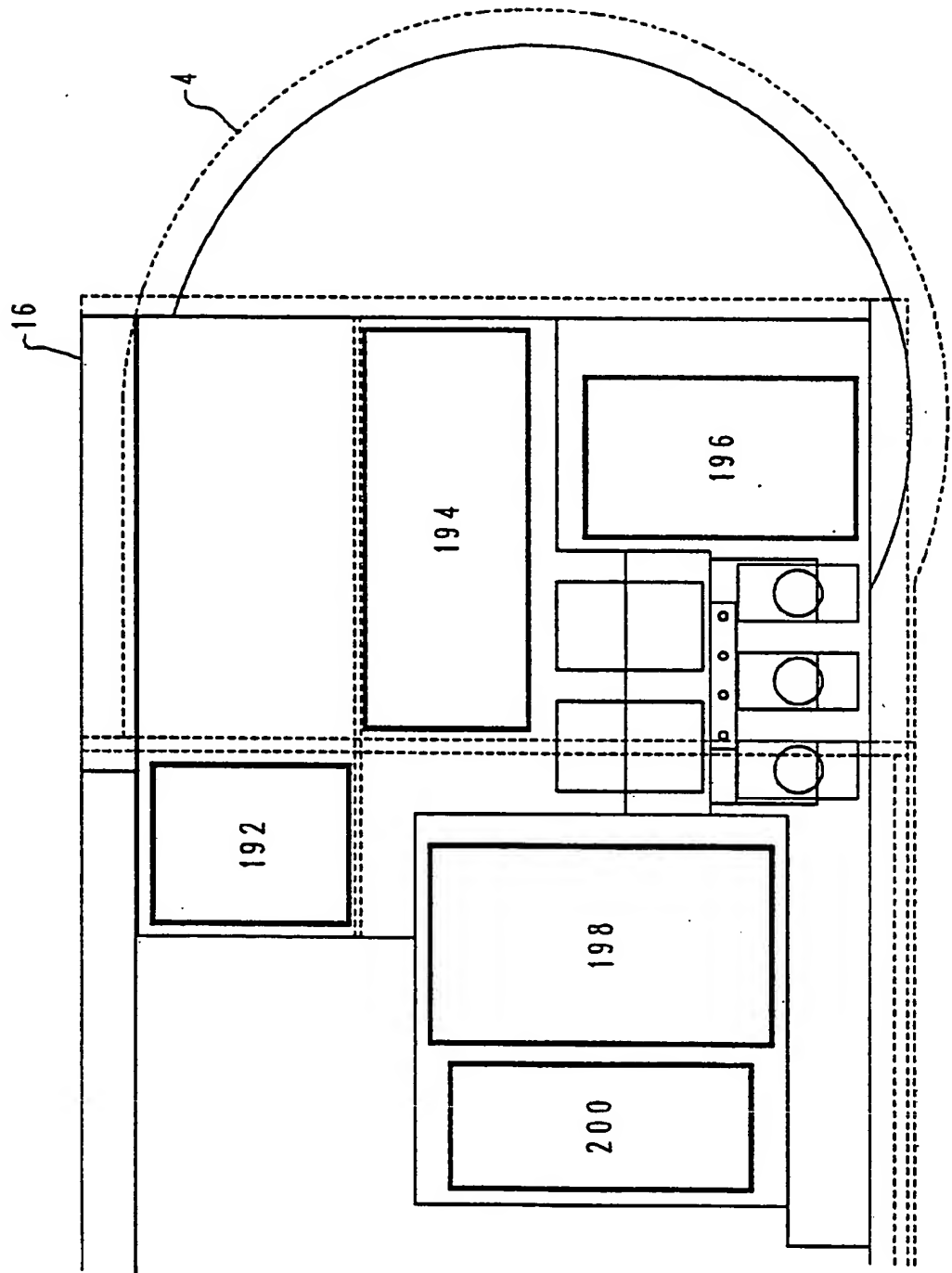


FIG. 13

15/38

*Fig. 14*


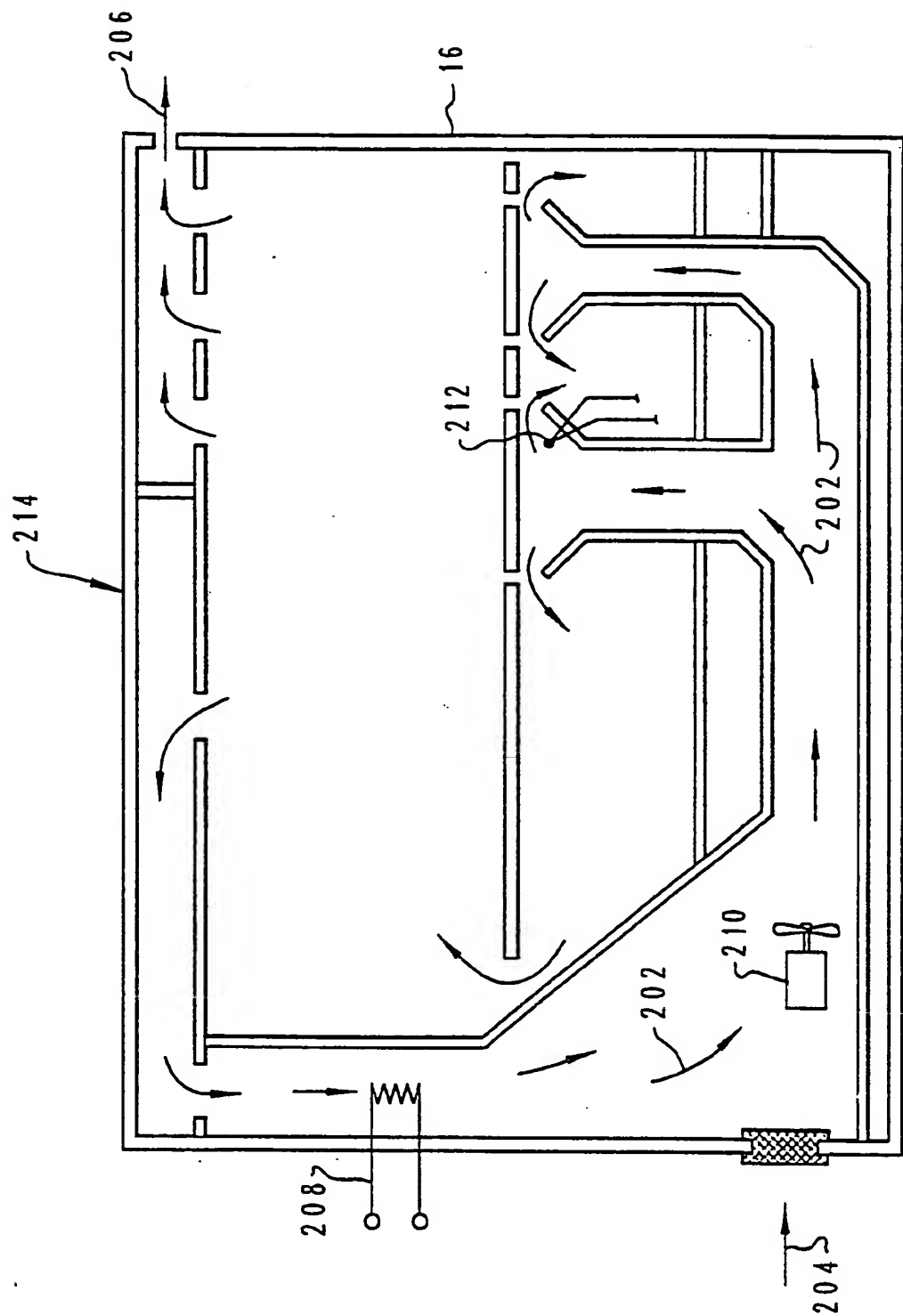


Fig. 15

17 / 38

FIG. 17

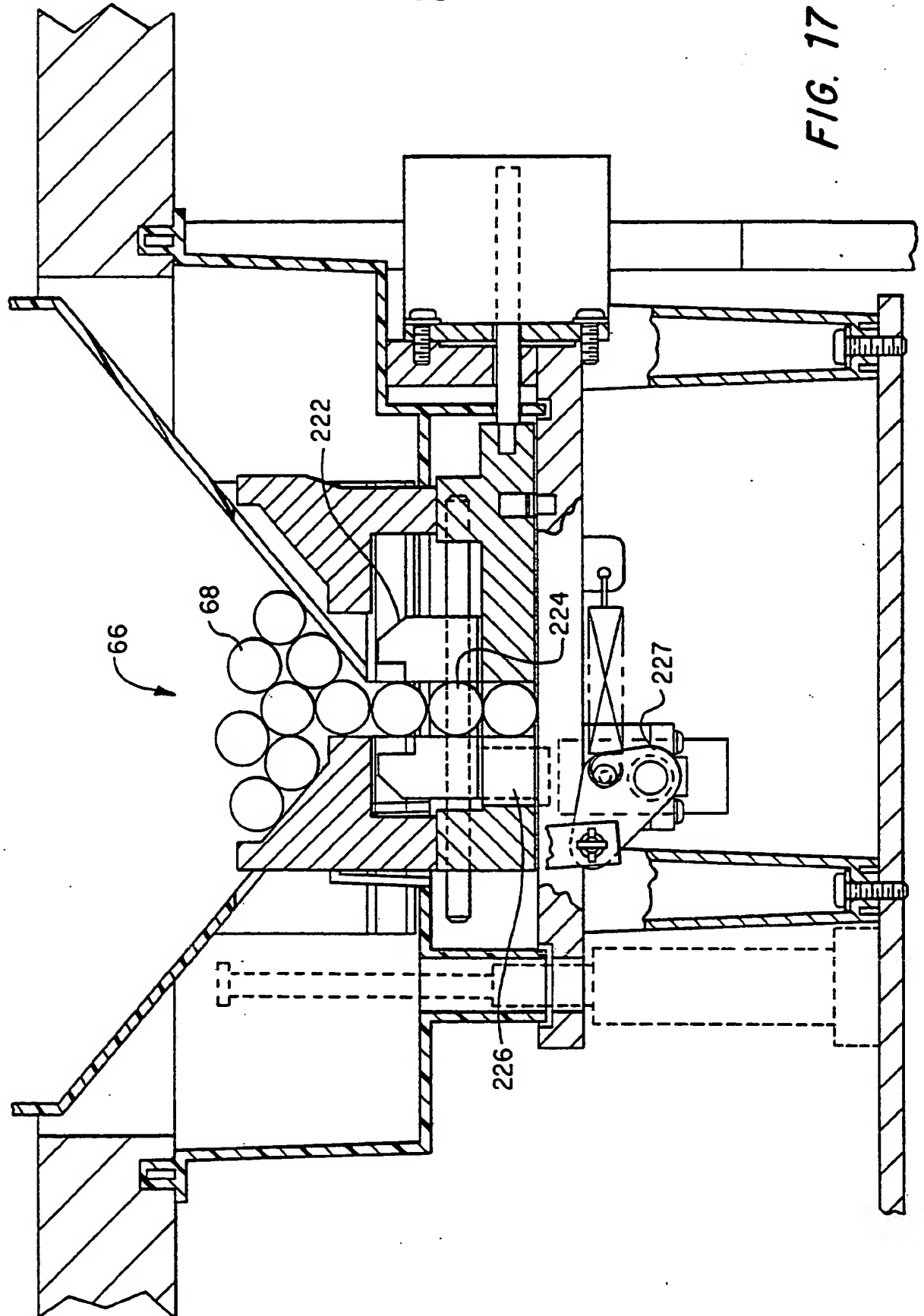

18 / 38

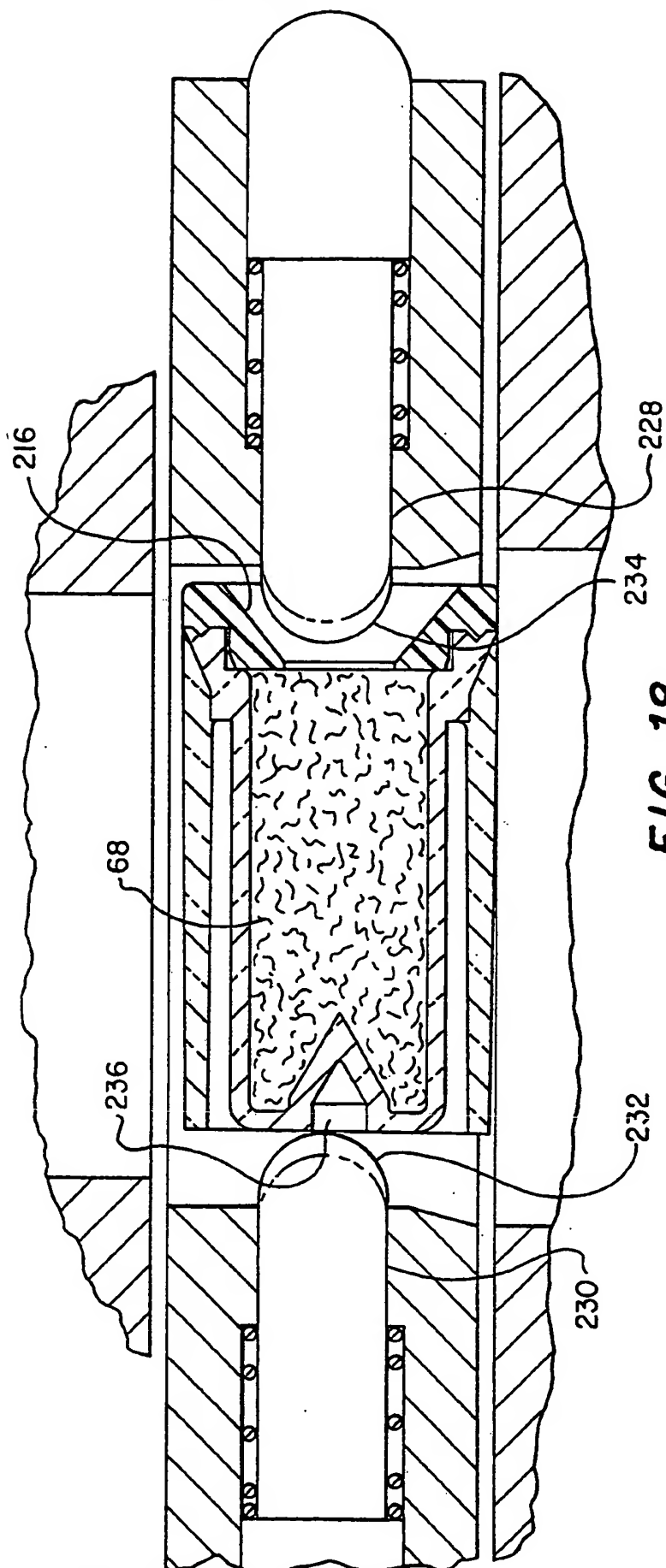


FIG. 18

19/38

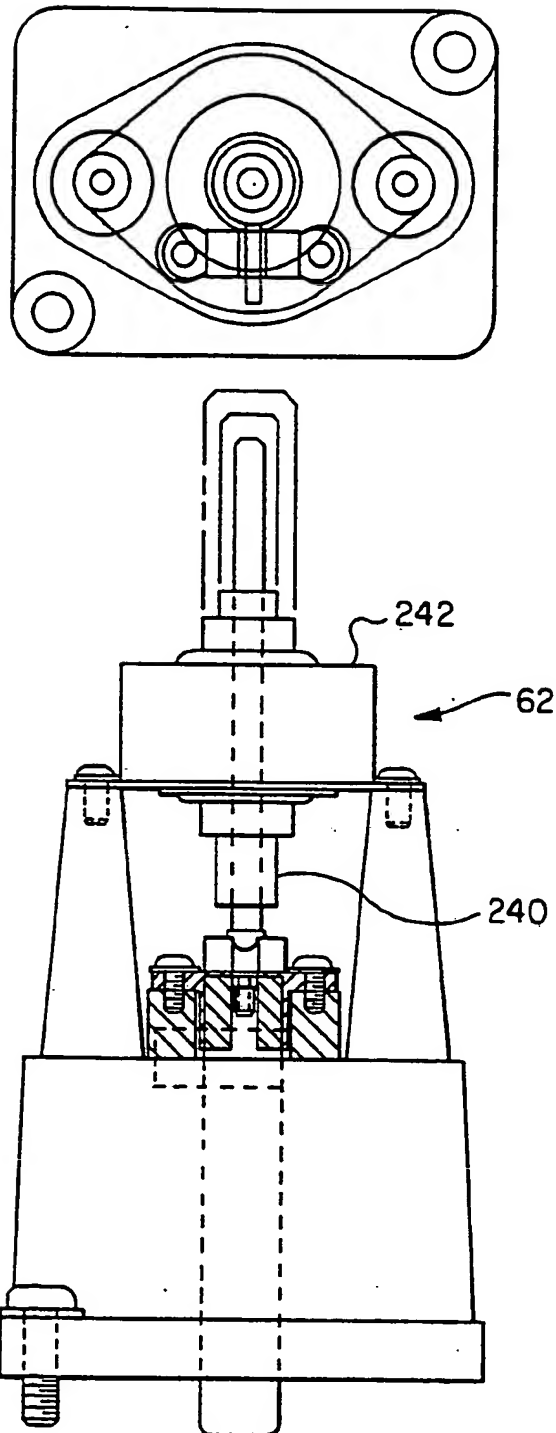


FIG. 19

20/38

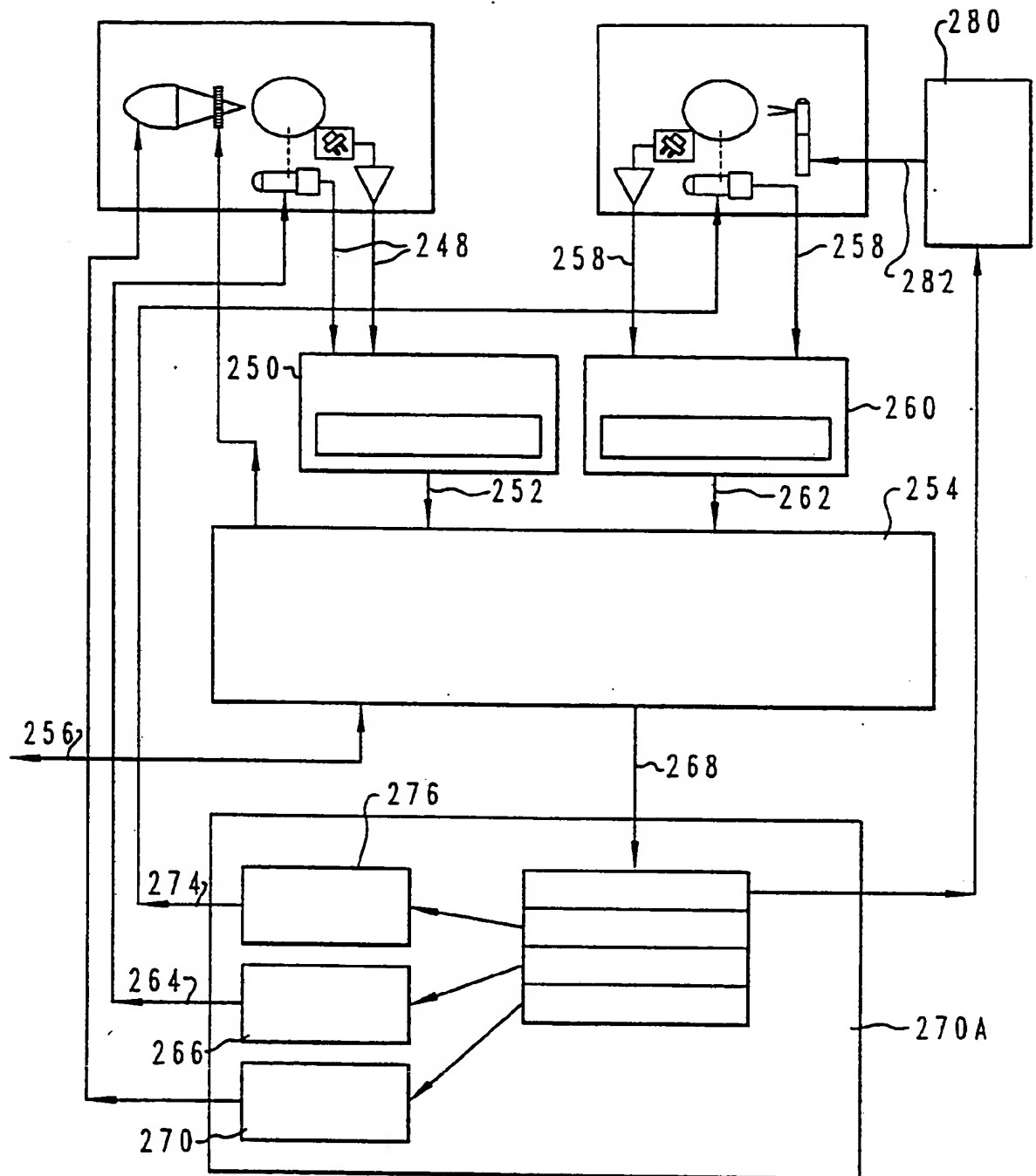


Fig. 20

21/38

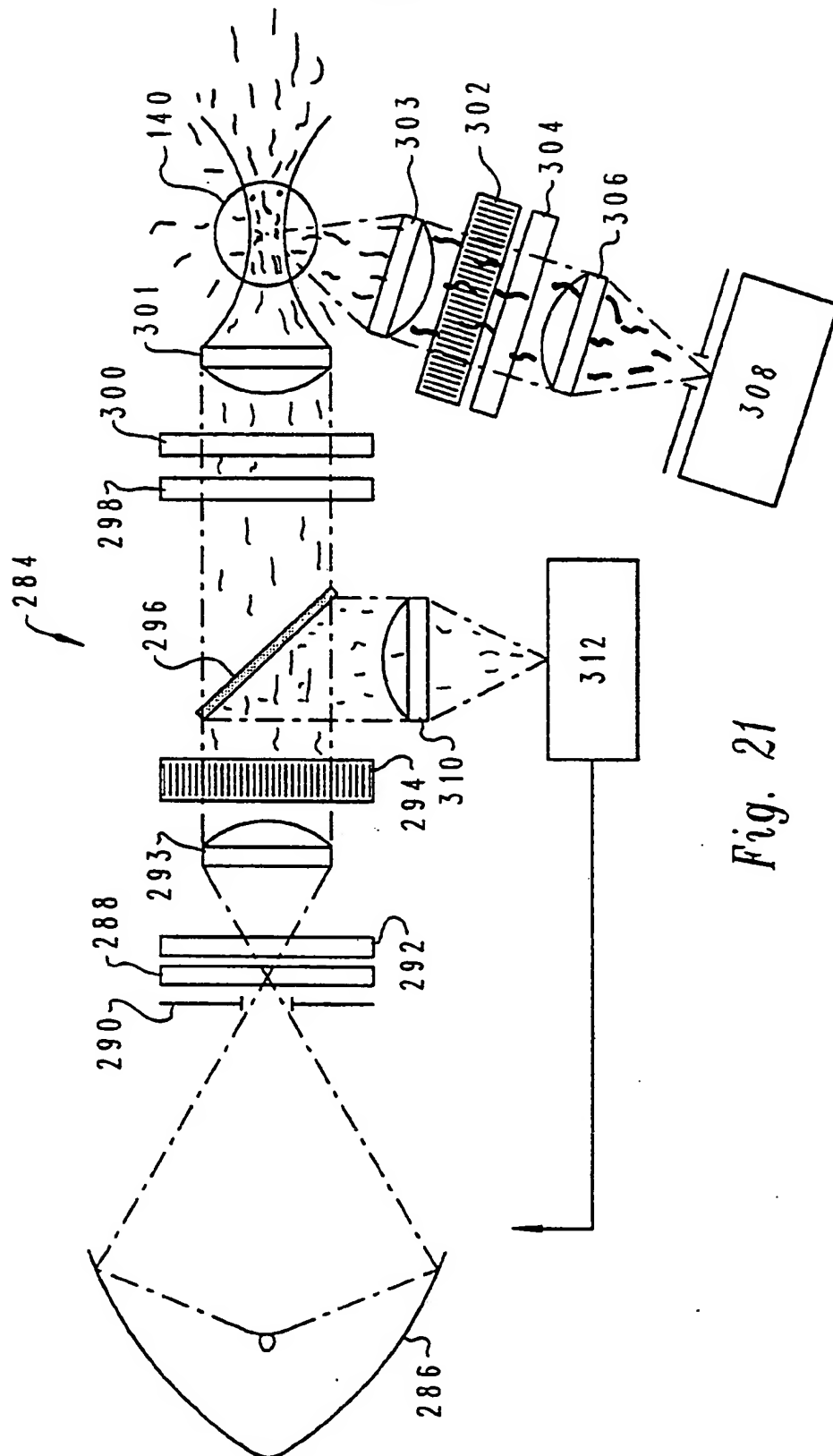


Fig. 21

22/38

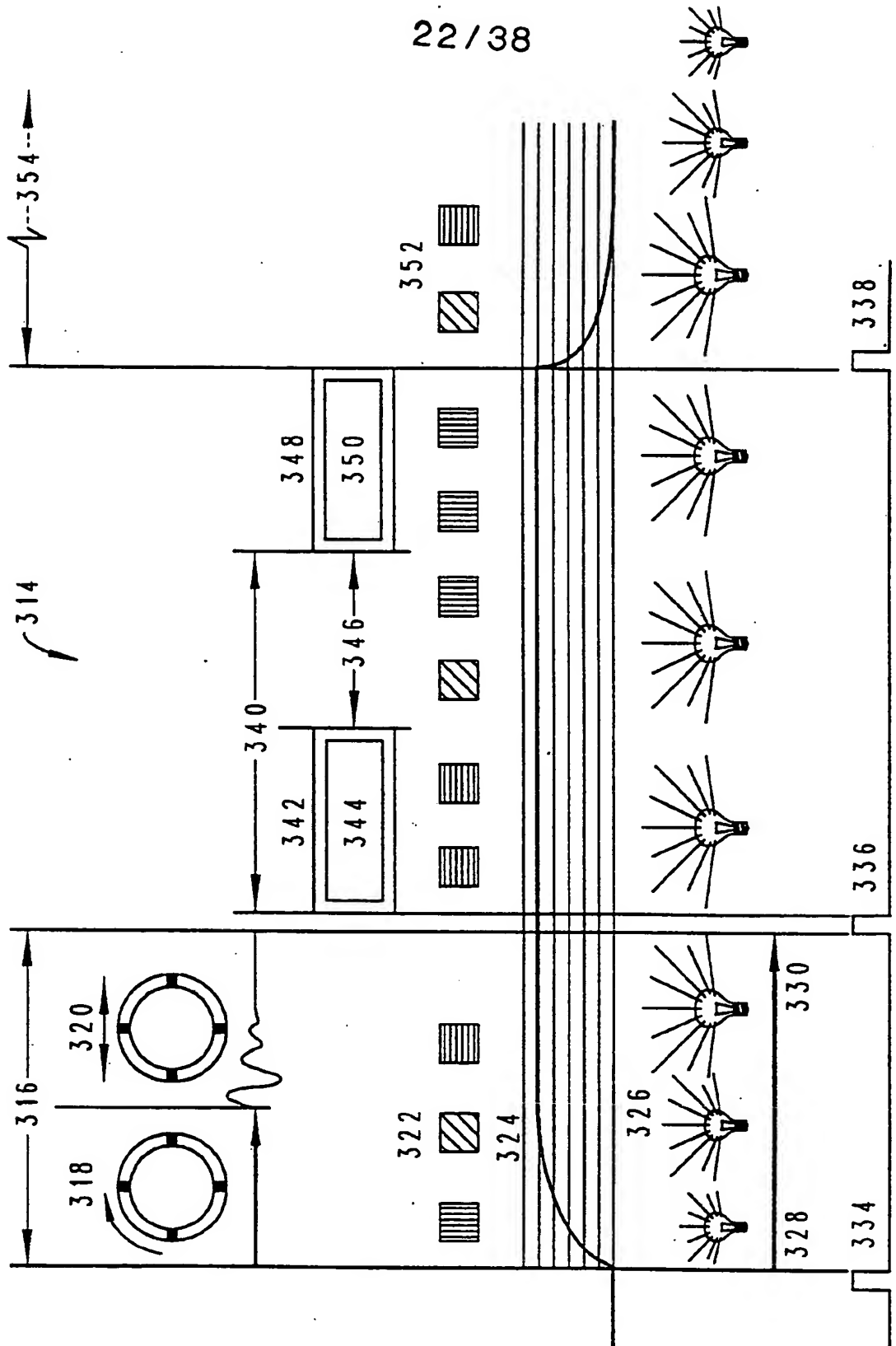


Fig. 22

23/38

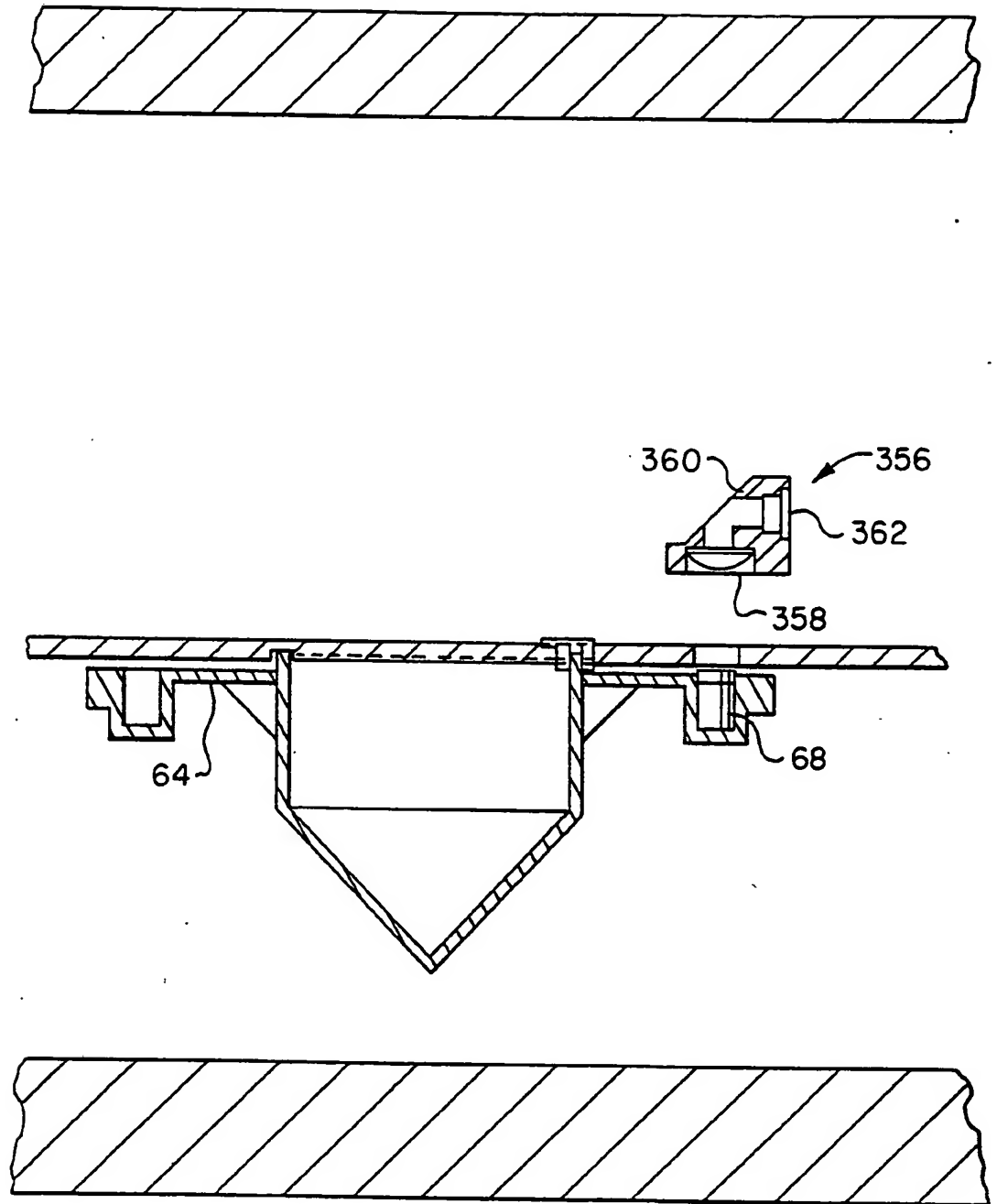
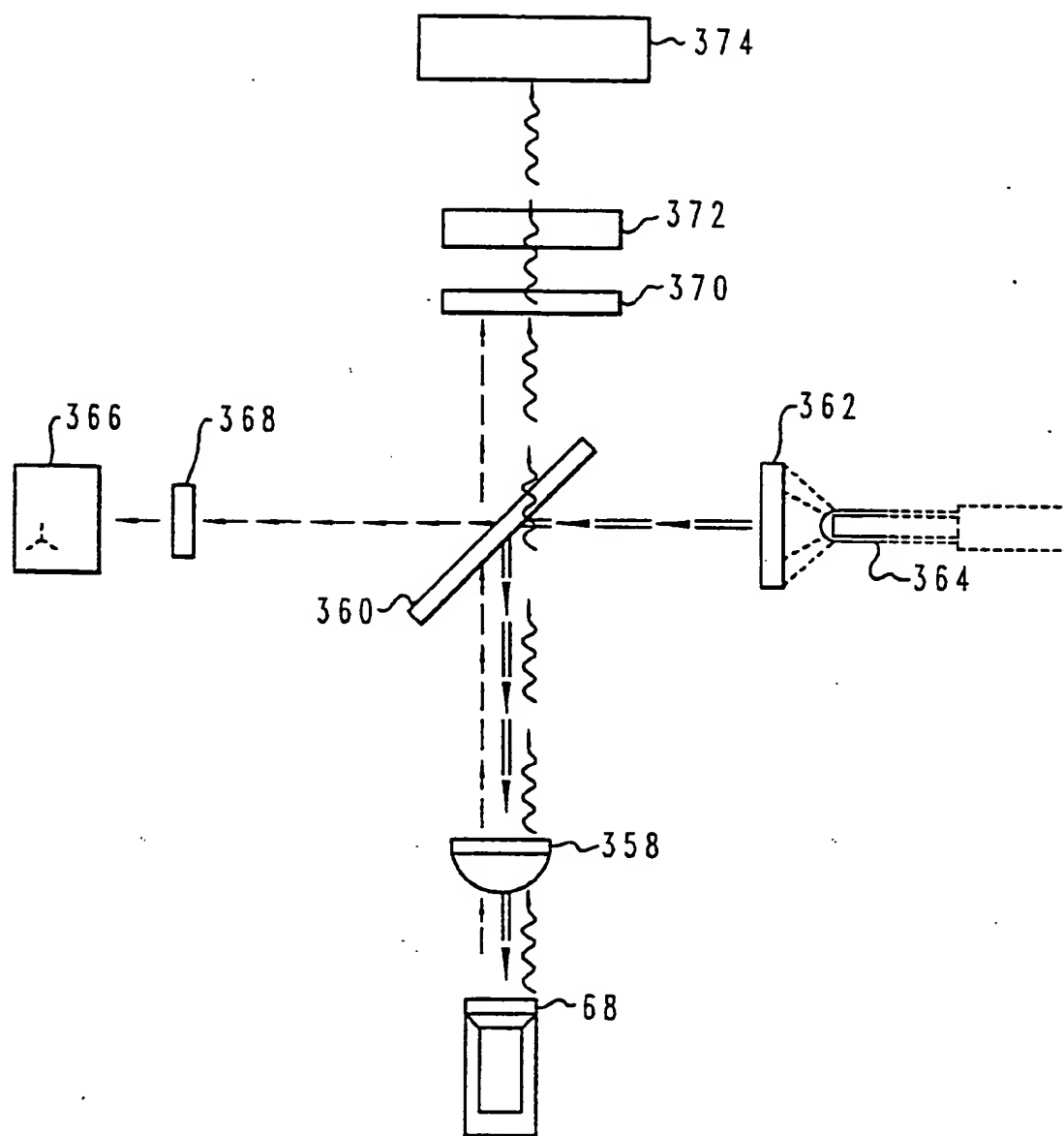


FIG. 23

24/38

*Fig. 24*

25 / 38

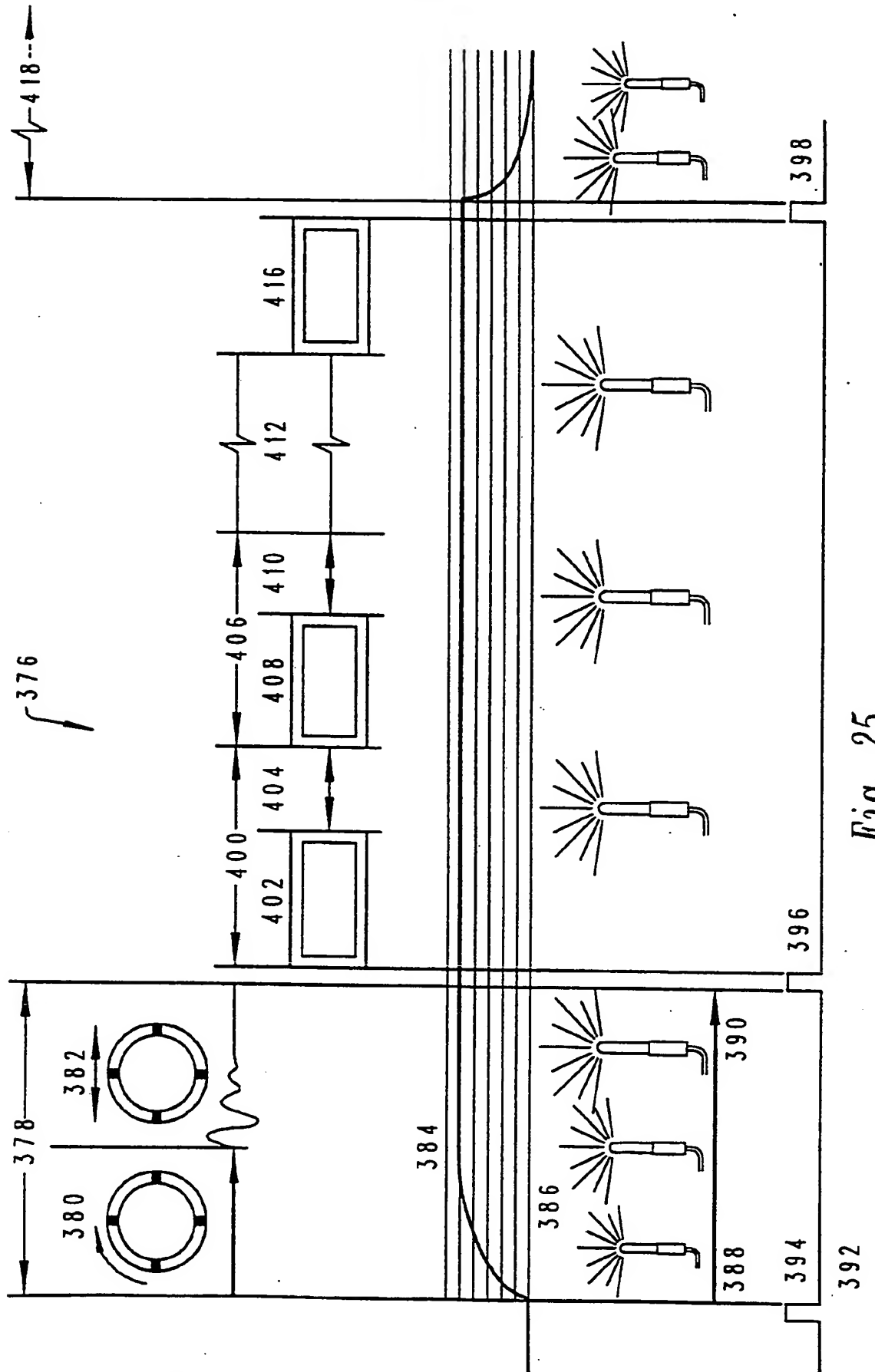


Fig. 25

26/38

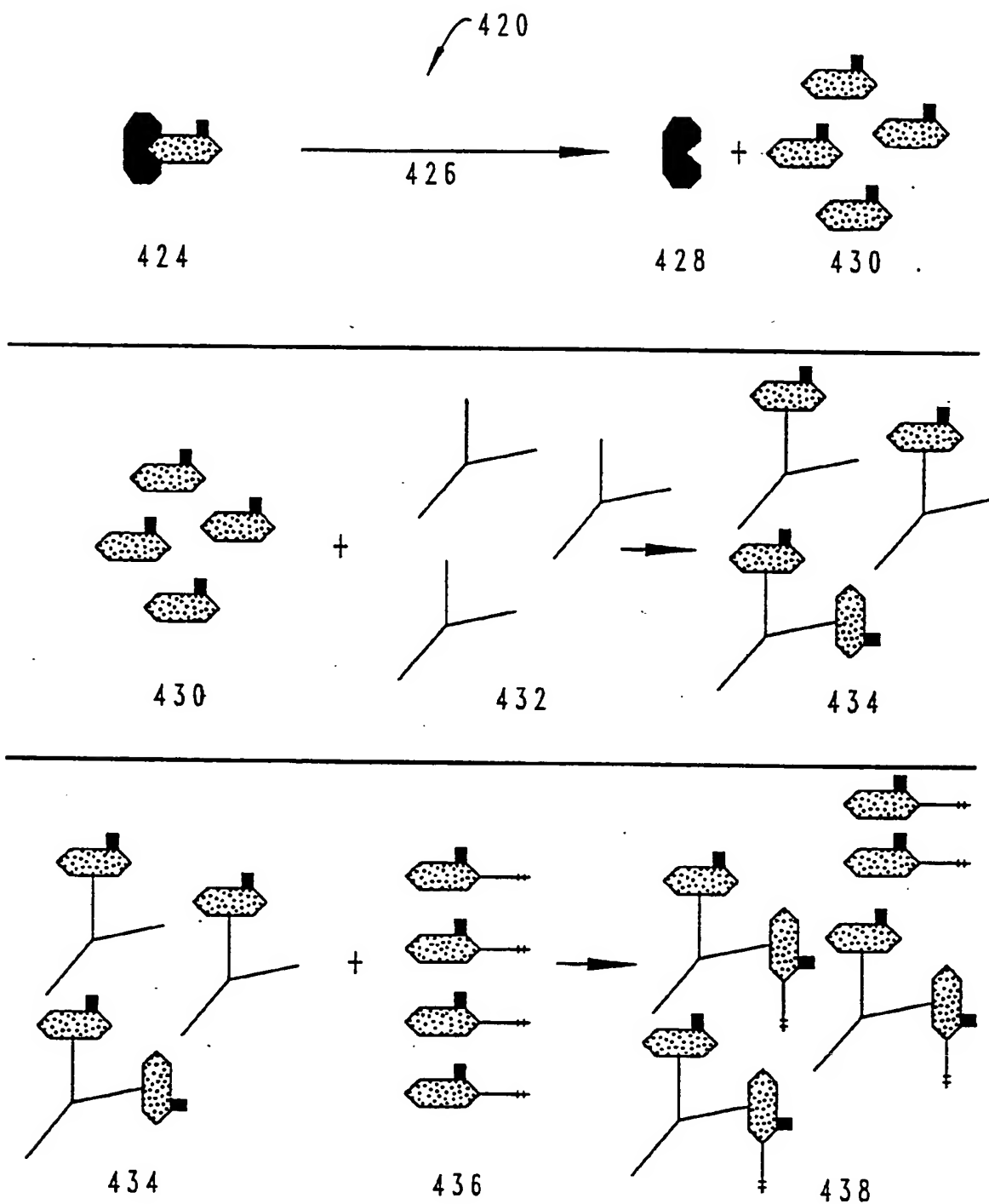


Fig. 26

27/38

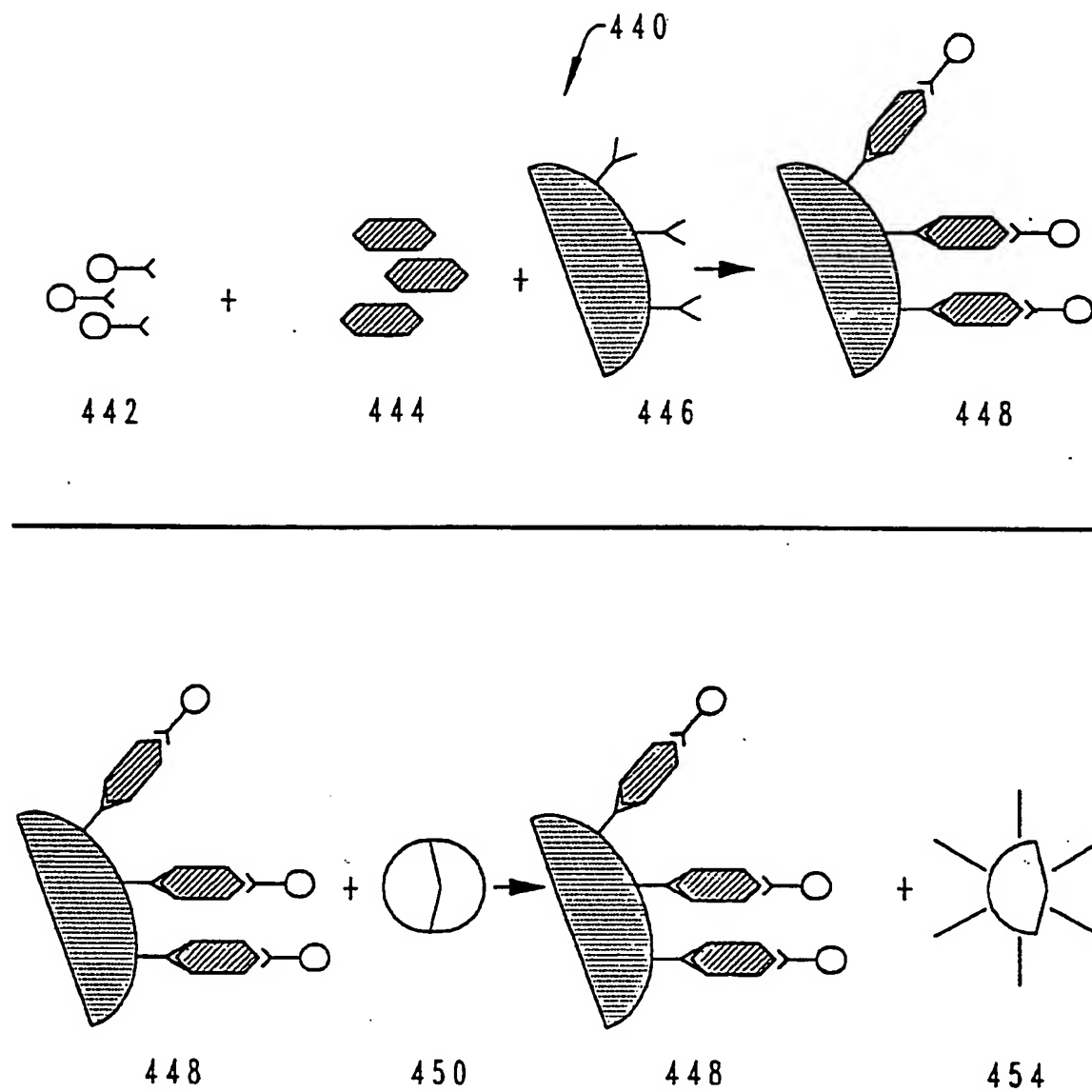


Fig. 27

28 / 38

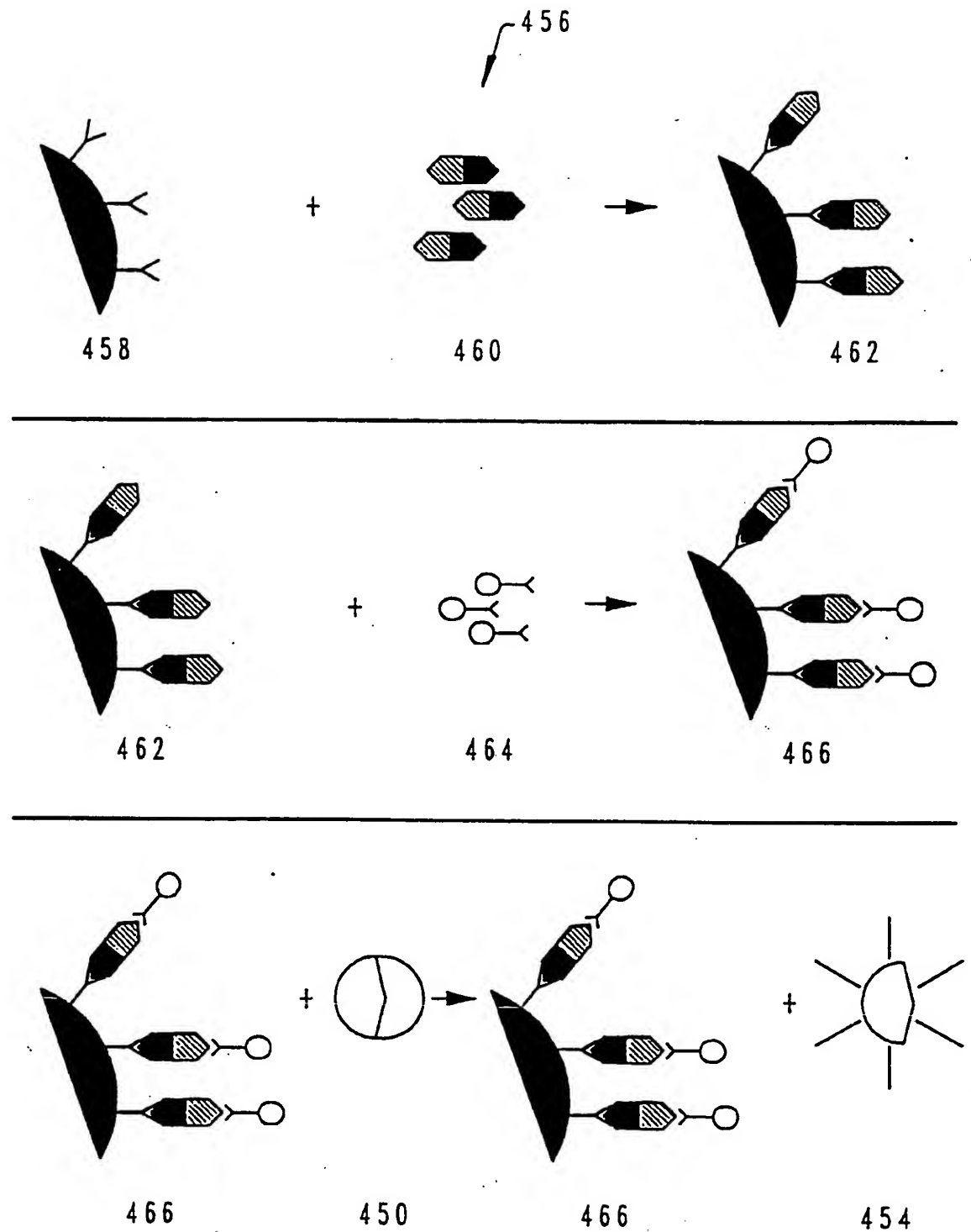
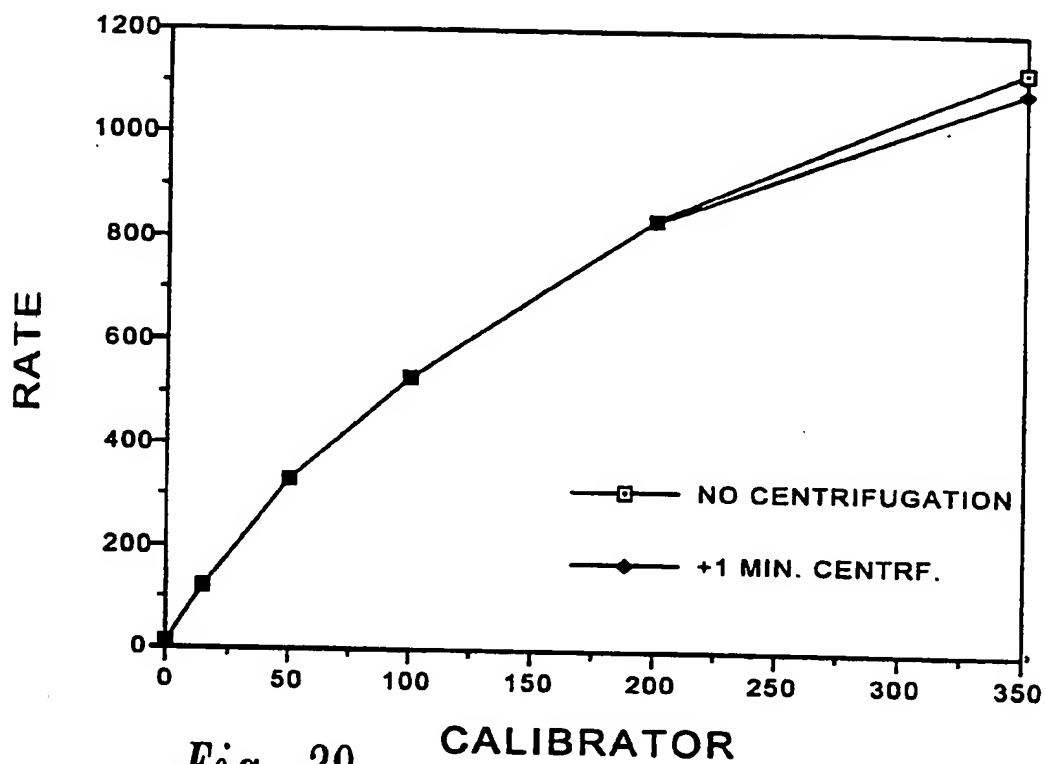



Fig. 28

29/38


*Fig. 29*

CALIBRATOR

30/38

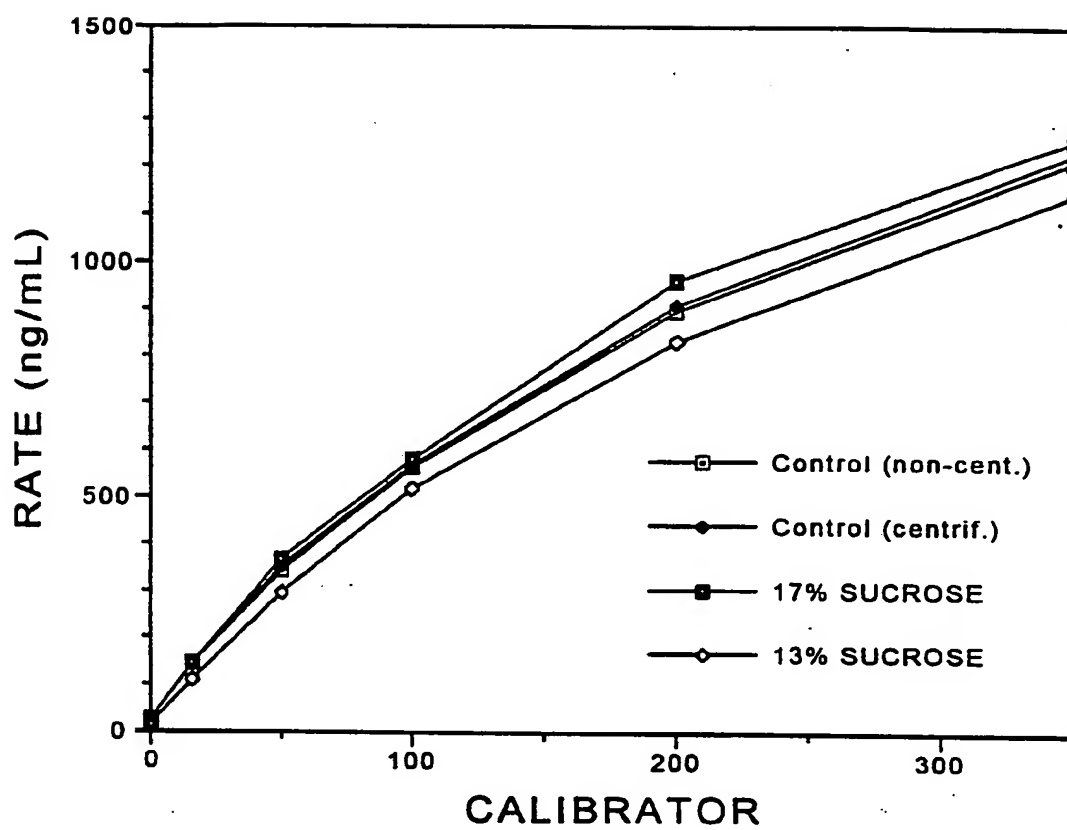*Fig. 30*

31/38

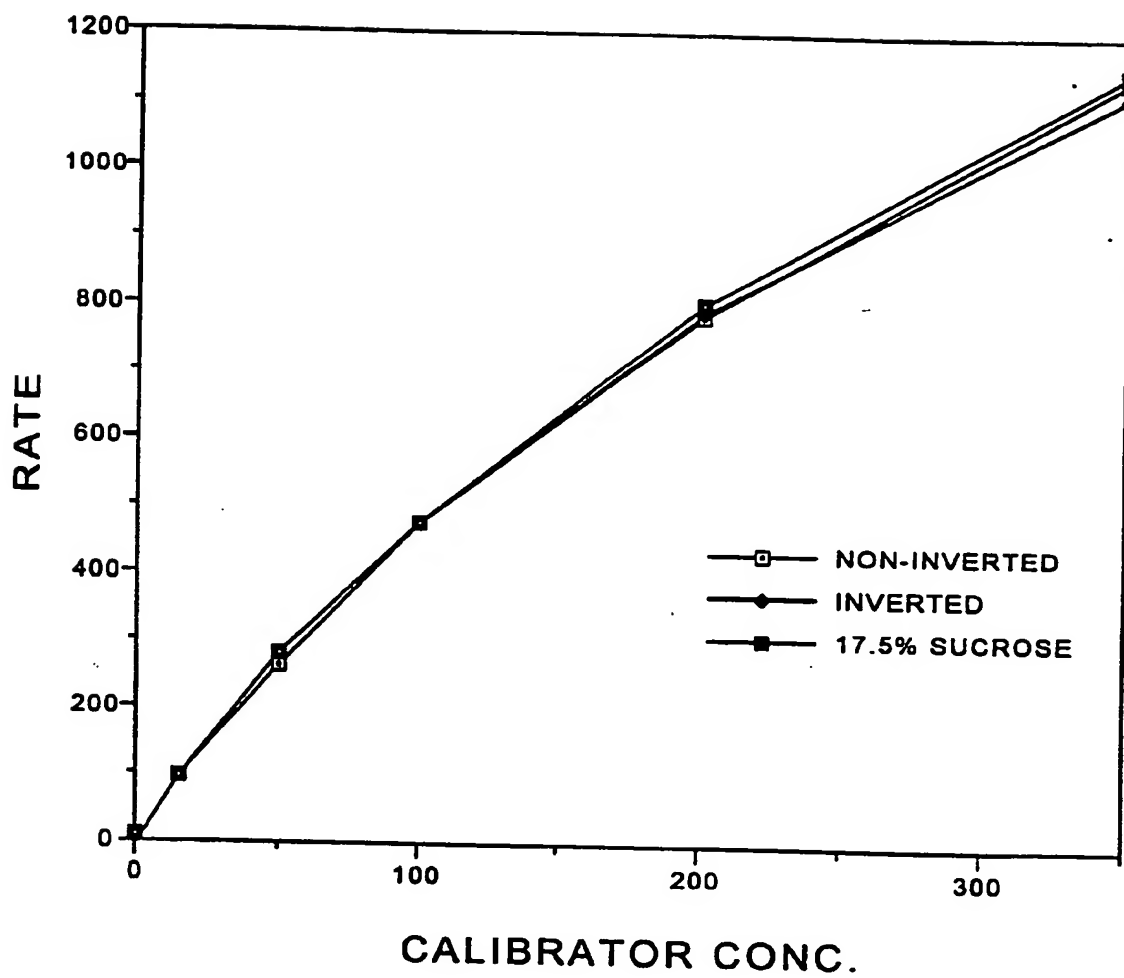*Fig. 31*

32/38

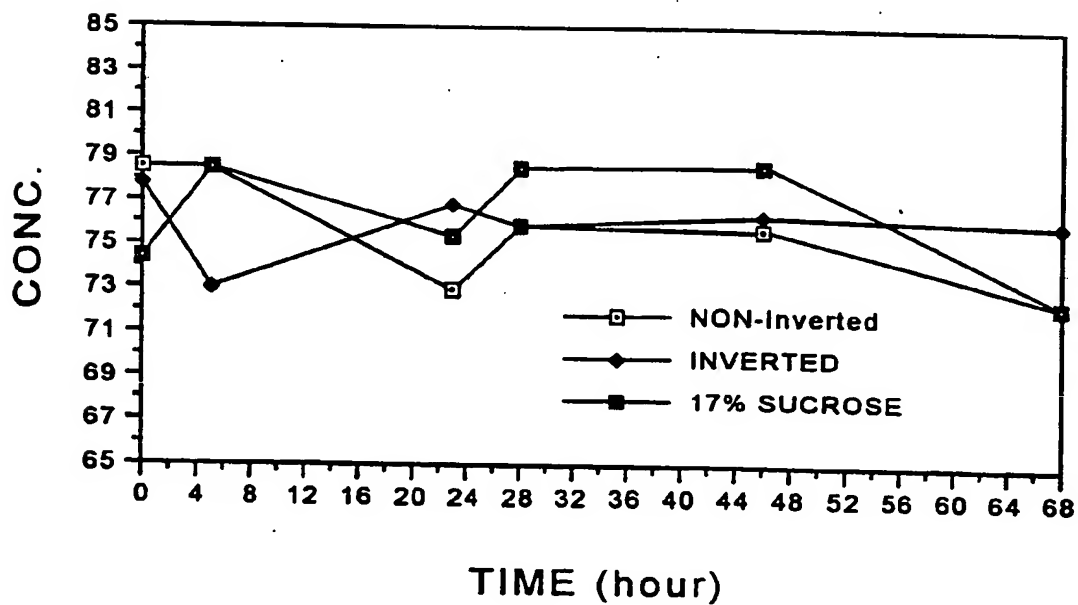
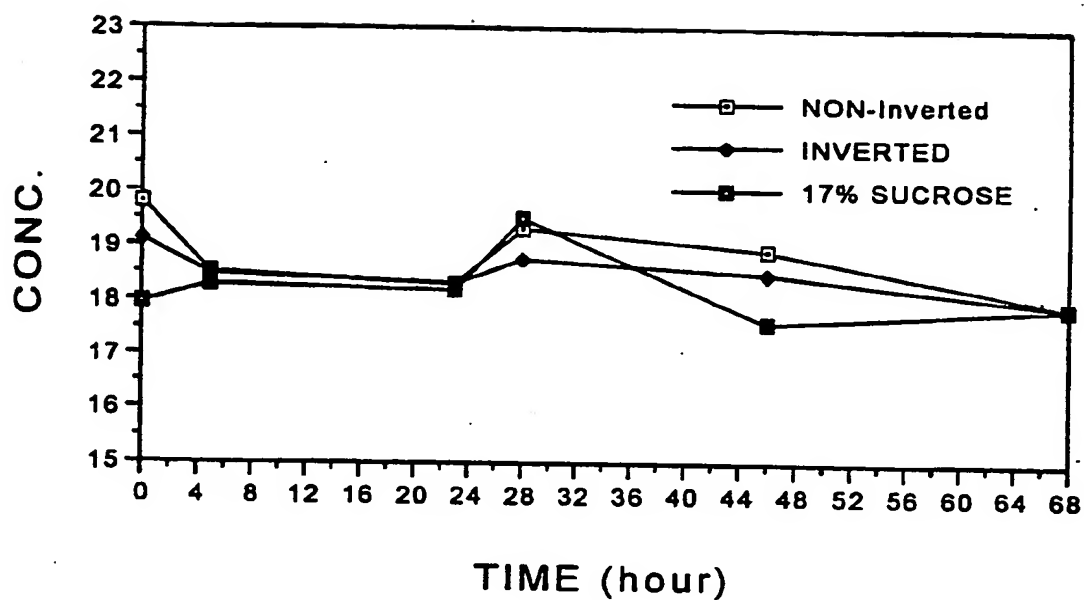
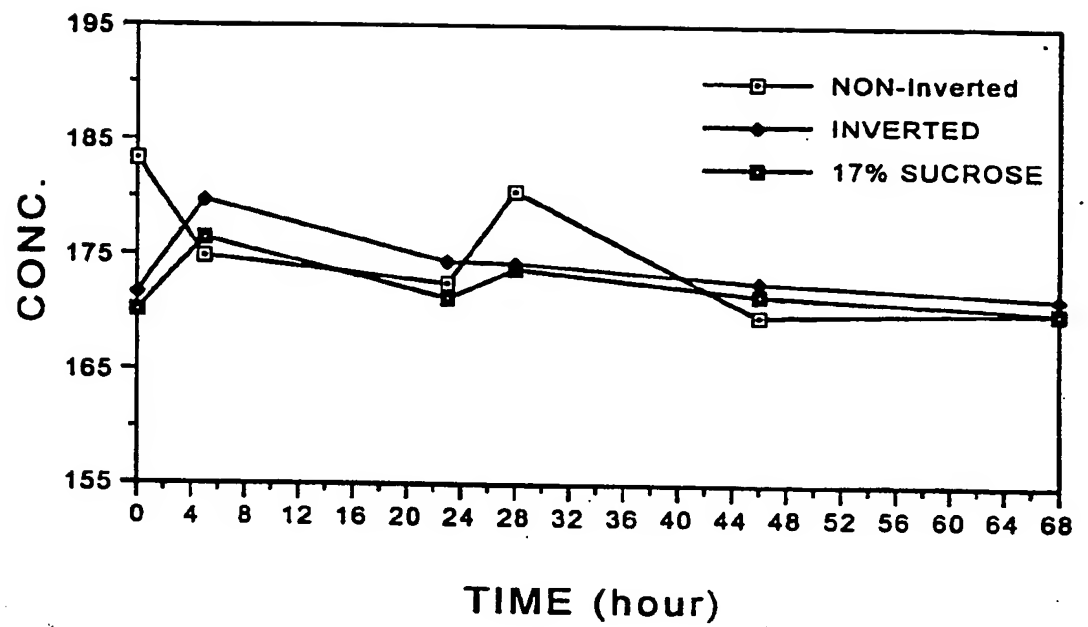
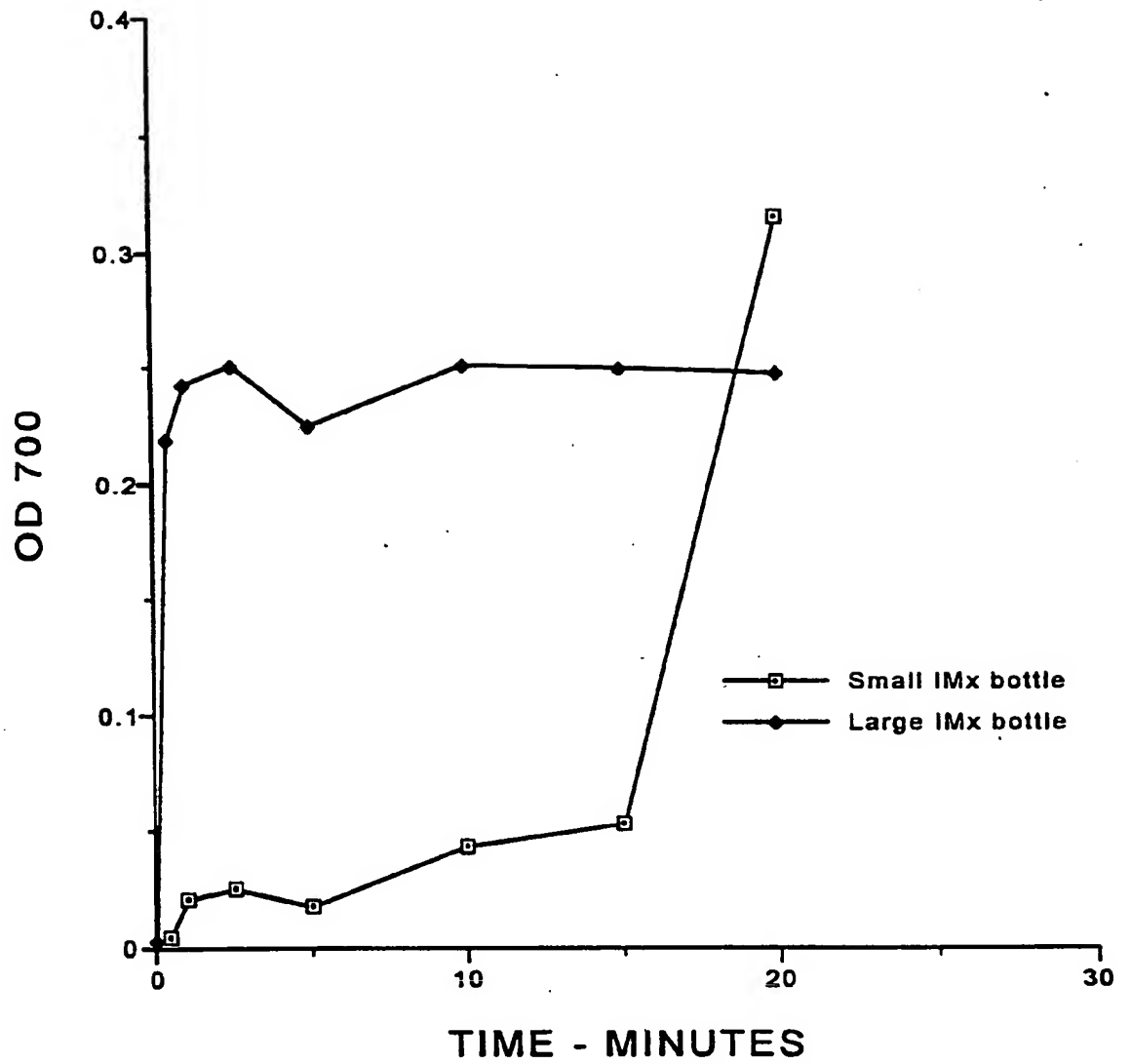


Fig. 32

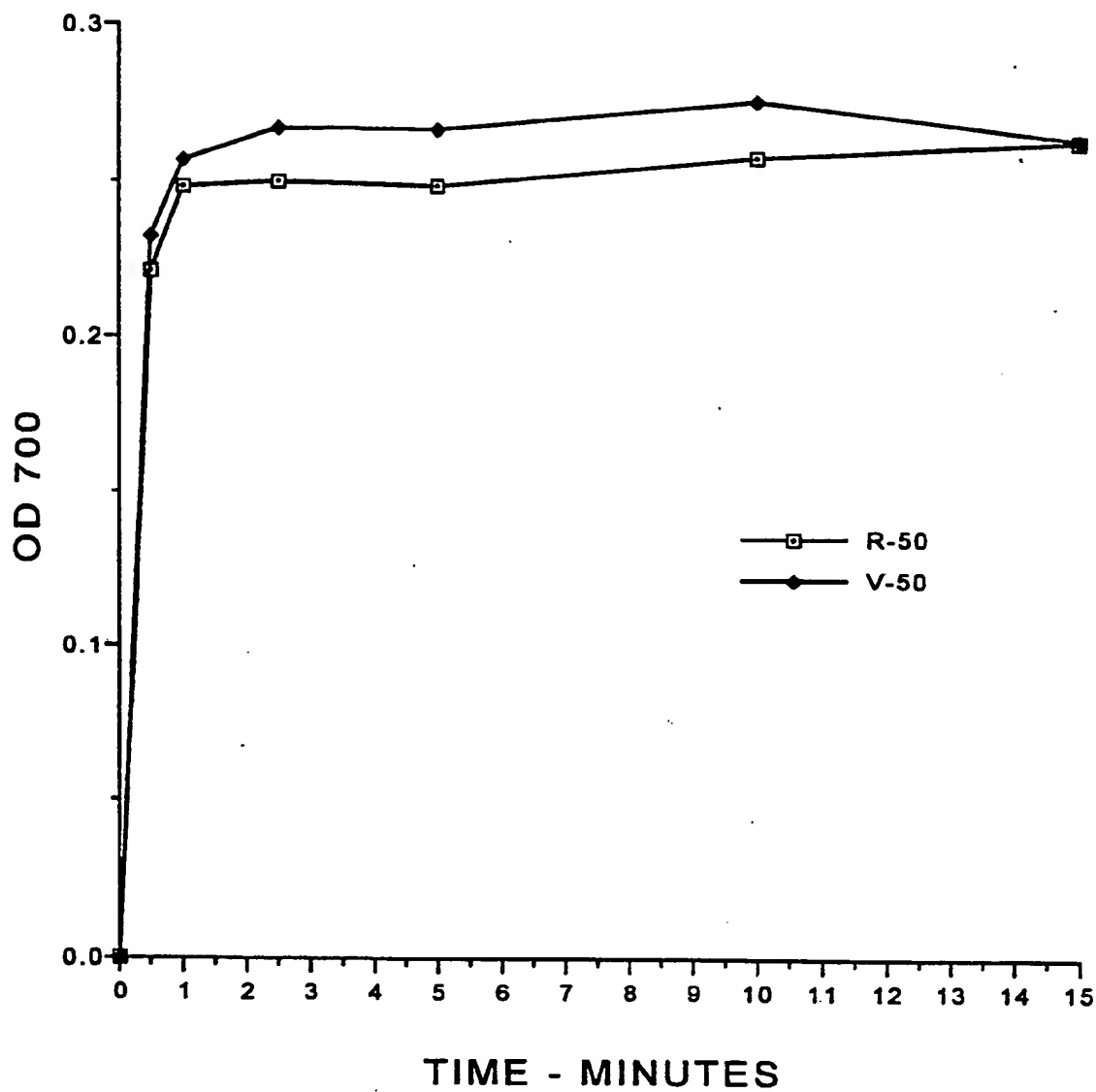
33/38

*Fig. 33*

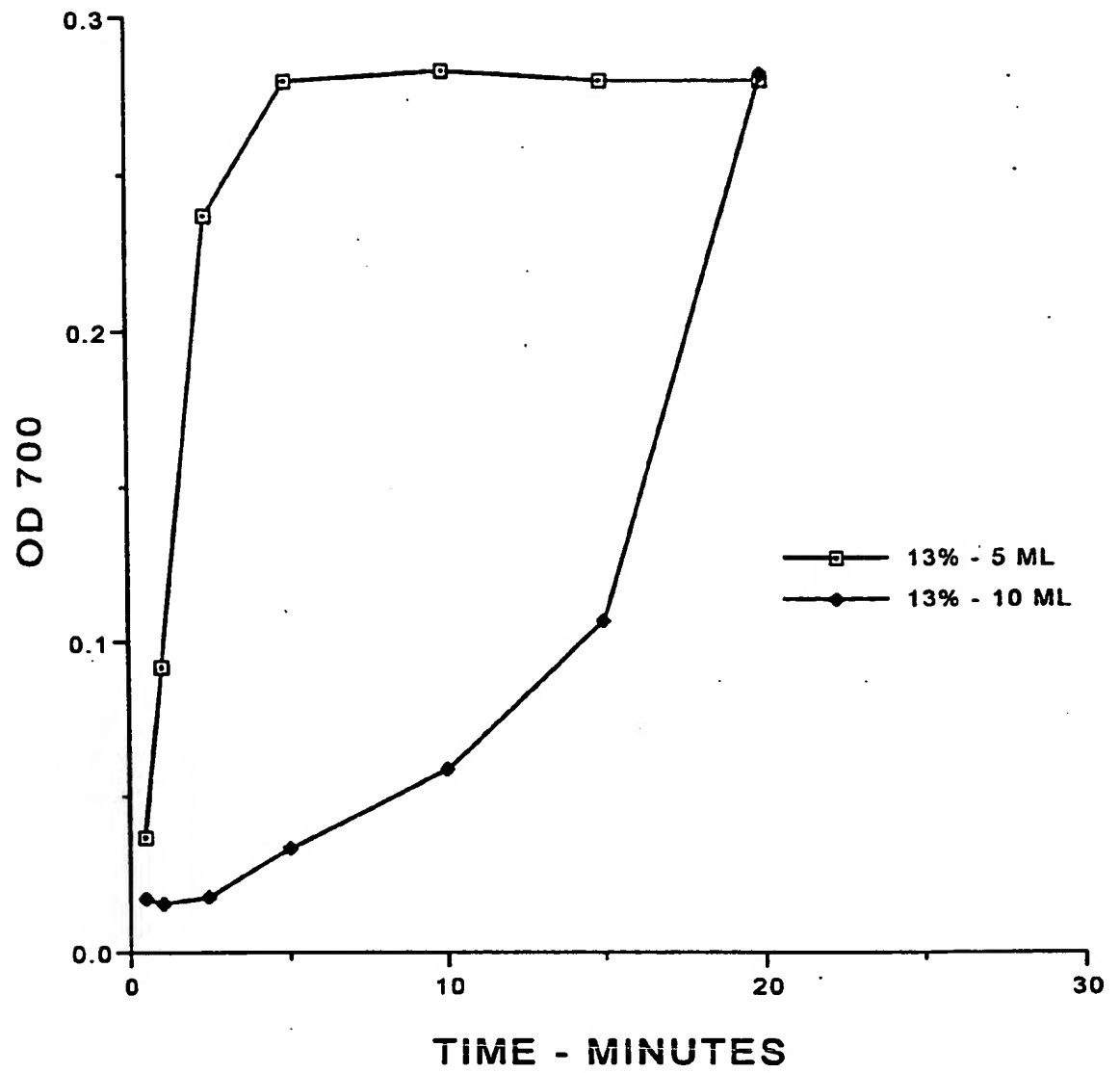
34/38

*Fig. 34*

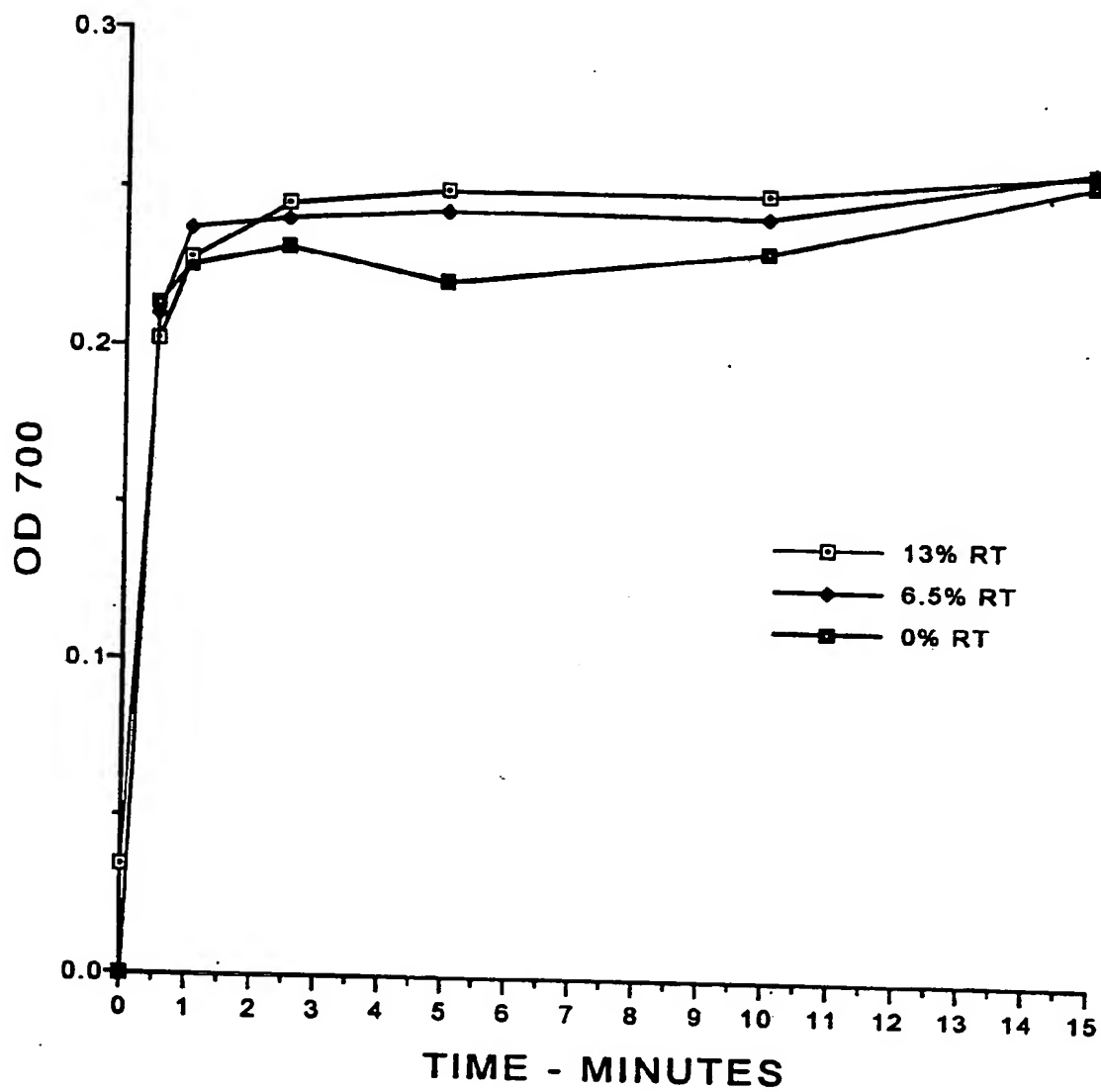
35/38

*Fig. 35*

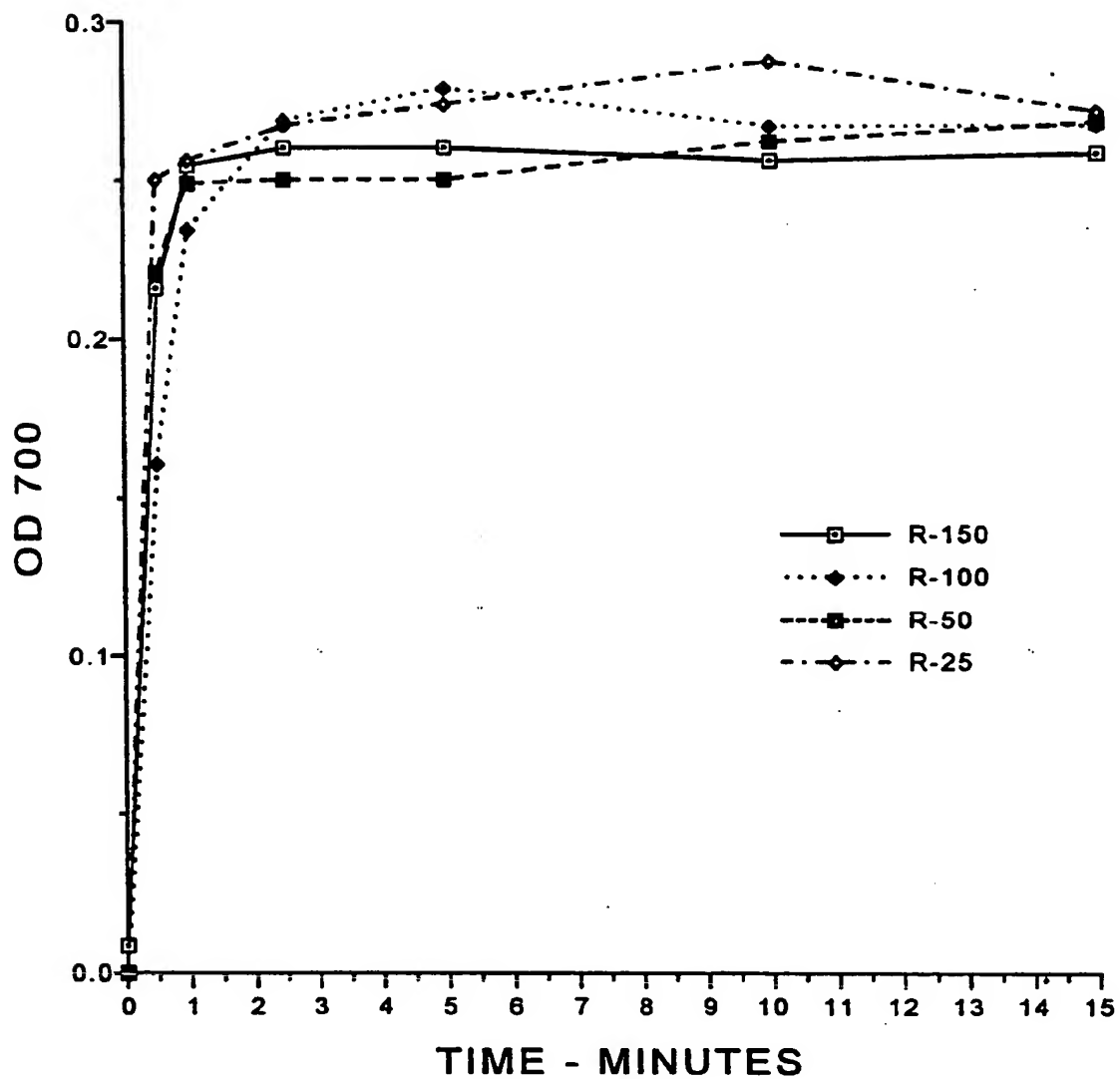
36/38

*Fig. 36*

37/38

*Fig. 37*

38/38

*Fig. 38*

INTERNATIONAL SEARCH REPORT

International application No.
PCT/US93/02776

A. CLASSIFICATION OF SUBJECT MATTER

IPC(5) : G01N 33/546, 35/02

US CL : 436/533, 47

According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

U.S. : 436/533, 534, 47, 50, 826; 422/64, 67

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practicable, search terms used)

APS, JPO

search terms: sucrose, metrizamide or metrioic acid or metroate, particle# or microparticle# or bead#, assay

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X Y	Seradyn, Inc., "Microparticle Immunoassay Techniques", 2nd Edition, Published 1988 by Seradyn, Inc. (Indianapolis, IN), pages 1-49, see page 5, Test Type Overview; page 6, ELISA; page 7, Instrument Based; page 35, Neutral Bouyancy.	<u>1-6,9-15,18-20</u> 7-8,16-17,21-43
X Y	US,A, 4,656,143 (Baker et al.) 07 April 1987, see especially Abstract; column 3, lines 1-30; column 4, lines 8-23; column 5, lines 10-36.	<u>1-5,9-11,13-15,18-20</u> 6-8,12,14,16-17,21-43
Y	L. Pauling, "General Chemistry", published 1953 by W.H. Freeman and Company (San Francisco), page 399, see second paragraph.	7-8,16-17,21-43

☒ Further documents are listed in the continuation of Box C. ☐ See patent family annex.

* Special categories of cited documents:	T	later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
"A" document defining the general state of the art which is not considered to be part of particular relevance	X	document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
"E" earlier document published on or after the international filing date	Y	document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art
"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)	A	document member of the same patent family
"O" document referring to an oral disclosure, use, exhibition or other means		
"P" document published prior to the international filing date but later than the priority date claimed		

Date of the actual completion of the international search

14 June 1993

Date of mailing of the international search report

23 JUN 1993

Name and mailing address of the ISA/US
Commissioner of Patents and Trademarks
Box PCT
Washington, D.C. 20231

Authorized officer

SUSAN C. WOLSKI

Facsimile No. NOT APPLICABLE

Telephone No. (703) 308-0196

Form PCT/ISA/210 (sec 3/21/2005, EAST Version: 2.0.1.4)

INTERNATIONAL SEARCH REPORT

International application No.
PCT/US93/02776

C (Continuation). DOCUMENTS CONSIDERED TO BE RELEVANT

Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
Y	EP,A, 0,410, 645 (Cohen et al.) 30 January 1991, see especially Abstract; page 4, lines 37-51; page 6, lines 36-47; page 10, lines 44-55; page 11, lines 16-22.	23-43

Form PCT/ISA/210 (continuation of second sheet)(July 1992)*

INTERNATIONAL SEARCH REPORT

International application No.
PCT/US93/02716

Box I Observations where certain claims were found unsearchable (Continuation of item 1 of first sheet)

This international report has not been established in respect of certain claims under Article 17(2)(a) for the following reasons:

1. ☐ Claims Nos.:
because they relate to subject matter not required to be searched by this Authority, namely:

2. ☐ Claims Nos.:
because they relate to parts of the international application that do not comply with the prescribed requirements to such an extent that no meaningful international search can be carried out, specifically:

3. ☐ Claims Nos.:
because they are dependent claims and are not drafted in accordance with the second and third sentences of Rule 6.4(a).

Box II Observations where unity of invention is lacking (Continuation of item 2 of first sheet)

This International Searching Authority found multiple inventions in this international application, as follows:
(Telephone Practice)

- I. Claims 1-22, drawn to a method of providing a liquid assay reagent and a method of performing a heterogeneous immunoassay, classified in Class 435, subclass 533.
 - II. Claims 23-43, drawn to a method of operating an automated analytical system, classified in Class 435, subclass 47.
- The claims of these two groups are drawn to distinct inventions which are not so linked as to form a single general inventive concept. PCT Rules 13.1 and 13.2 do not provide for multiple distinct methods.

1. ☒ As all required additional search fees were timely paid by the applicant, this international search report covers all searchable claims.
2. ☐ As all searchable claims could be searched without effort justifying an additional fee, this Authority did not invite payment of any additional fee.
3. ☐ As only some of the required additional search fees were timely paid by the applicant, this international search report covers only those claims for which fees were paid, specifically claims Nos.:

4. ☐ No required additional search fees were timely paid by the applicant. Consequently, this international search report is restricted to the invention first mentioned in the claims; it is covered by claims Nos.:

Remark on Protest

- ☐ The additional search fees were accompanied by the applicant's protest.
☒ No protest accompanied the payment of additional search fees.